



Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Ifenprodil (Vasodilator)
 イフェンプロジル (末梢血管拡張薬)

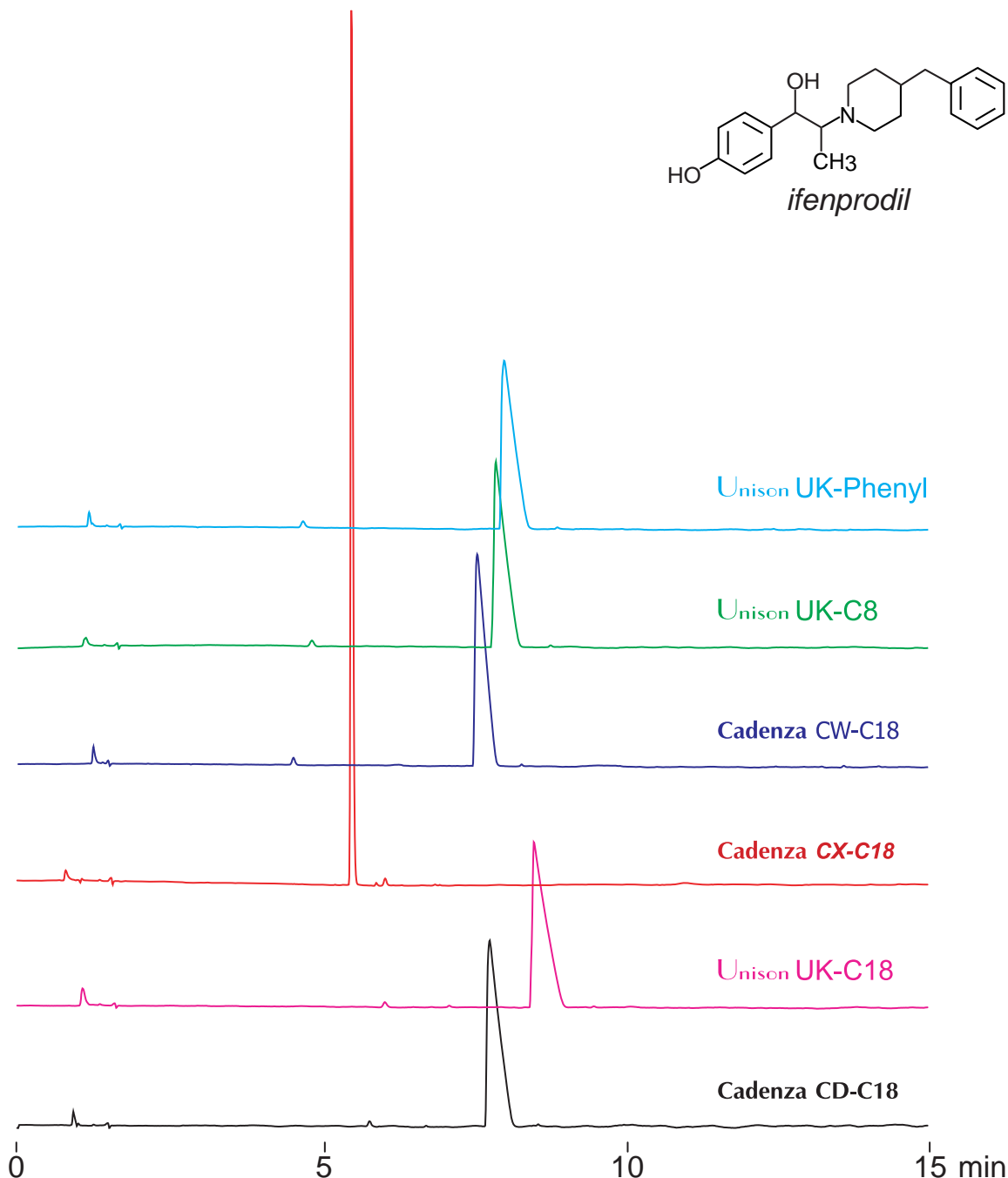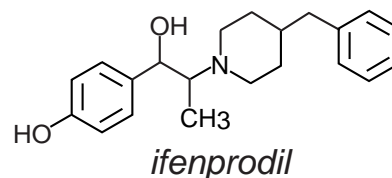

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

10-100%B (0-15min)

1.0mL/min (12MPa), 37°C, 280nm

10µL (1.0mg/mL, water / formic acid = 100 / 1)