

Unison UK-C8
Unison UK-C18



Unison UK-Phenyl
Cadenza CD-C18

150 x 4.6 mm

HPLC COLUMN APPLICATION

Chloronicotinic acid isomers

クロロニコチン酸異性体

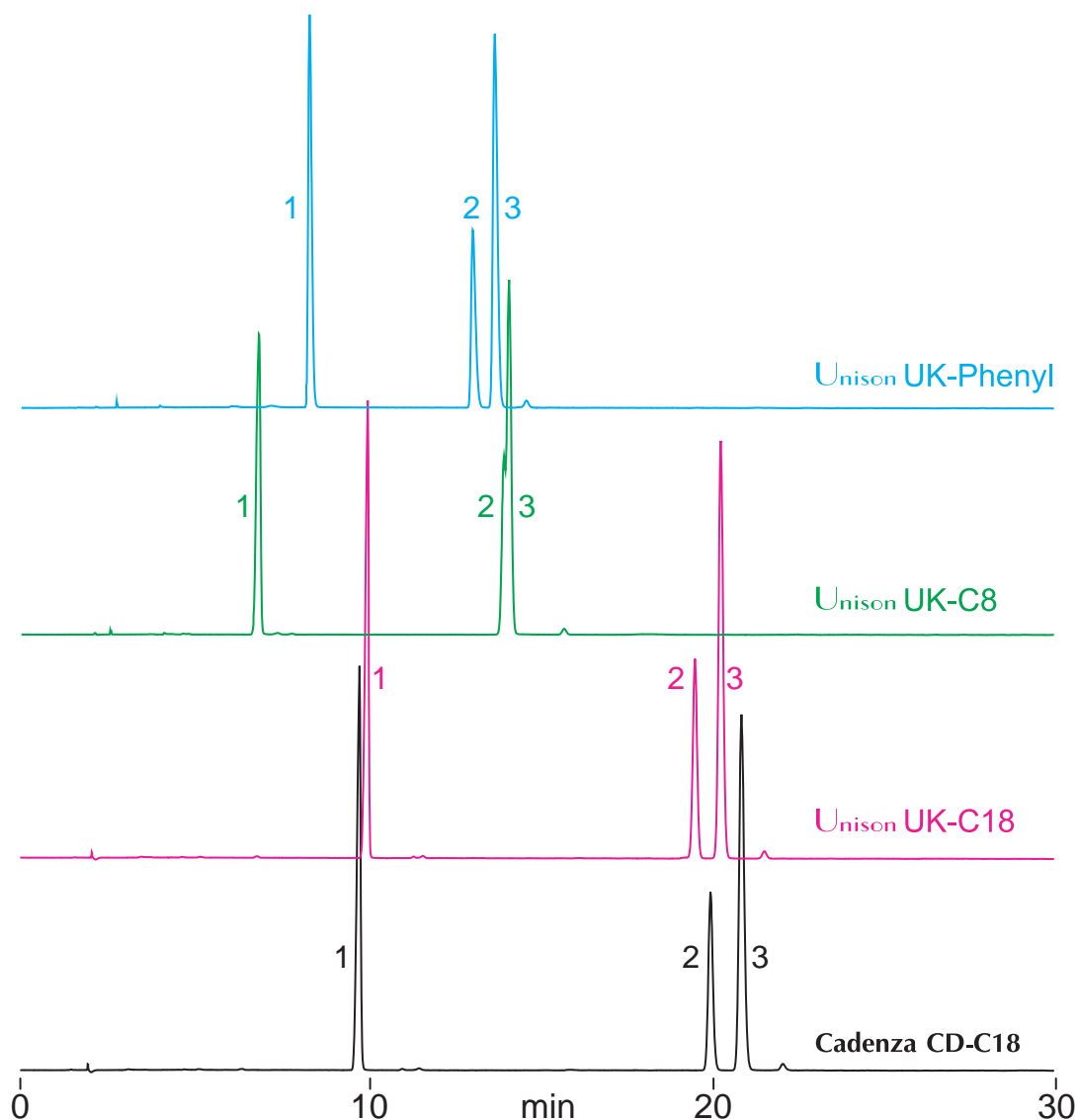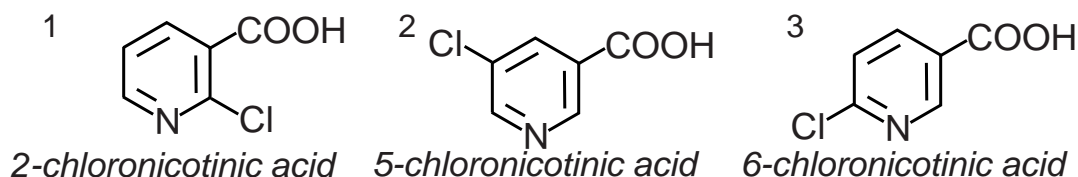

150 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-30%B (0-30min)

1.0mL/min (11MPa), 37°C, 260nm

10µL (0.3mg/mL, water / methanol / formic acid = 50 / 50 / 1)