



Unison UK-C8

Unison UK-Phenyl

100 x 4.6 mm

HPLC COLUMN APPLICATION

Monoterpeneketone
モノテルペンケトン

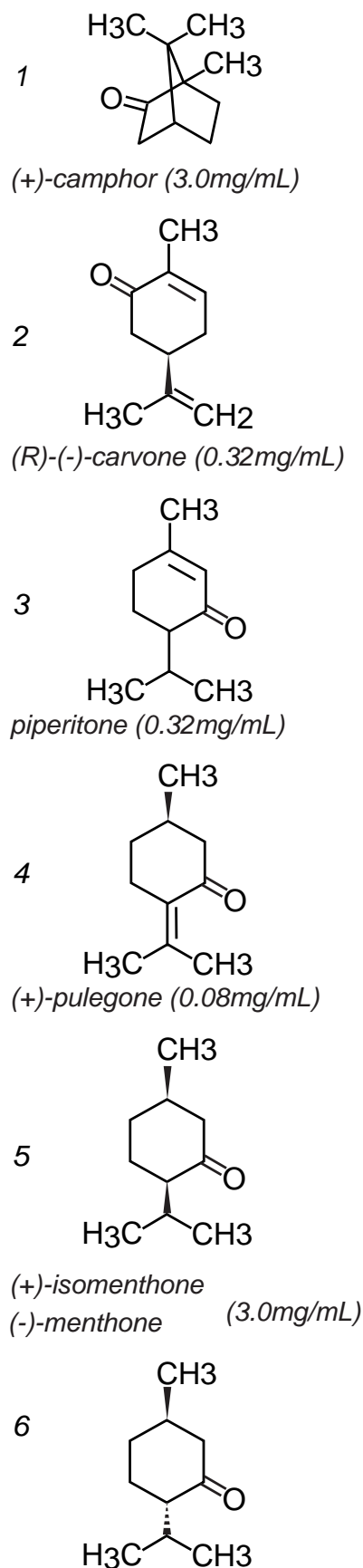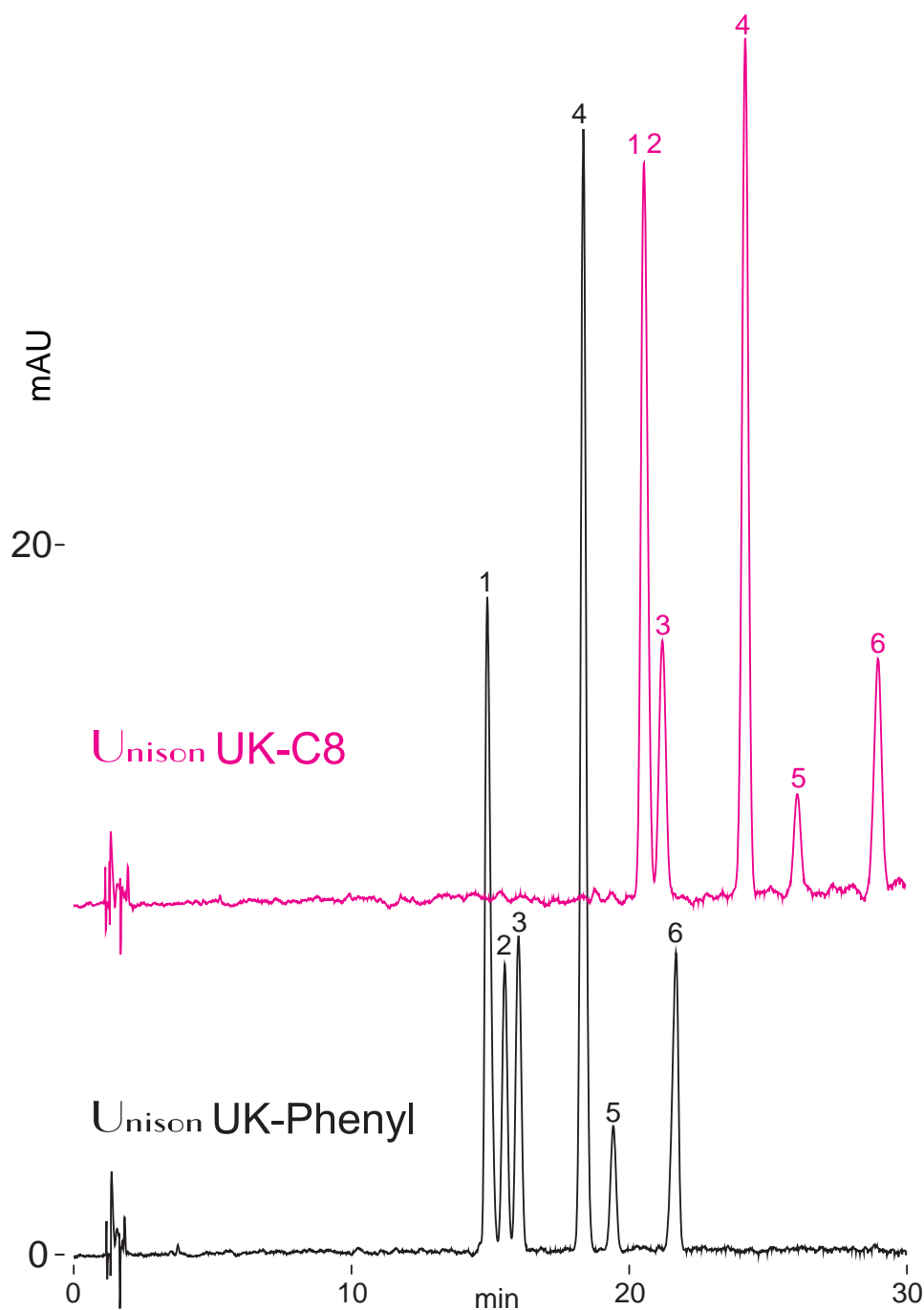

100 x 4.6 mm

A: 100mMHCOONH₄ / ACN = 90 / 10

B: 100mMHCOONH₄ / ACN = 50 / 50

0-75%B(0-30min)

1.0mL/min(9MPa), 280nm, 60°C

20uL (0.08-3.0mg/mL, water / MeOH / HCOOH = 50 / 50 / 1)