

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Bisacodyl ( Laxative )

ビサコジル ( 便秘薬 )

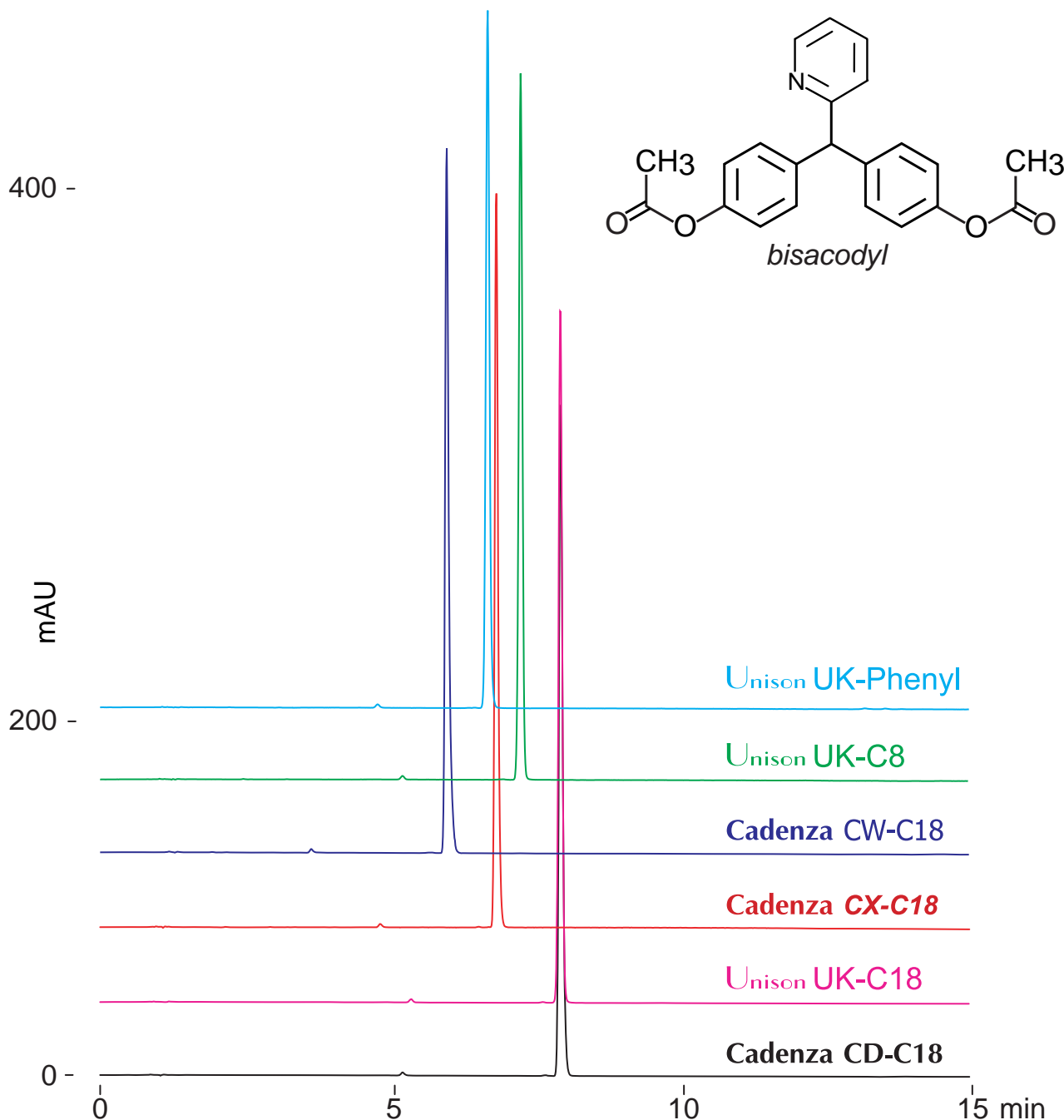

100 x 4.6 mm

A: 100mM HCOONH<sub>4</sub>

B: 100mM HCOONH<sub>4</sub> / acetonitrile = 10 / 90

40-100%B (0-15min)

1.0mL/min (13MPa), 37°C, 260nm

10μL (0.25mg/mL, 100mM HCOONH<sub>4</sub>/ acetonitrile= 50 / 50)