
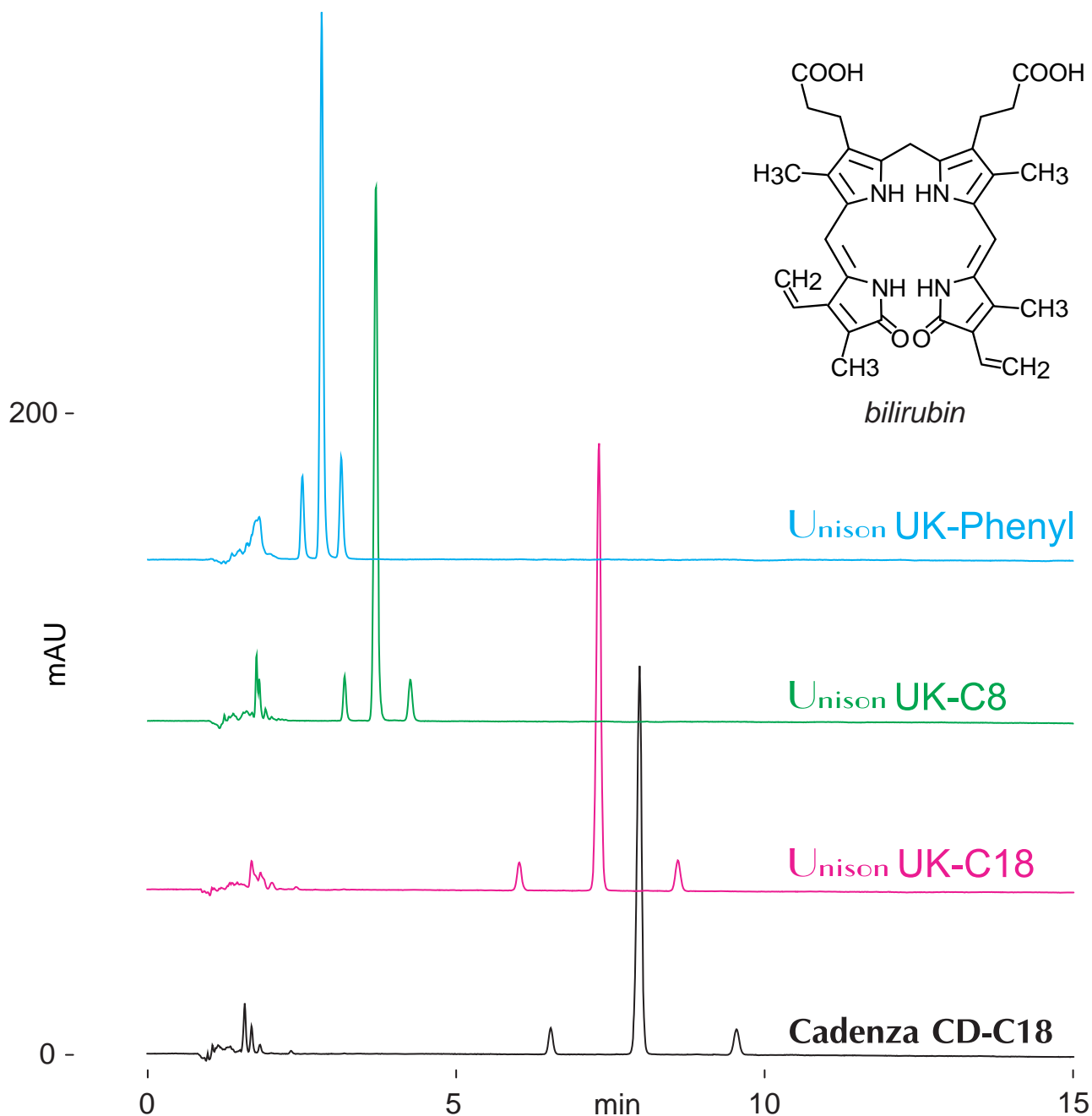


Unison UK-C18 Unison UK-Phenyl  
Cadenza CD-C18 Unison UK-C8 100 x 4.6 mm

HPLC COLUMN APPLICATION

**Bilirubin (Heme metabolite)**  
ビリルビン (ヘム代謝物)

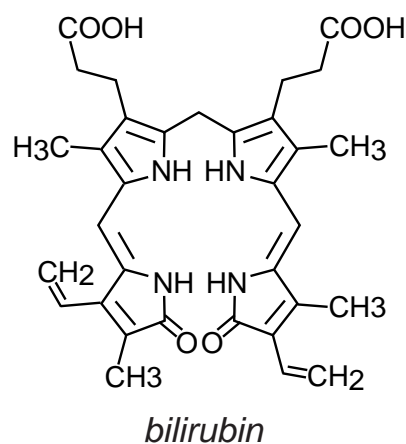

100 x 4.6 mm  
 A: 100mM ammonium acetate  
 B: 100mM ammonium acetate / acetonitrile = 10 / 90  
 40-100%B (0-13min)  
 1.0mL/min (13MPa), 37°C, 450nm  
 10µL (0.025mg/mL, DMSO / HCOOH = 100 / 1)