

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Desloratadine (Histamine receptor antagonist)

デスロラタジン (ヒスタミン受容体拮抗薬)

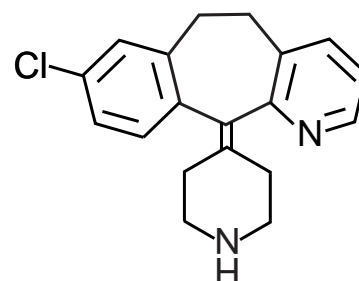

desloratadine

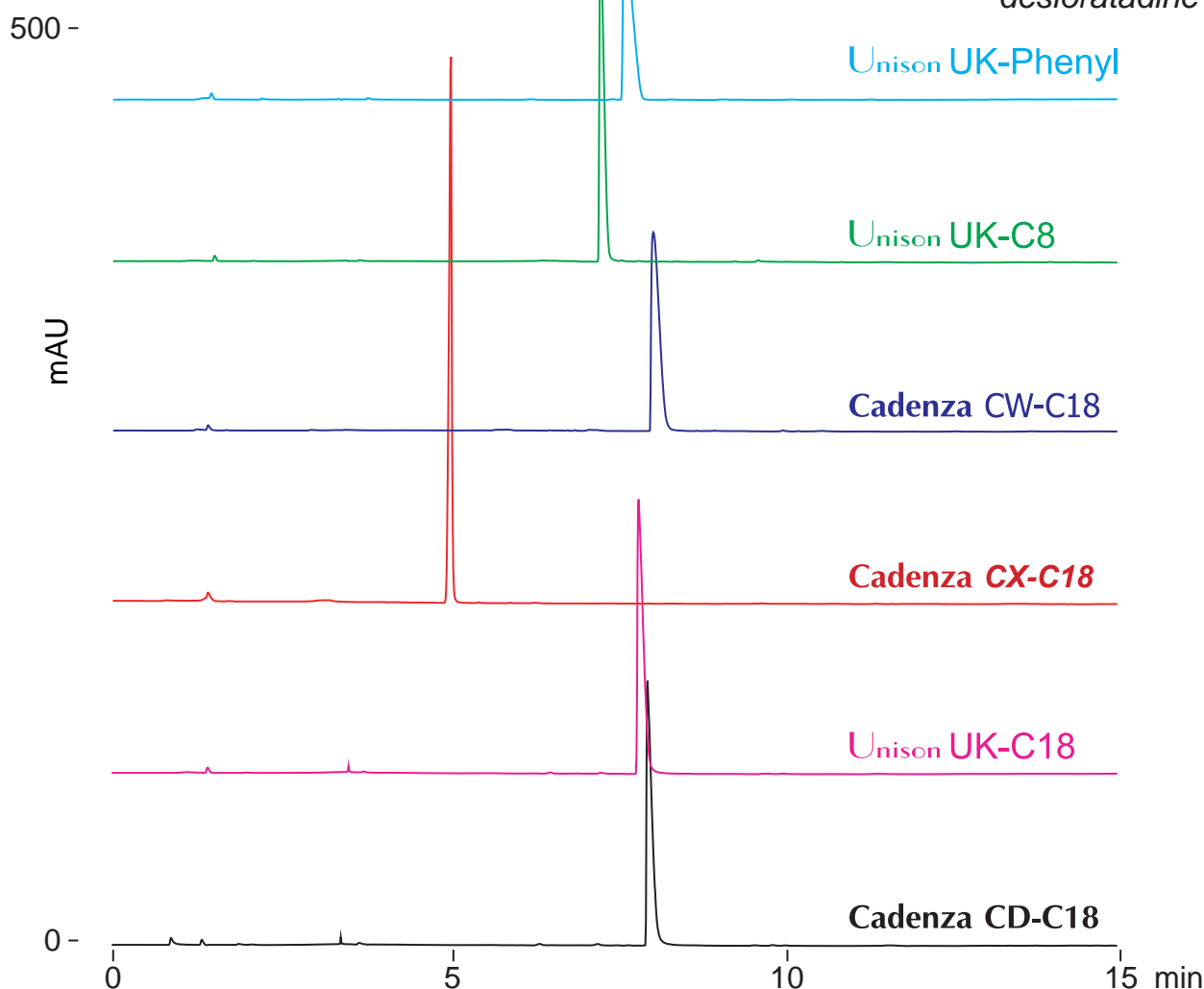

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-80%B (0-15min)

1.0mL/min (14MPa), 37°C, 260nm

10µL (0.05mg/mL, water / formic acid = 100 / 1)