

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18

Cadenza CW-C18


Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

3,3'-Diaminobenzidine

3,3'-ジアミノベンジジン

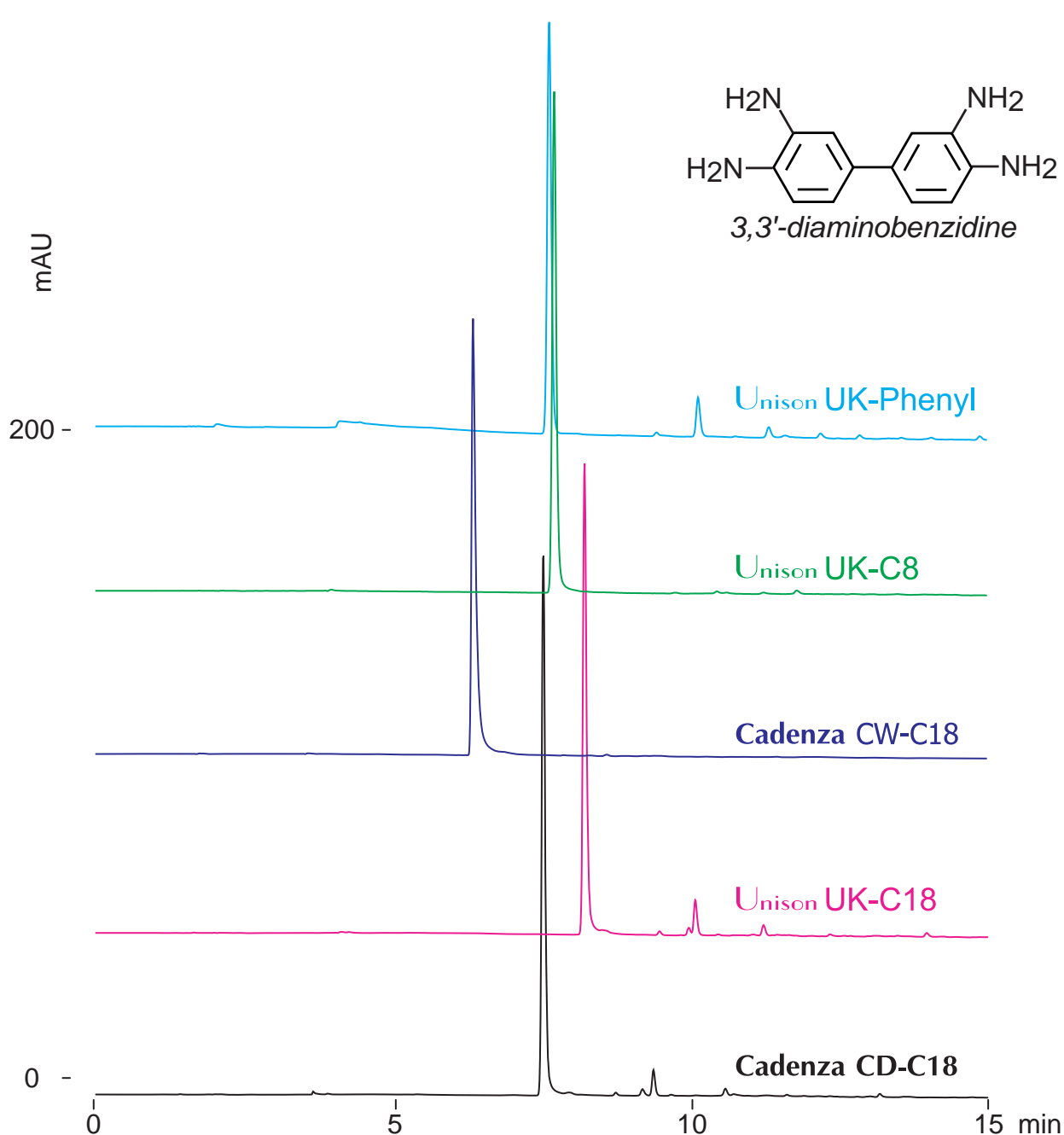

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-40%B (0-15min)

1.0mL/min (14MPa), 37°C, 260nm

10µL (0.05mg/mL, 100mM ammonium formate / acetonitrile = 95 / 5)