

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18


Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Loxoprofen (Nonsteroidal anti-inflammatory drug)

ロキソプロフェン (非ステロイド性抗炎症薬)

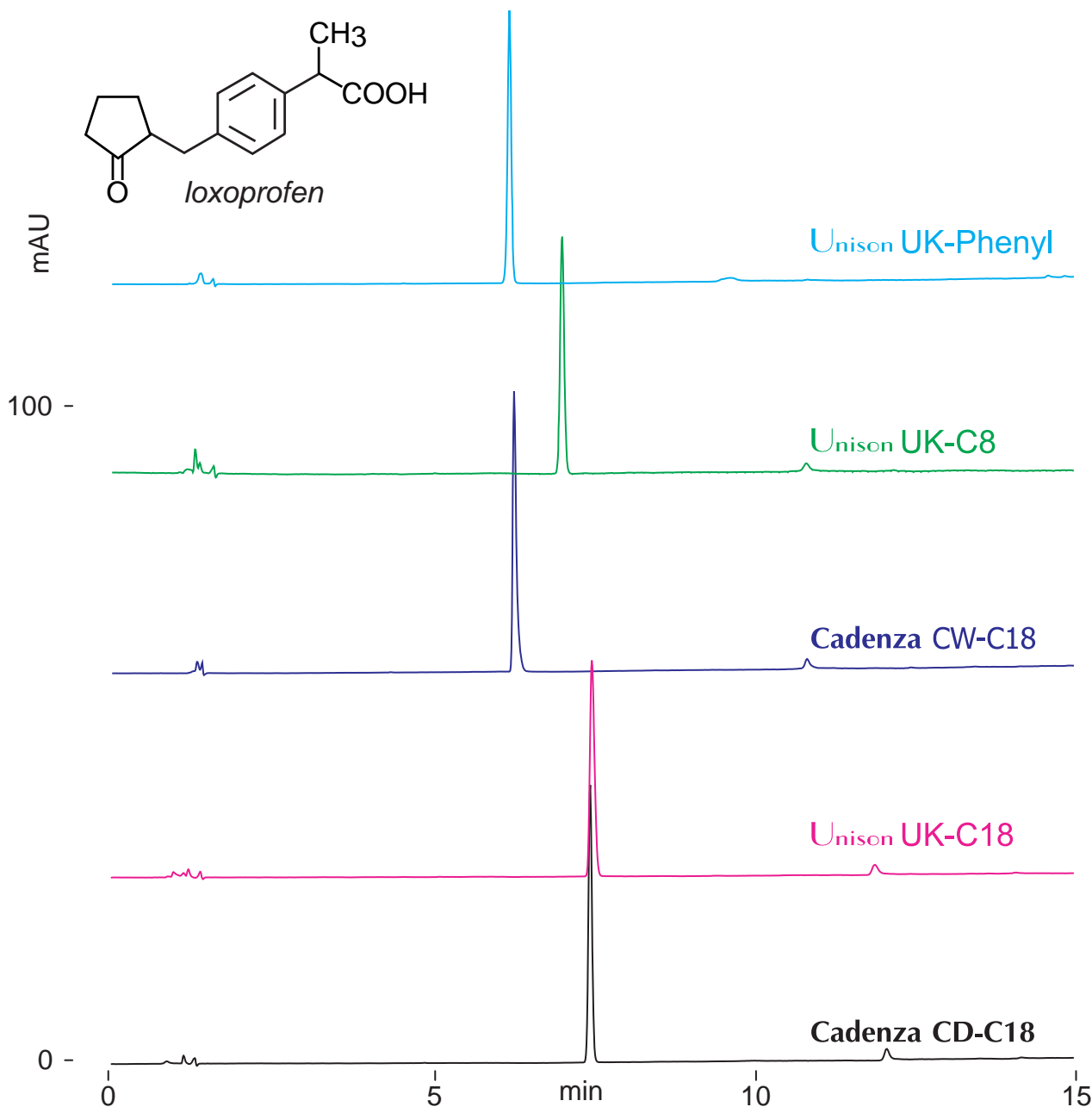

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

25-100%B (0-15min)

1.0mL/min (14MPa), 37°C, 260nm

10µL (1.0mg/mL, water / formic acid = 100 / 1)