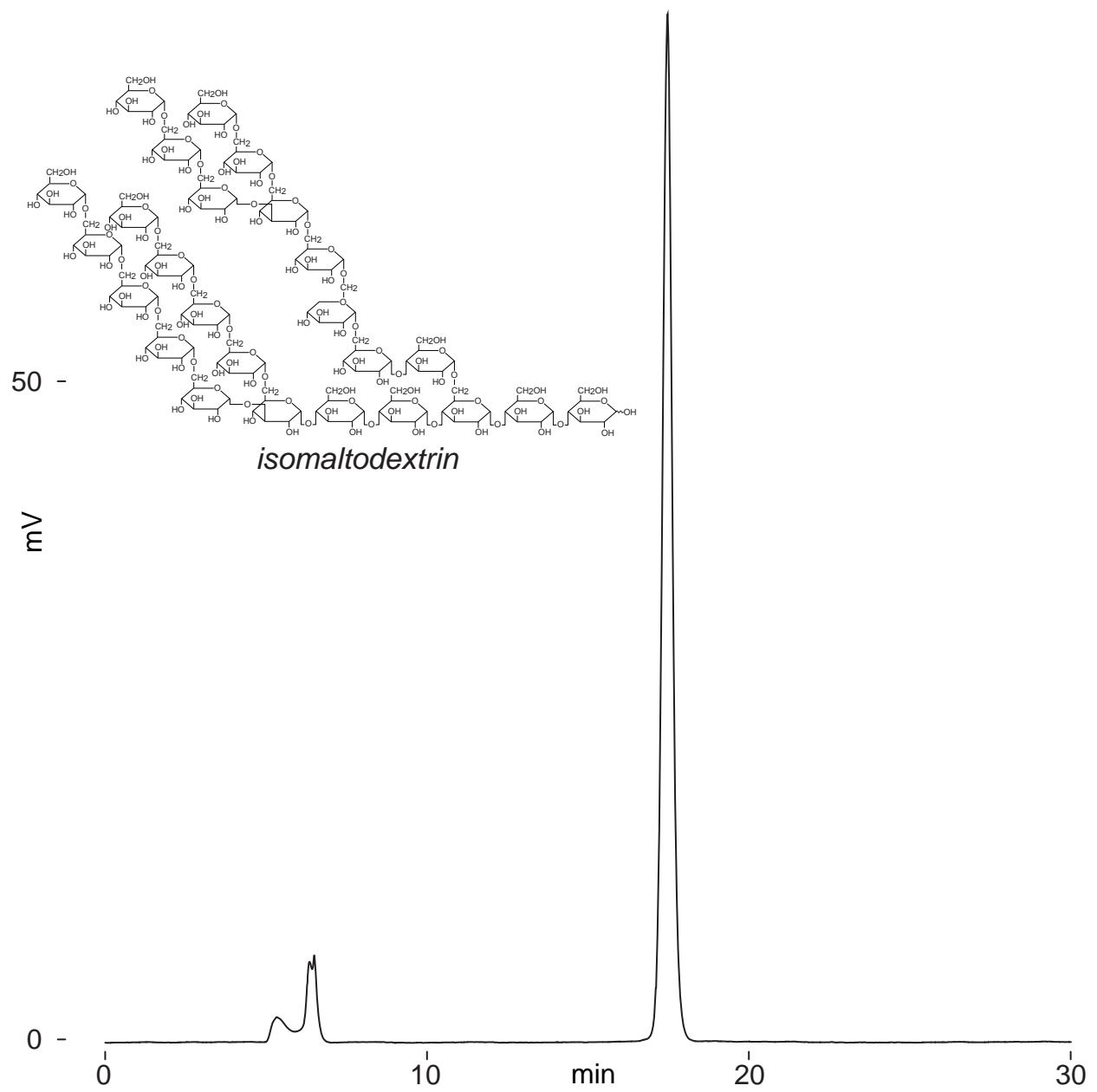


Presto FF-C18

250 x 4.6 mm

HPLC COLUMN APPLICATION

Isomaltodextrin
イソマルトデキストリン

Presto FF-C18, 250 x 4.6 mm
A: water / formic acid = 100 / 0.1
B: water / acetonitrile / formic acid = 10 / 90 / 0.1
0-50% B (0-30min)
0.3mL/min (20MPa), 37°C,
ELSD (spray chamber 30°C, drift chamber 60°C)
5μL (5mg/mL, water / formic acid = 100 / 1)