

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18

Cadenza CW-C18


Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Styrene urinary metabolites

尿中スチレン代謝物

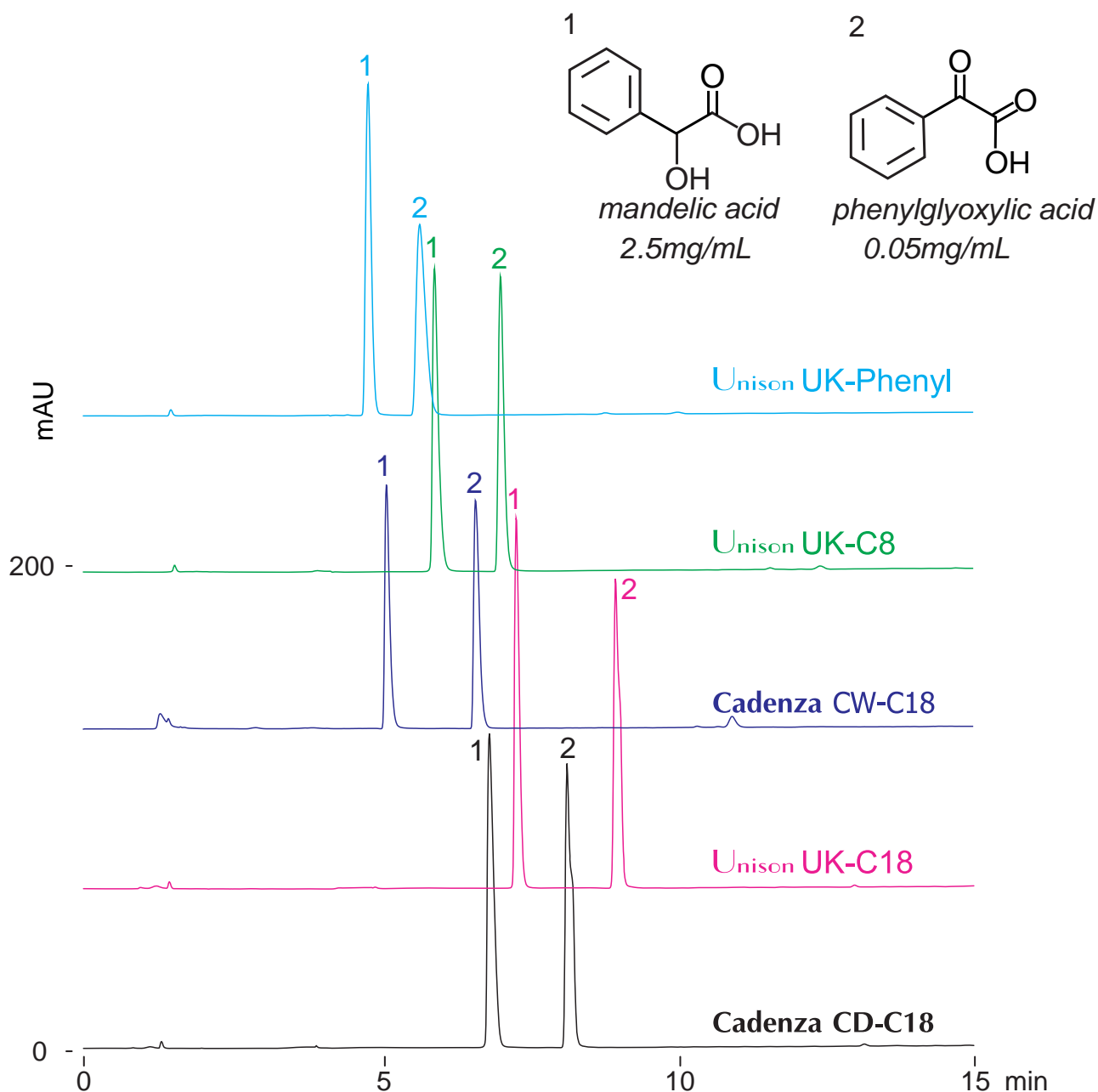

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-40%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.05-2.5mg/mL, water / formic acid = 100 / 1)