

Unison UK-Phenyl Cadenza CW-C18



Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Rutin (flavonoid glycoside)

ルチン (フラボノイド配糖体)

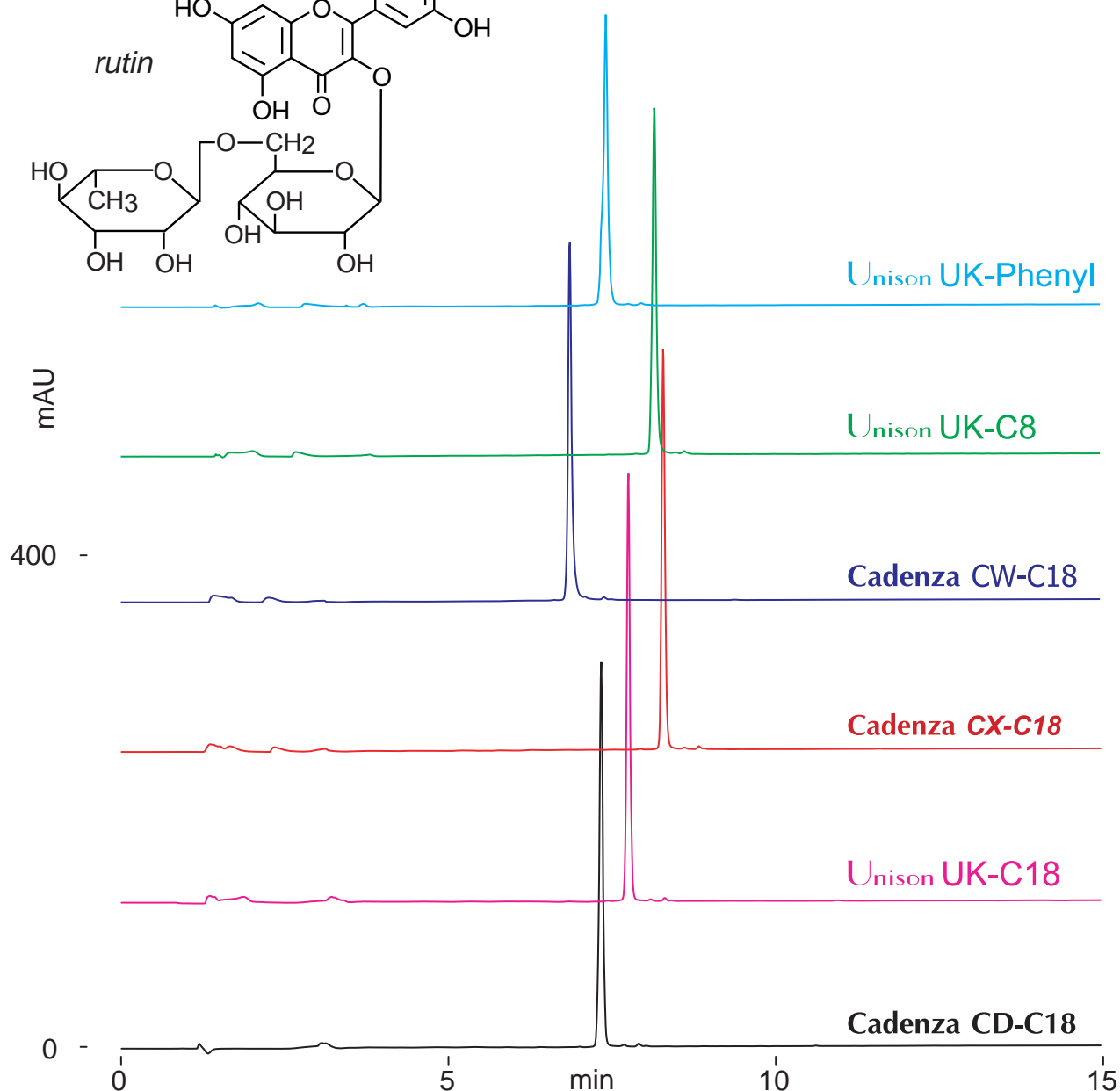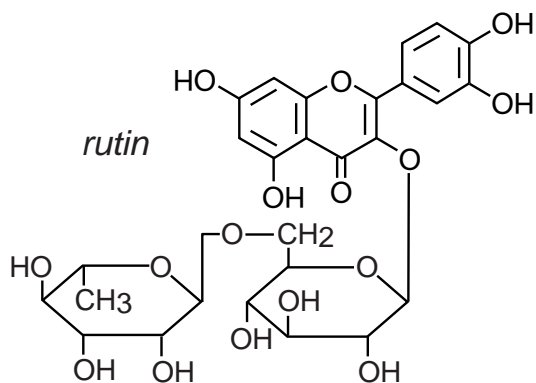

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-65%B (0-15min)

1.0mL/min (9MPa), 37°C, 360nm

10µL (0.15mg/mL, methanol / formic acid = 100 / 1)