


Unison UK-C8
 Unison UK-C18
 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Ceftriaxone (Cephem antibiotics)

セフトリアキソン (セフェム系抗生物質)

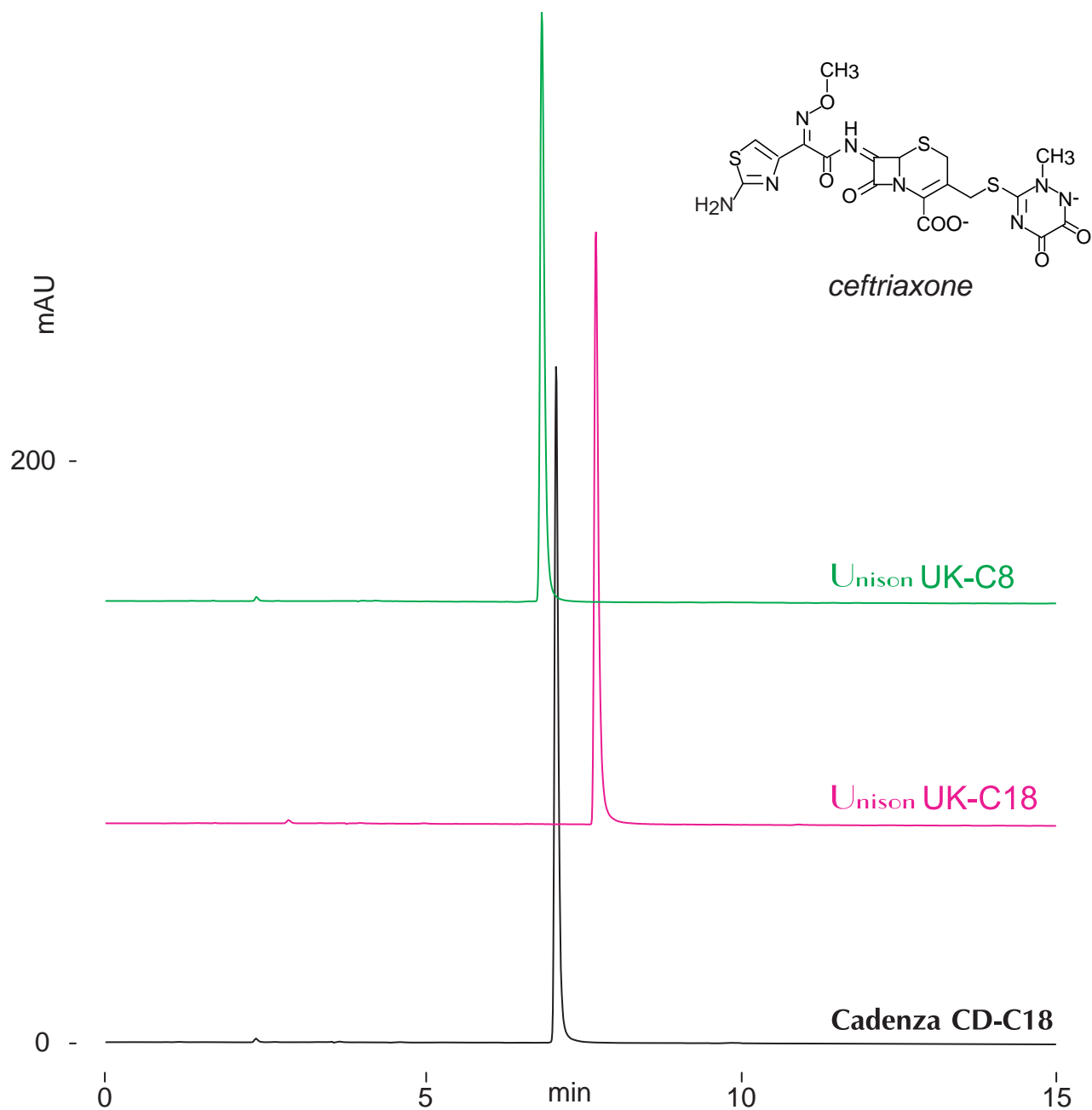

100 x 4.6 mm
 A: 100mM ammonium formate
 B: 100mM ammonium formate / acetonitrile = 50 / 50
 0-40%B (0-15min)
 1.0mL/min (9MPa), 37°C, 260nm
 10μL (0.075mg/mL, 100mM ammonium formate)