


Unison UK-C18 MF

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Cephem antibiotics

セフェム系抗生物質の LC-MS 分析

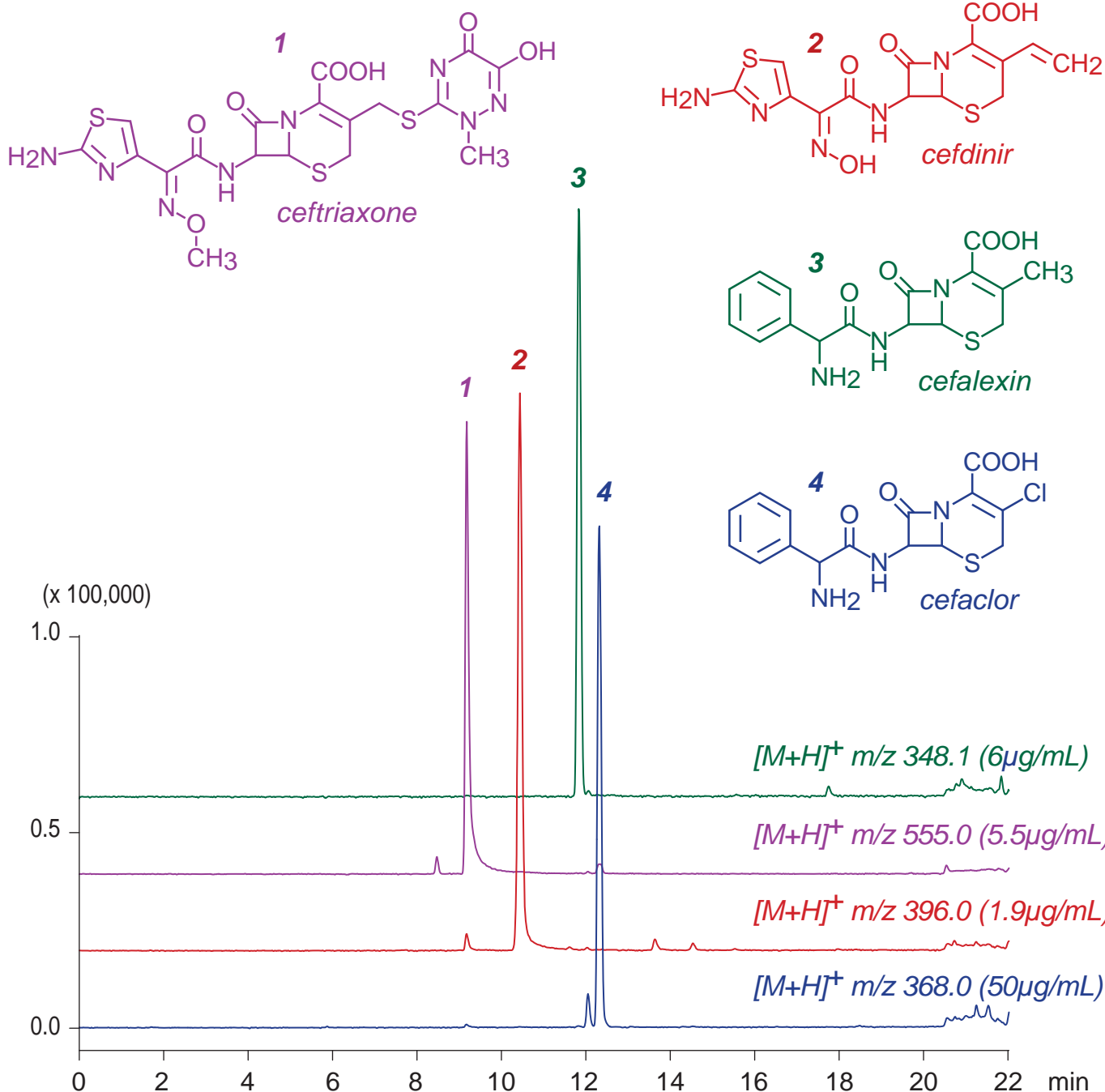

Unison UK-C18 MF, 150 x 2 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-30%B (0-18min), 100%B (18-22min)

0.2mL/min (11MPa), 37°C, Single Quad. MS (ESI)

10 μL (1.9-50 $\mu\text{g/mL}$, 100mM ammonium formate)