

Unison UK-Phenyl Cadenza CW-C18



Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Bezafibrate (Hypolipidemic agent)

ベザフィブラート (高脂血症治療薬)

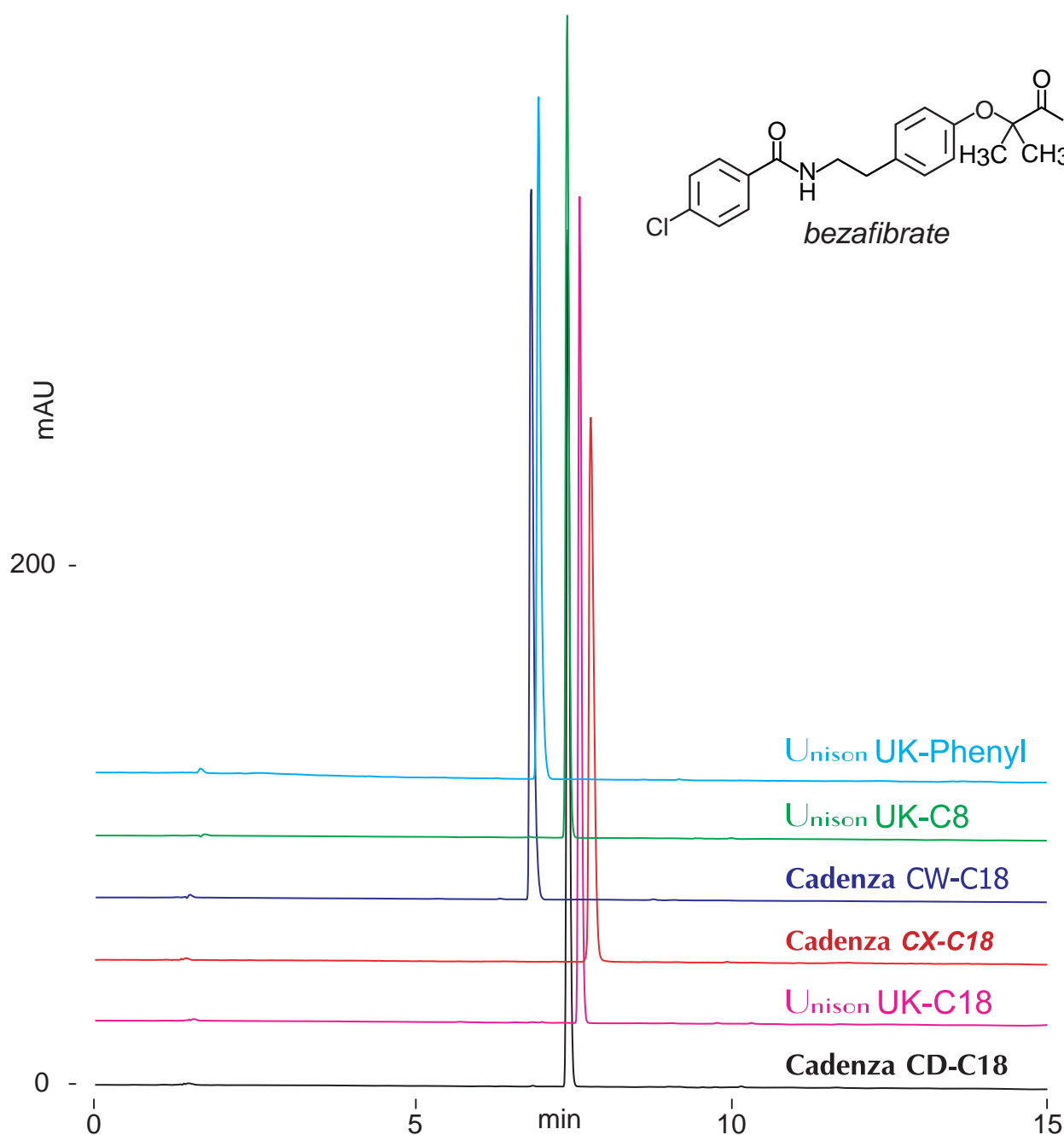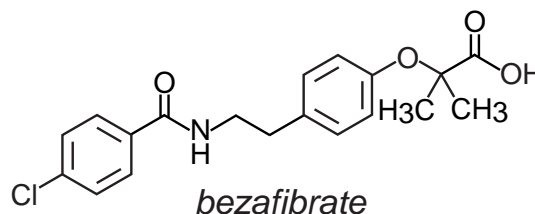

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

10-80%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.25mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)