


Unison UK-Phenyl

250 x 4.6 mm


HPLC COLUMN APPLICATION

Chlorophenylboronic acid


クロロフェニルボロン酸

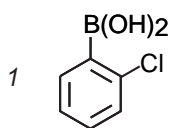


2-chlorophenyl boronic acid

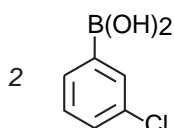


3-chlorophenyl boronic acid

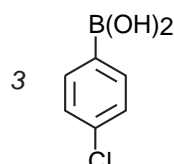

4-chlorophenyl boronic acid

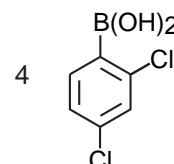

2,4-Dichlorophenyl boronic Acid

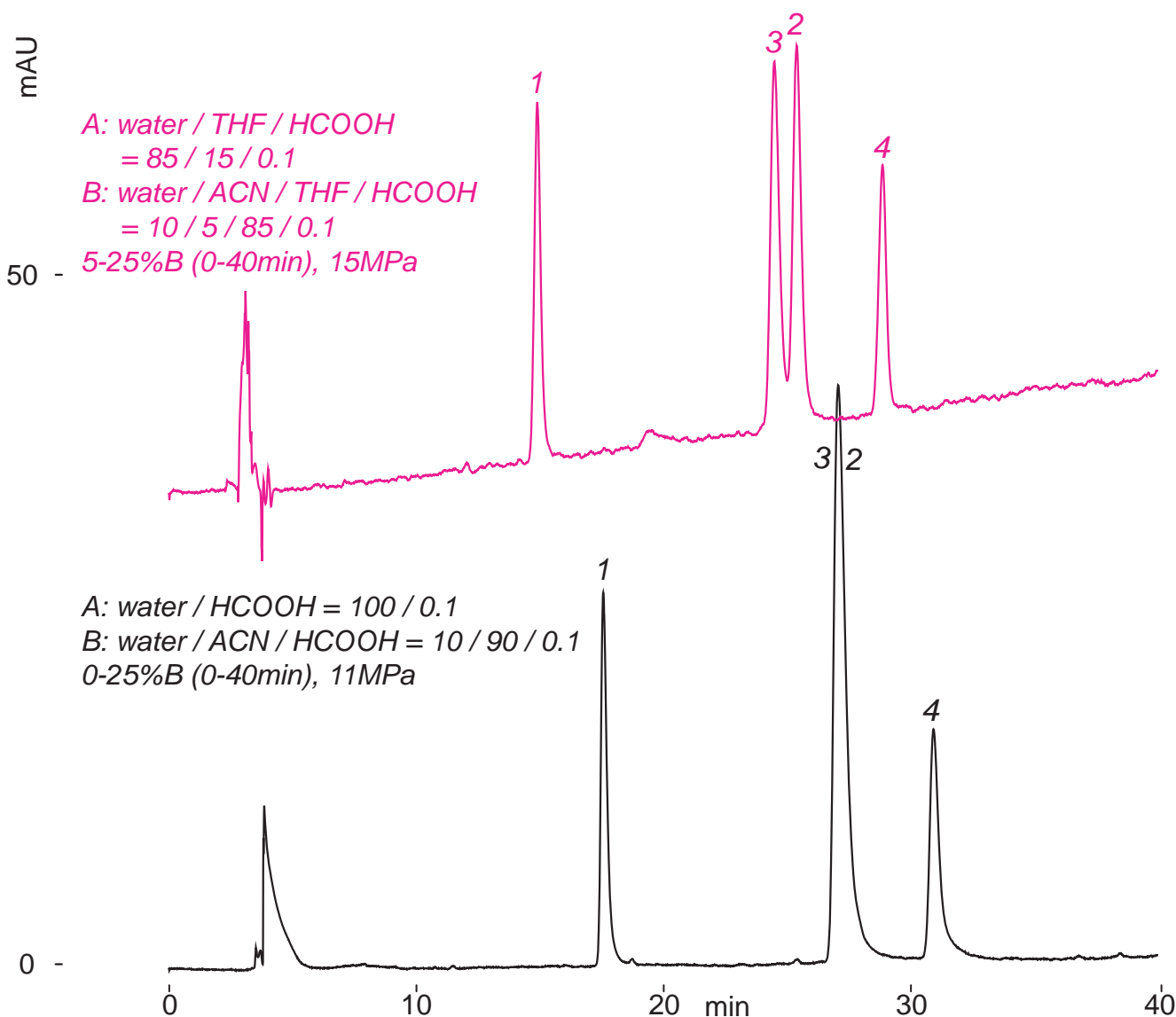

Unison UK-Phenyl, 250 x 4.6 mm
 1mL/min, 60°C, 260nm
 10µL (0.75mg/mL, DMSO/ formic acid = 100 / 1)