

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Tamsulosin (Antidysuria)

タムスロシン (排尿障害改善薬)

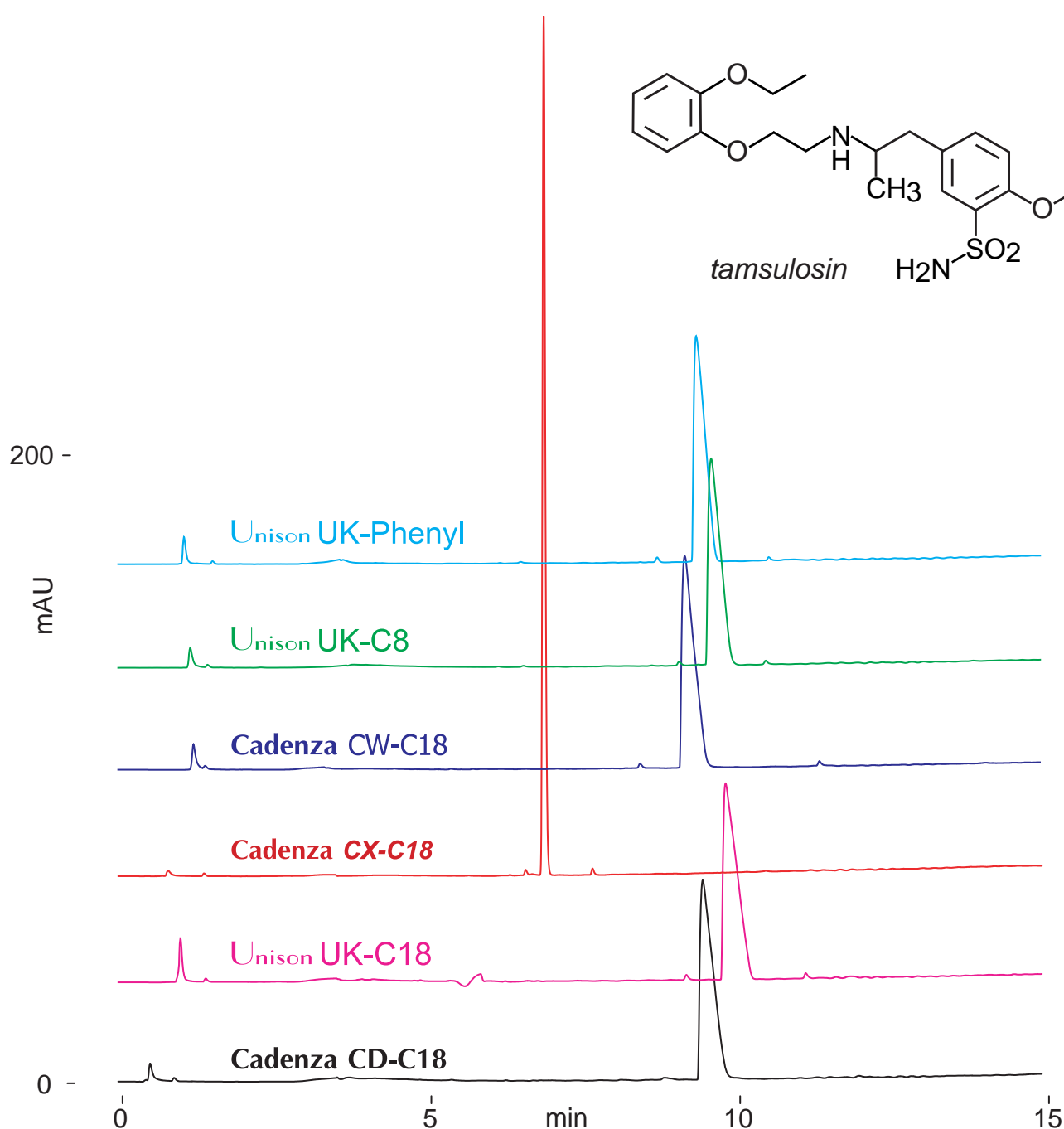

100 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

0-60%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, water / formic acid = 100 / 1)