

Nardis ND-RX


Unison UK-Phenyl

150 x 4.6 mm

HPLC COLUMN APPLICATION

### Catecholamines

カテコールアミン類の分析

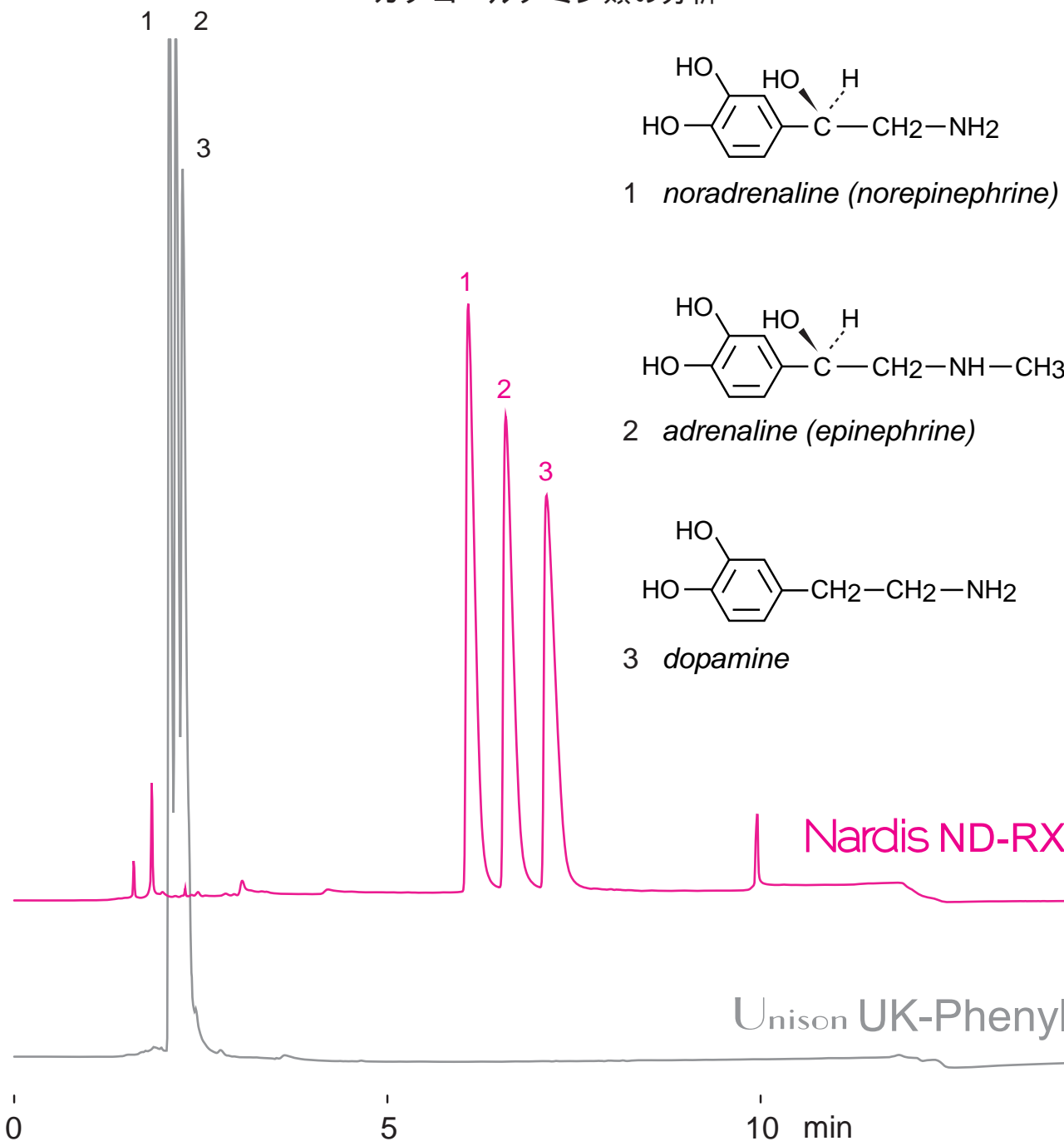

150 x 4.6 mm

A: 10mM ammonium acetate / acetonitrile = 90 / 10

B: 50mM ammonium acetate-acetic acid (5:5, pH 4.7) / acetonitrile = 70 / 30

0 -100%B (0 -10 min), 1.0 mL/min, 37°C, 260nm

10 µL (1.0mg/mL, water/HCOOH=100/1)