

Nardis ND-NX


Unison UK-Amino

150 x 4.6 mm

HPLC COLUMN APPLICATION

Acyclovir, Ganciclovir (Antiviral)

アシクロビル, ガンシクロビル (抗ウイルス薬)

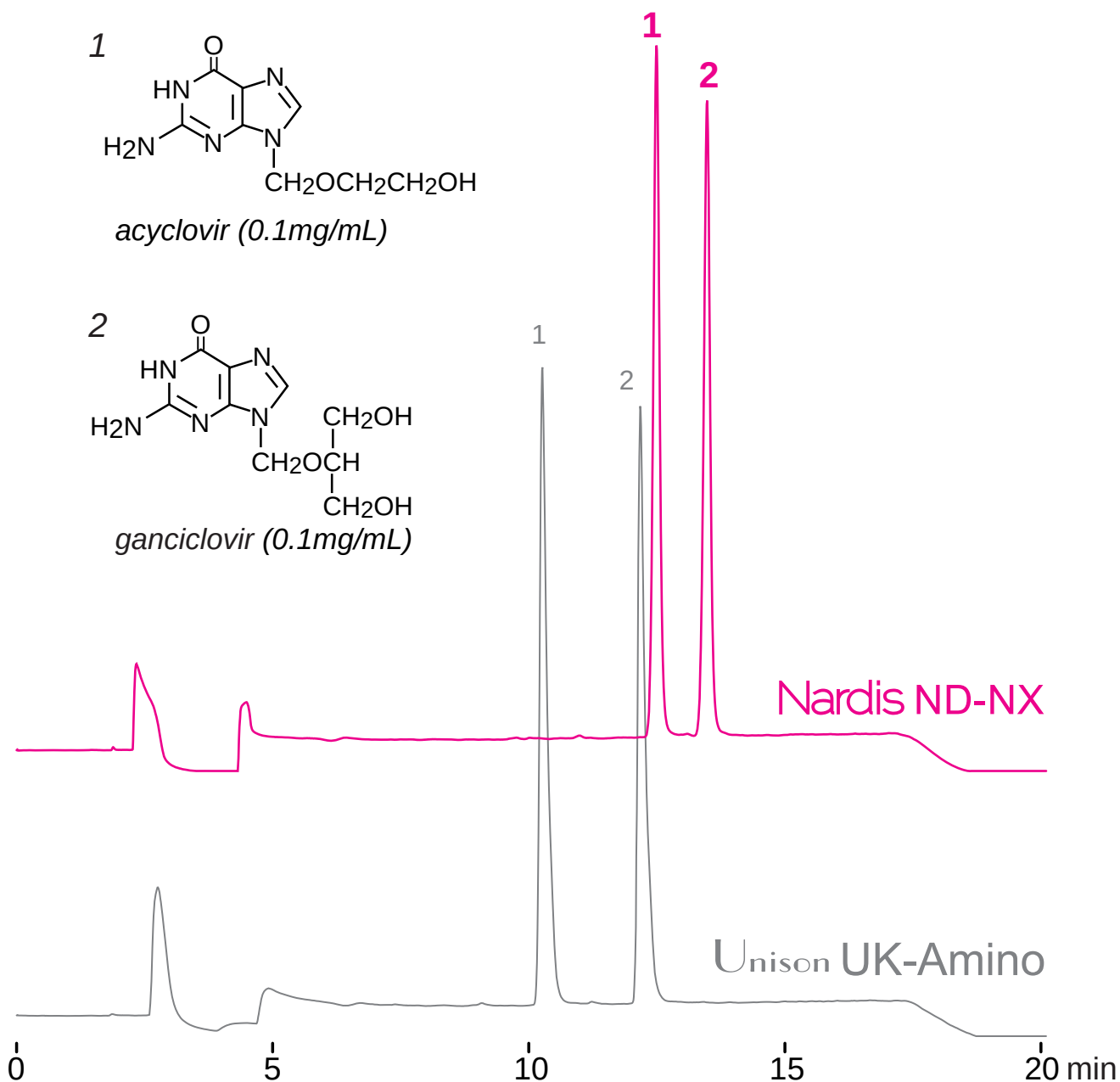

150 x 4.6 mm

A: acetonitrile / formic acid = 100 / 0.1

B: acetonitrile / water / formic acid = 50 / 50 / 0.5

0-50%B (0-15min)

1.0mL/min (4-6 MPa) , 37°C, 260nm

10μL (0.1mg/mL, water/formic acid=100/1)