

Nardis ND-RX


Unison UK-Phenyl

150 x 4.6 mm

HPLC COLUMN APPLICATION

Indole derivatives

インドール誘導体

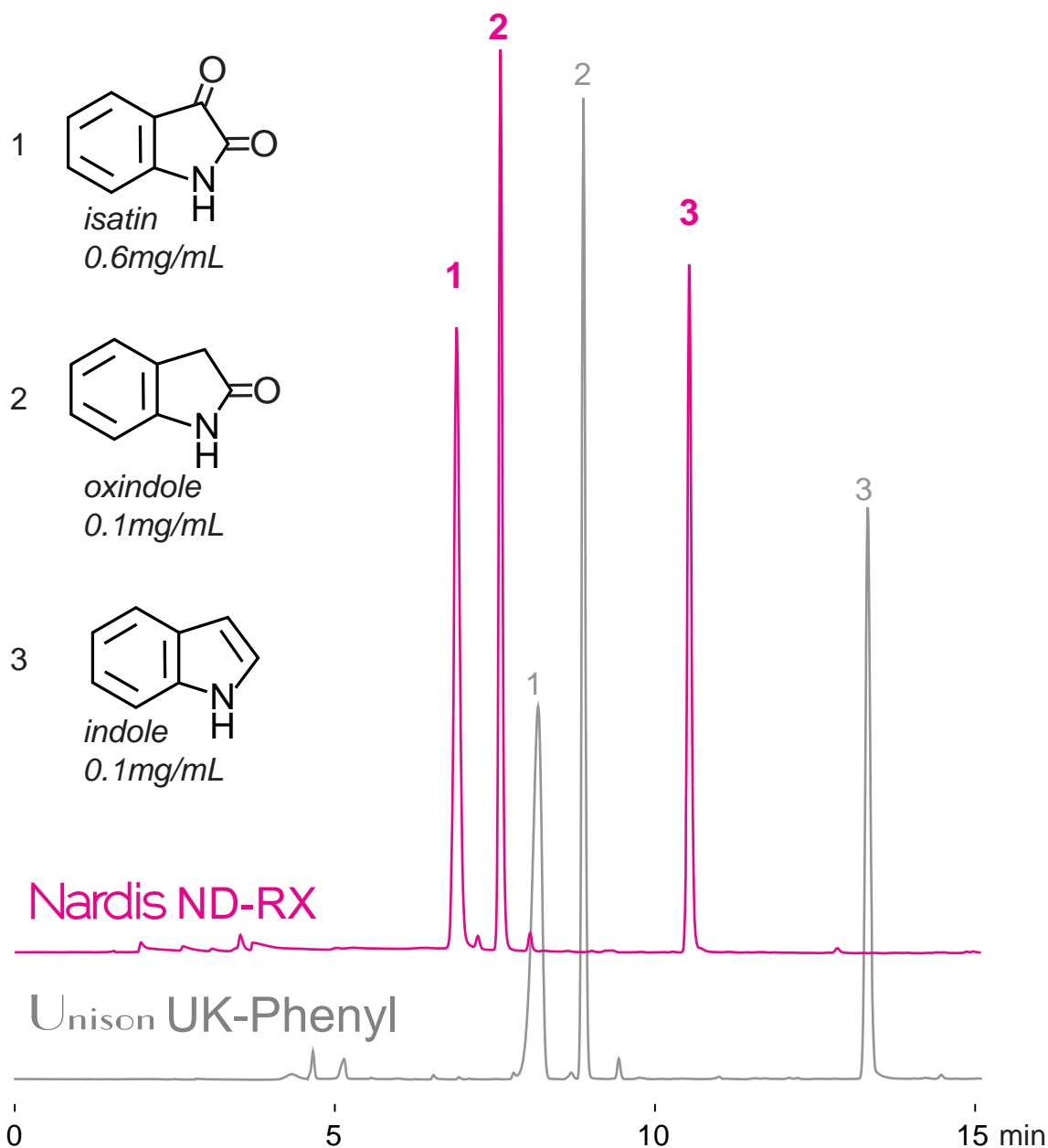

150 x 4.6 mm

A: 10mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-60%B (0-15min)

1.0mL/min (10-12MPa), 37°C, 260nm

10µL (0.1-0.6mg/mL, 100mM HCOONH₄-HCOOH(1:9)/ ACN=70/30)