


Nardis ND-NX **MF**

150 x 2 mm


**HPLC COLUMN APPLICATION**

Coenzymes (Nicotinamide adenine dinucleotide phosphate)


補酵素 (ニコチンアミドアデニンジヌクレオチドリン酸)

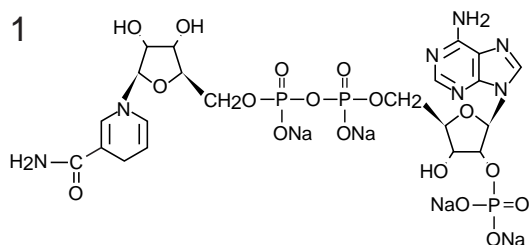

NADPH (0.6mg/mL)

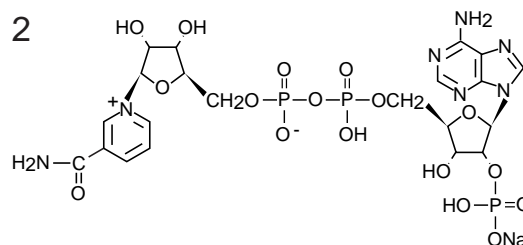

NADP (0.2mg/mL)

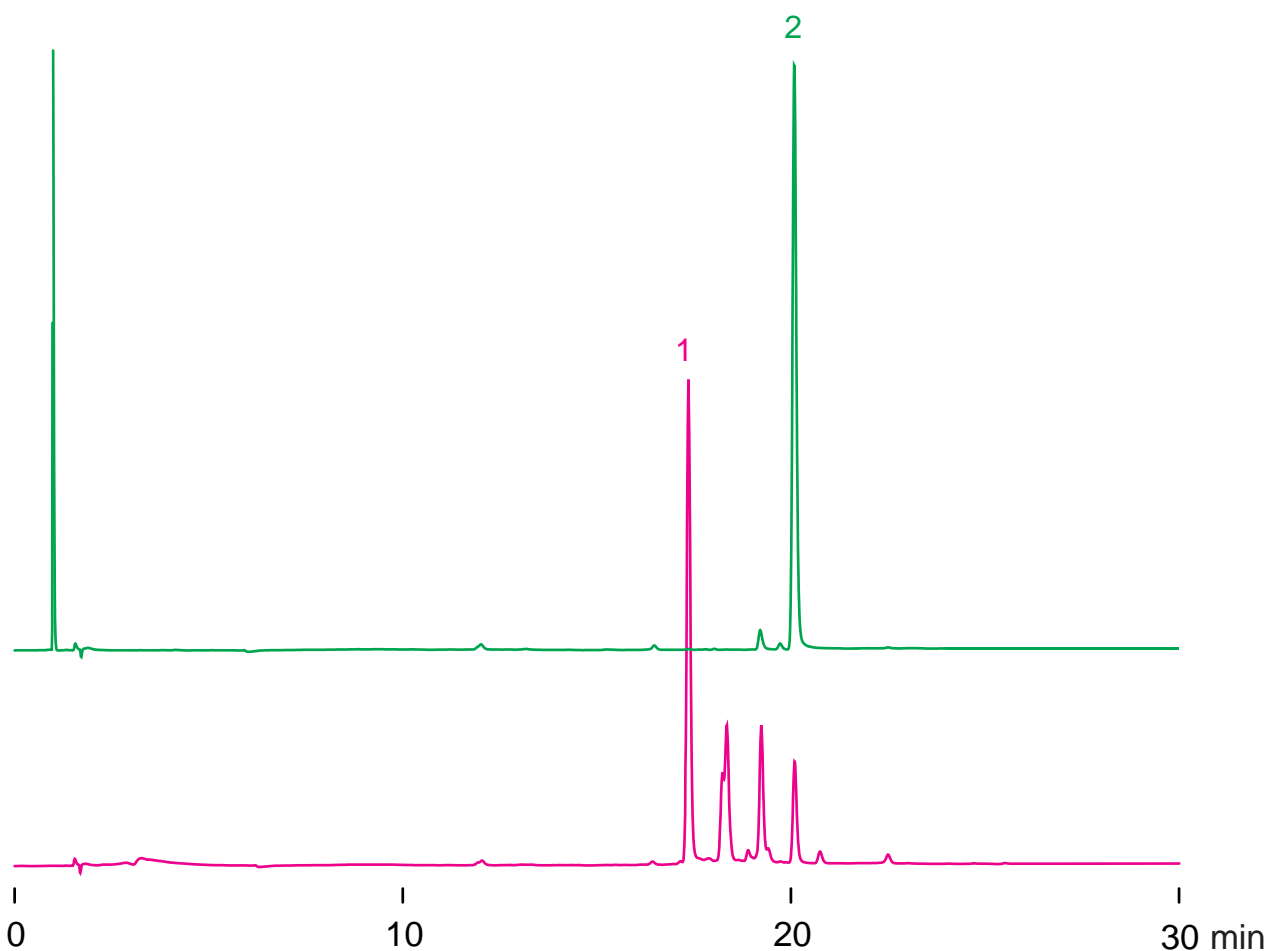

Nardis ND-NX MF, 150 x 2 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium acetate = 10 / 90

20-70%B (0-30min)

0.2mL/min (5MPa), 37°C, 260nm

2µL (0.2-0.6mg/mL, 100mM ammonium acetate)