

Nardis ND-NX


Unison UK-Amino

250 x 4.6 mm


HPLC COLUMN APPLICATION

Xylulose, xylitol, xylose


キシロース, キシリトール, キシロース

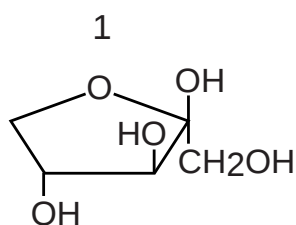


xylulose
(2mg/mL)

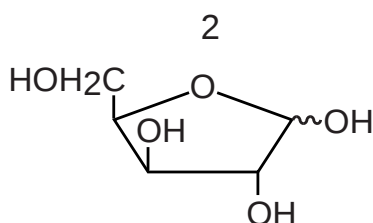

xylose
(1mg/mL)

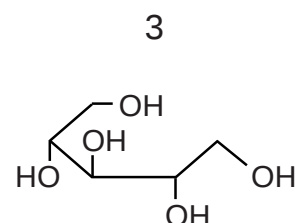

xylitol
(1mg/mL)

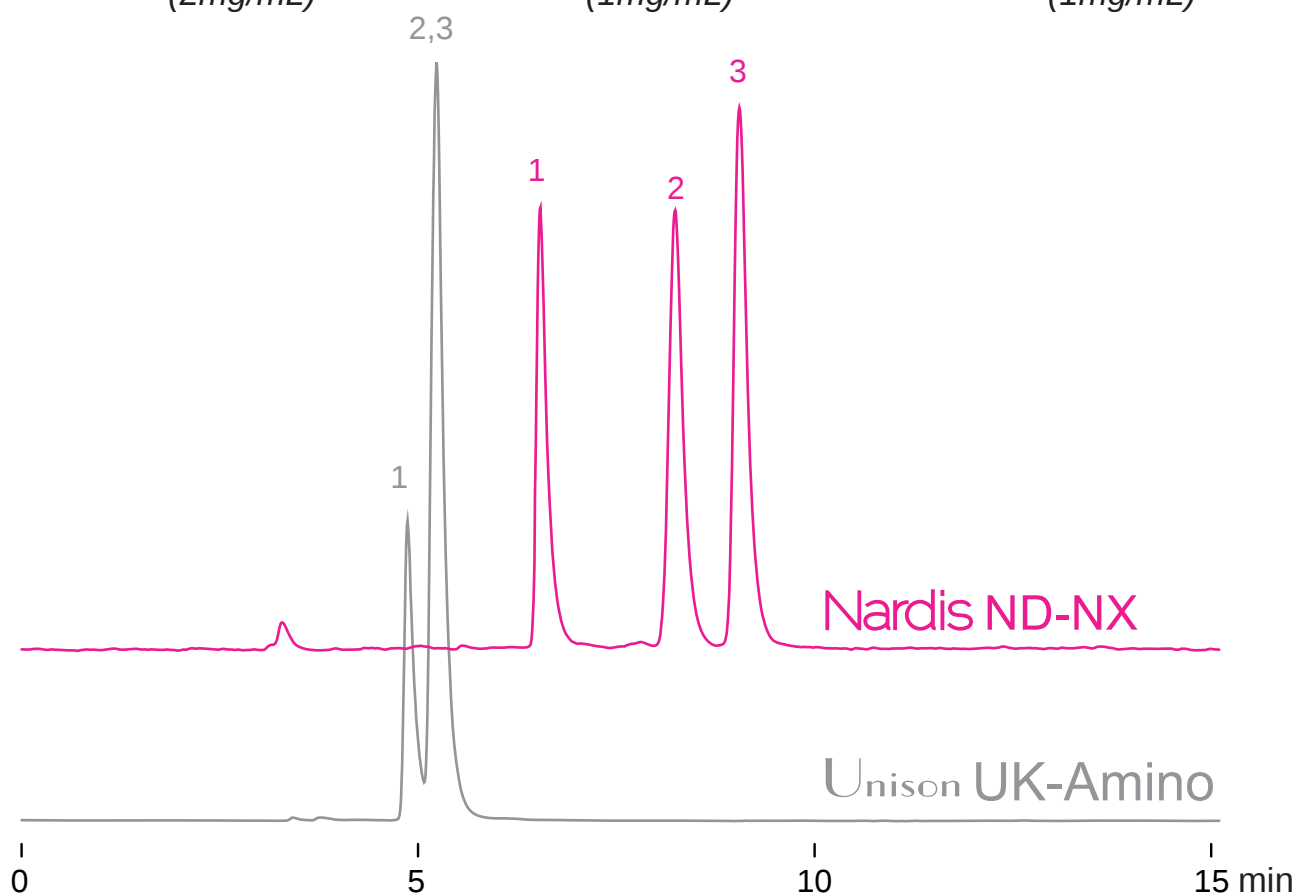

250 x 4.6 mm

A: acetonitrile

B: acetonitrile / 10mM ammonium acetate = 50 / 50

0 -100%B (0 -10 min), 100%B (10 -15 min)

0.8 mL/min (8 MPa) , 60°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10 μ L (1-2mg/mL, water)