


Unison UK-Silica
Nardis ND-NX

150 x 4.6 mm

HPLC COLUMN APPLICATION

Gemcitabine (Antineoplastic)

ゲムシタビン (抗悪性腫瘍薬)

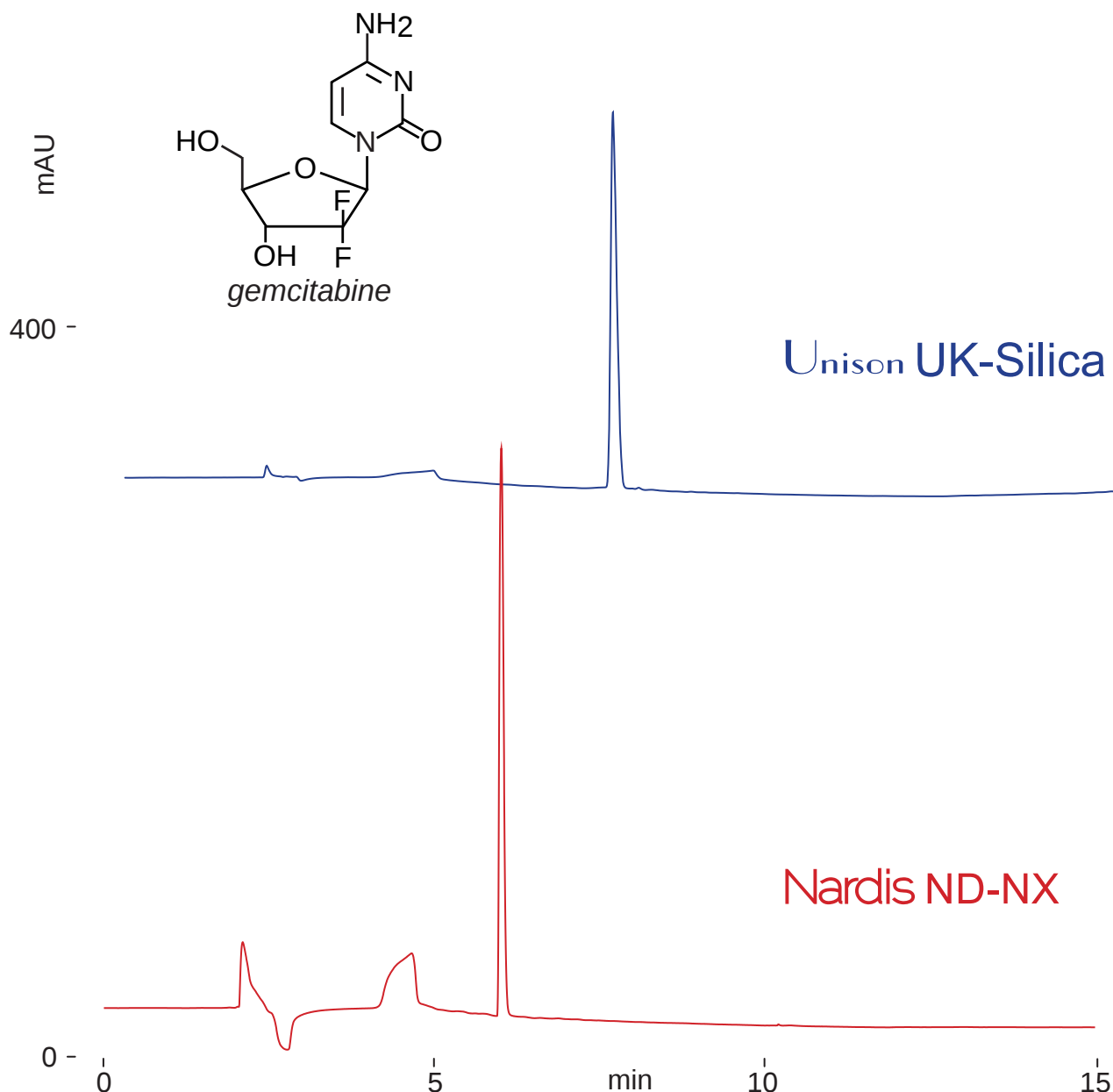

150 x 4.6 mm

A : acetonitrile / formic acid = 100 / 0.5

B : water / acetonitrile / formic acid = 50 / 50 / 0.5

0%B (0-1min), 0-100%B (1-10min), 100%B (10-15min)

1mL/min (4-6MPa), 37°C, 260nm

10μL (0.1mg/ml, water / formic acid = 100 / 1)