

Nardis ND-RP


Unison UK-Phenyl

150 x 4.6 mm


HPLC COLUMN APPLICATION

Neuroblastoma related compounds


神経芽細胞腫関連化合物

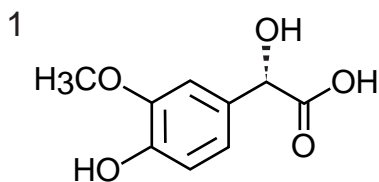


L-vanillylmandelic acid
(VMA, 0.5mg/mL)

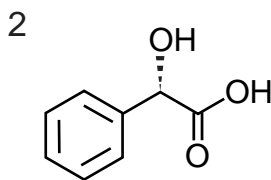

L-mandelic acid
(MA, 2.5mg/mL)

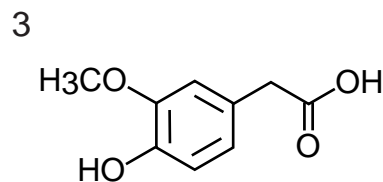

homovanillic acid
(HVA, 0.5mg/mL)

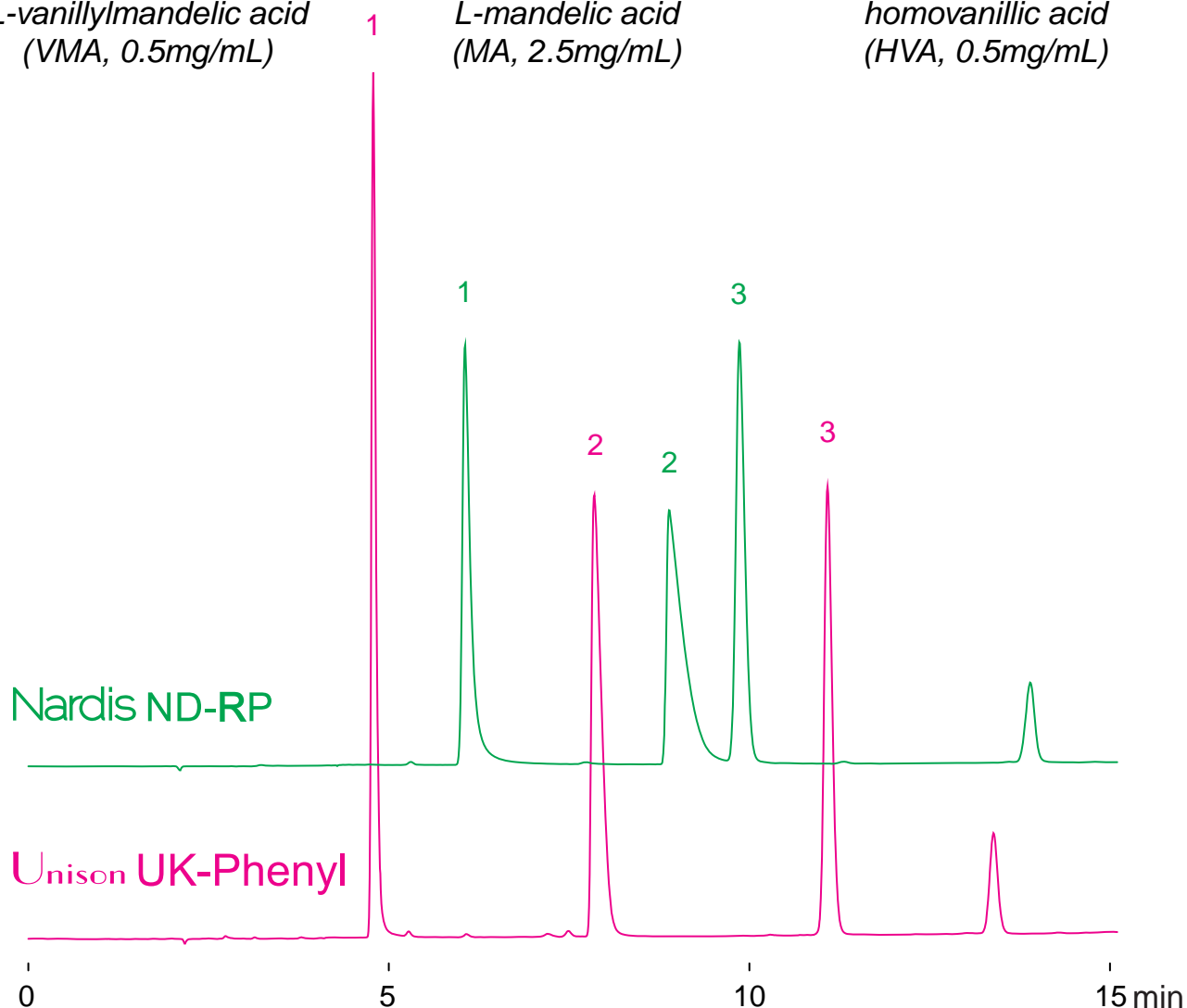

150 x 4.6 mm

A: water / formic acid = 100 / 0.5

B: water / acetonitrile / formic acid = 50 / 50 / 0.5

0-25%B (0-15min)

1.0mL/min (11-12MPa), 37°C, 260nm

10µL (0.5-2.5mg/mL, water)