


Unison UK-C18

150 x 3 mm


HPLC COLUMN APPLICATION

LC-MS: Aminonaphtholsulfonic acid isomers


アミノナフトールスルホン酸異性体の LC-MS 分析

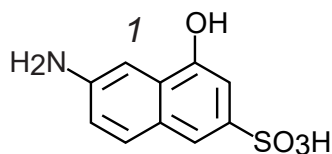


1
7-amino-1-naphthol
-3-sulfonic acid

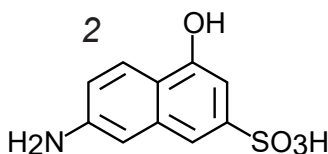



2
6-amino-1-naphthol
-3-sulfonic acid

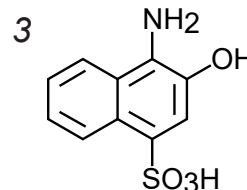

3
1-amino-2-naphthol
-4-sulfonic acid

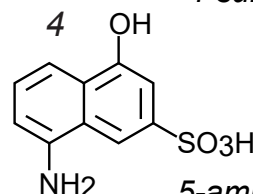

4
5-amino-1-naphthol
-3-sulfonic acid

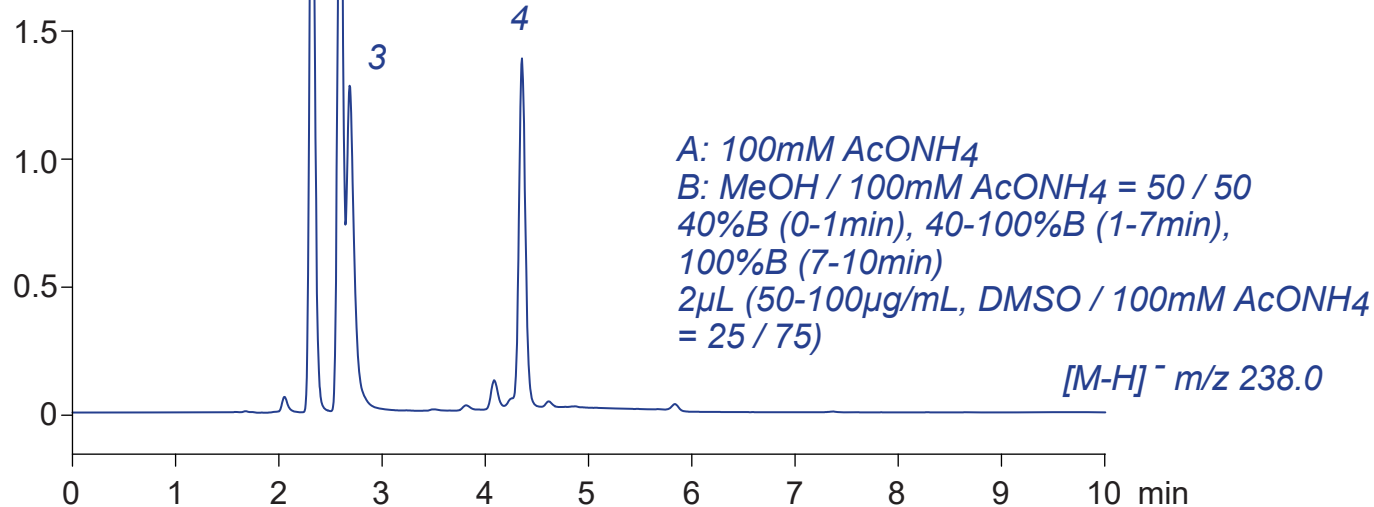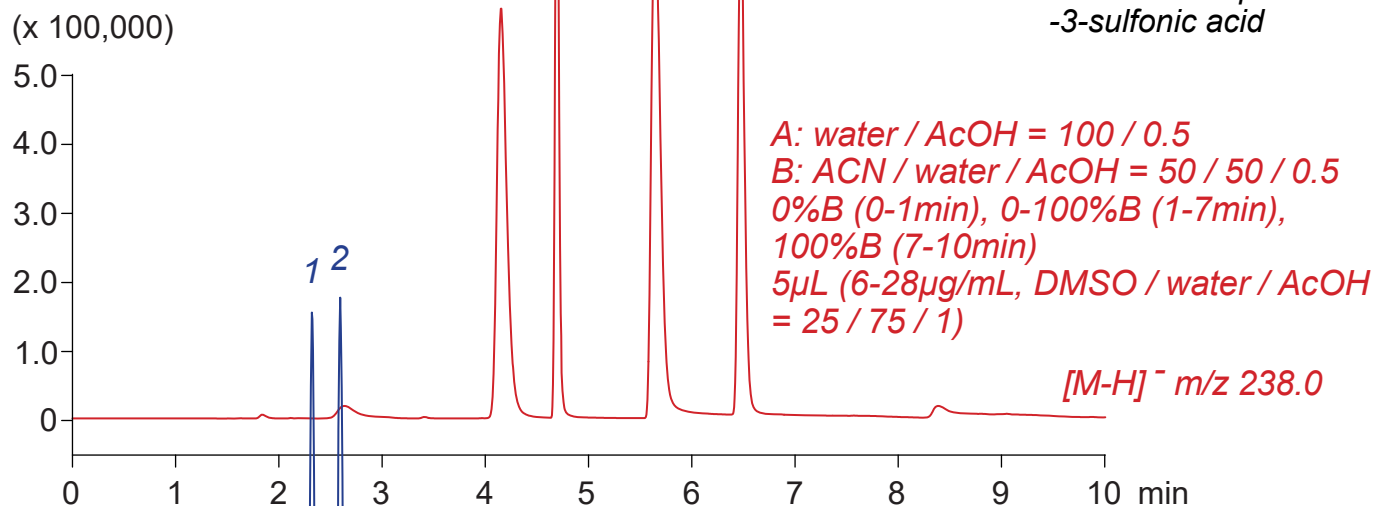

Unison UK-C18, 150 x 3 mm
0.4mL/min (8-12MPa), 37°C, Single Quad. MS (ESI)