

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Rotenone (Botanical insecticide)

ロテノン (植物性殺虫成分)

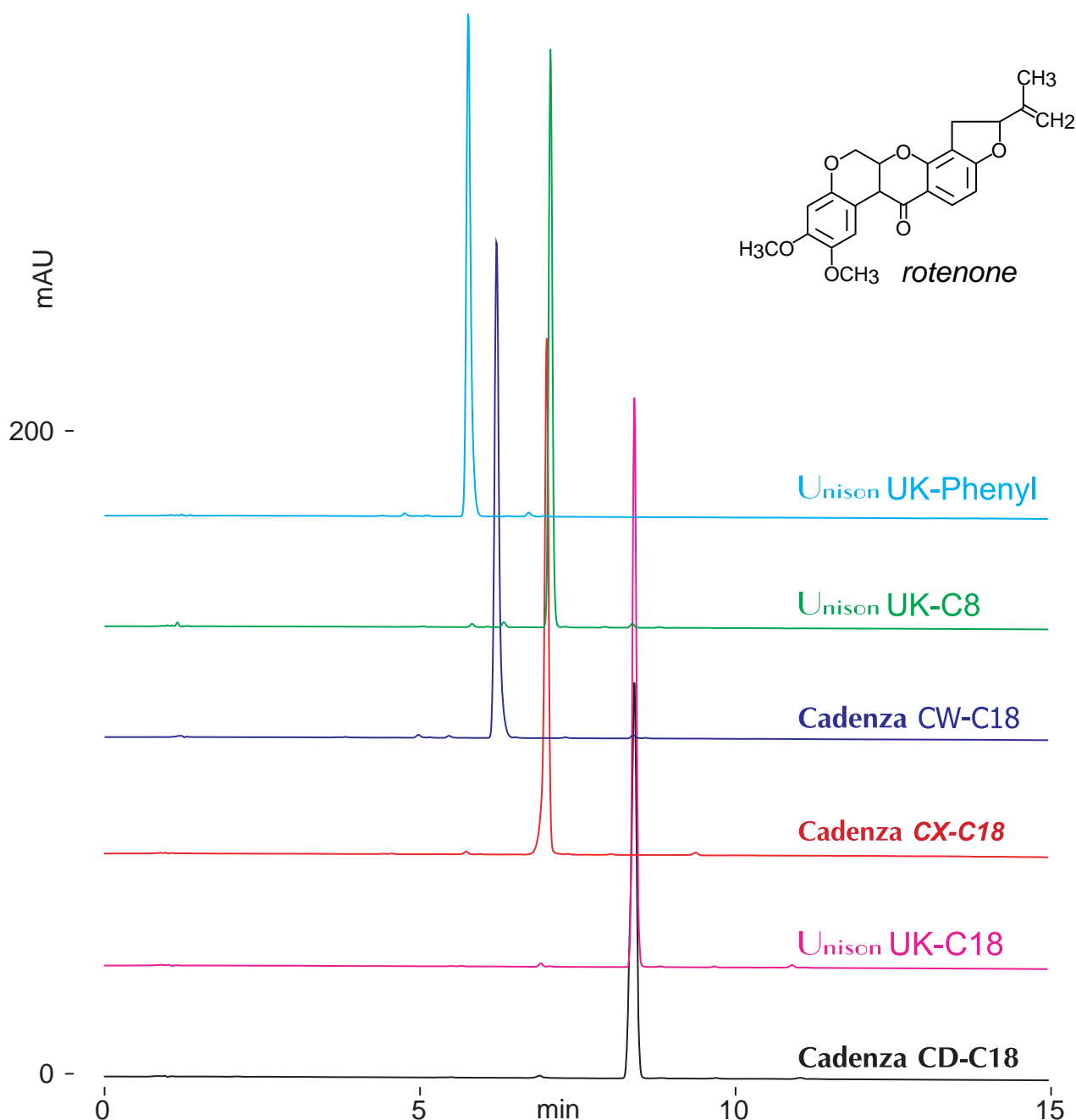

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

45-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 290nm

10µL (0.075mg/mL, water / acetonitrile / formic acid = 50 / 50 / 1)