

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Sesamin (sesame lignan)

セサミン (ゴマリグナン)

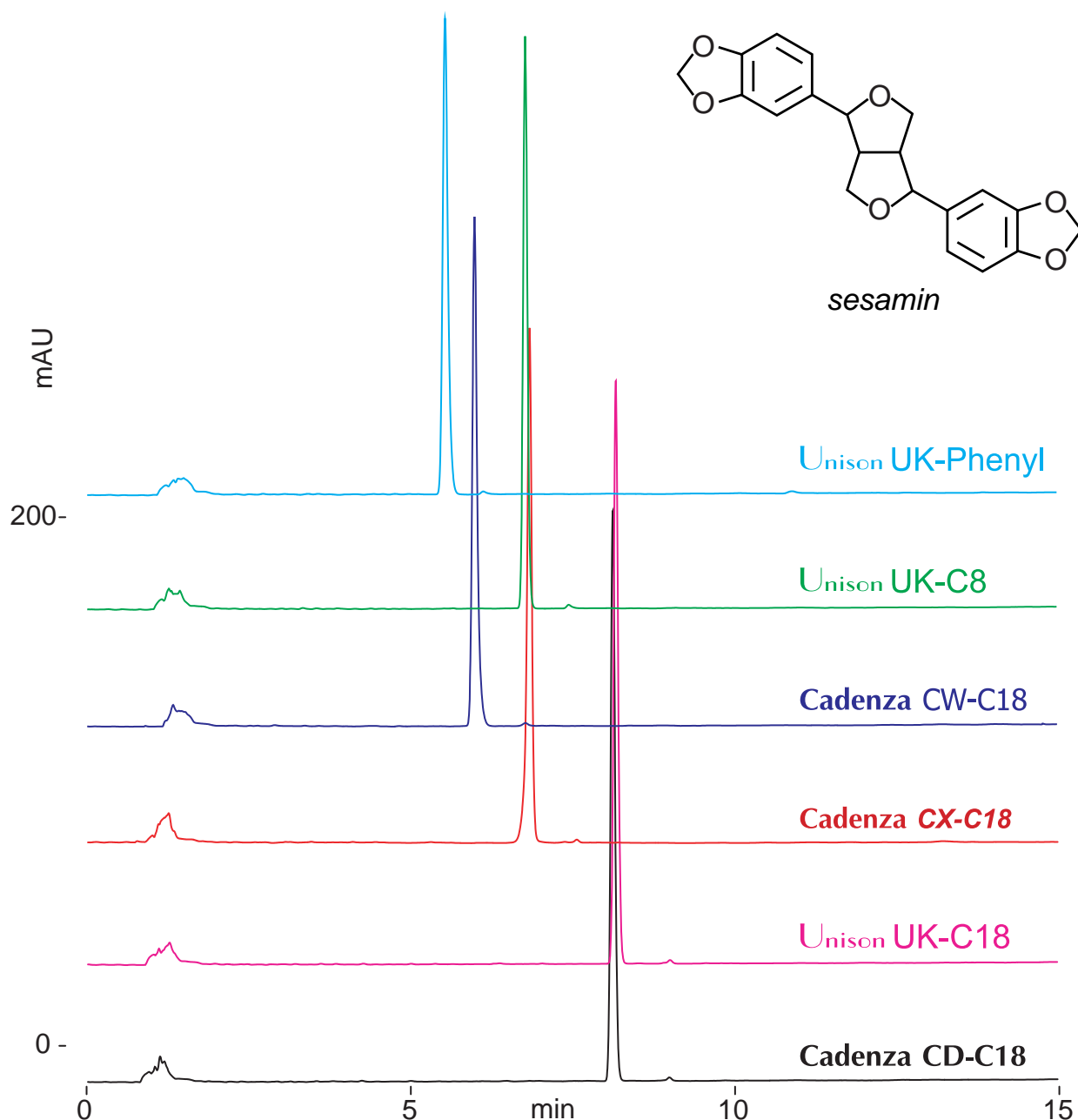

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

40-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 290nm

10µL (0.2mg/mL, dimethyl sulfoxide / formic acid = 100 / 1)