

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Phloretin (Polyphenol)

フロレチン (ポリフェノール)

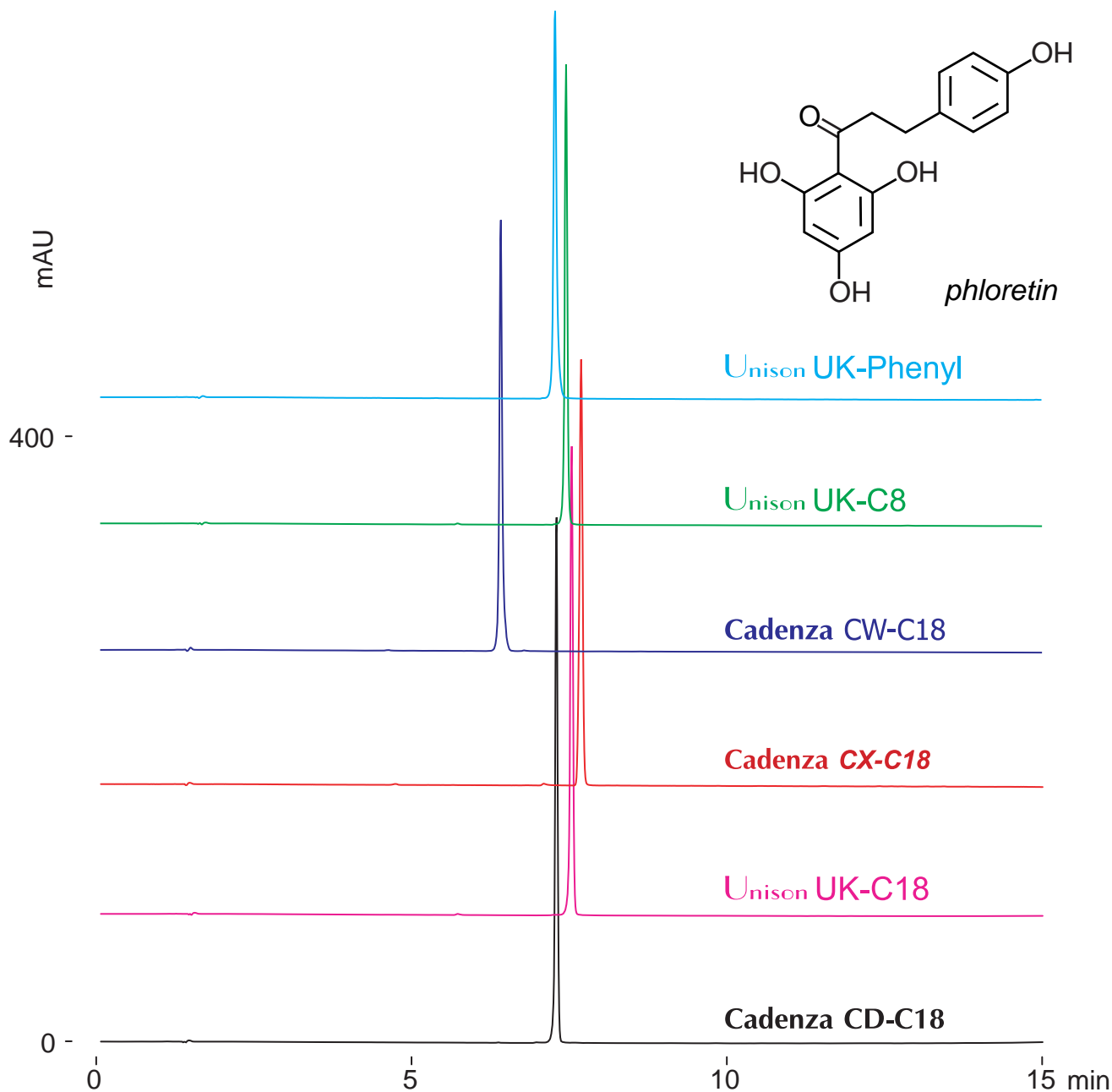

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

15-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.2mg/mL, 100mM ammonium formate / acetonitrile = 50 / 50)