

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Indocyanine Green (Near infrared fluorescent dye)

インドシアニングリーン (近赤外蛍光色素)

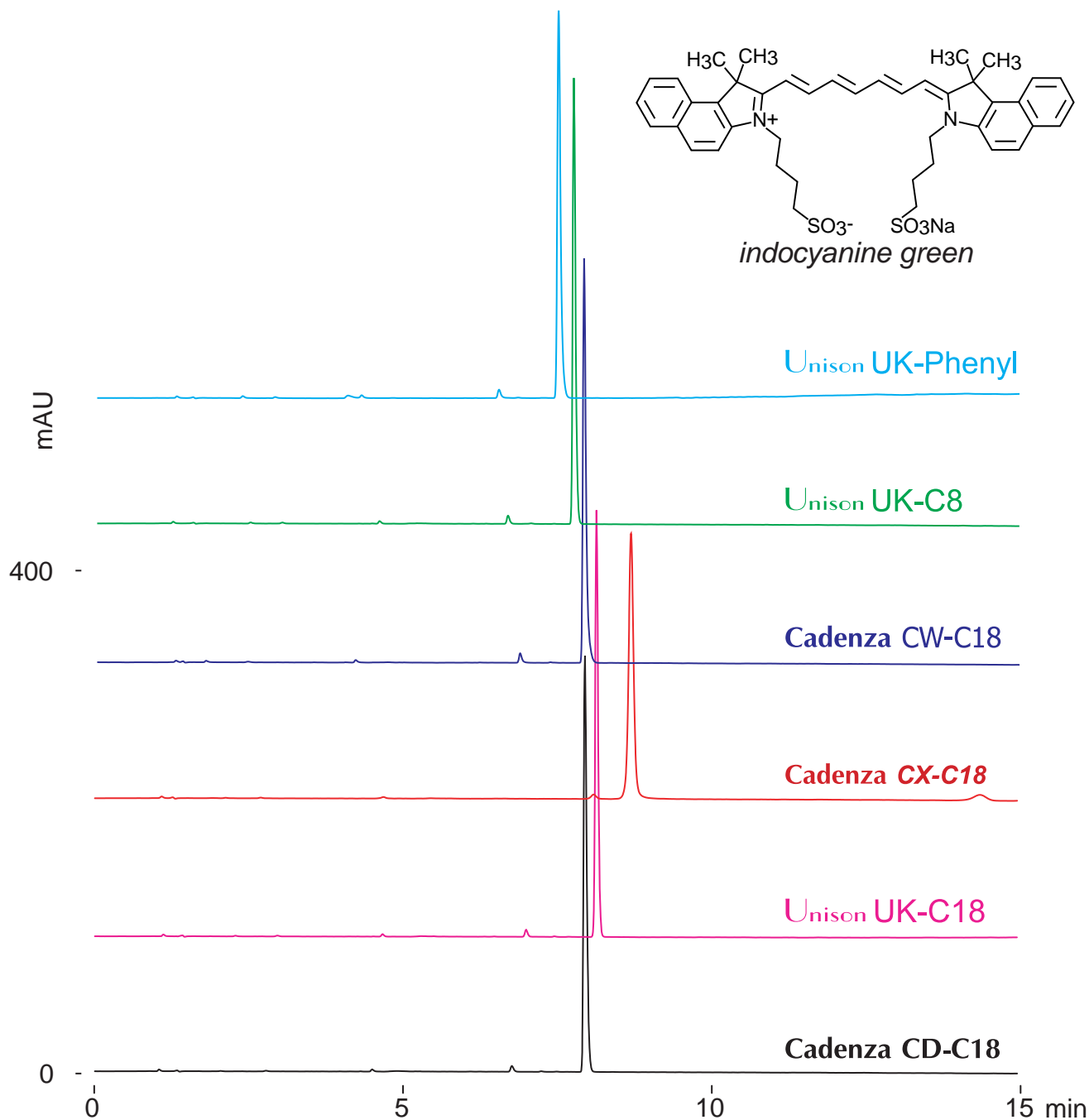

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

25-100%B (0-15min)

1.0mL/min (11MPa), 37°C, 260nm

10µL (0.15mg/mL, 100mM ammonium formate)