

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Gramine (Indole alkaloid)

グラミン (インドールアルカロイド)

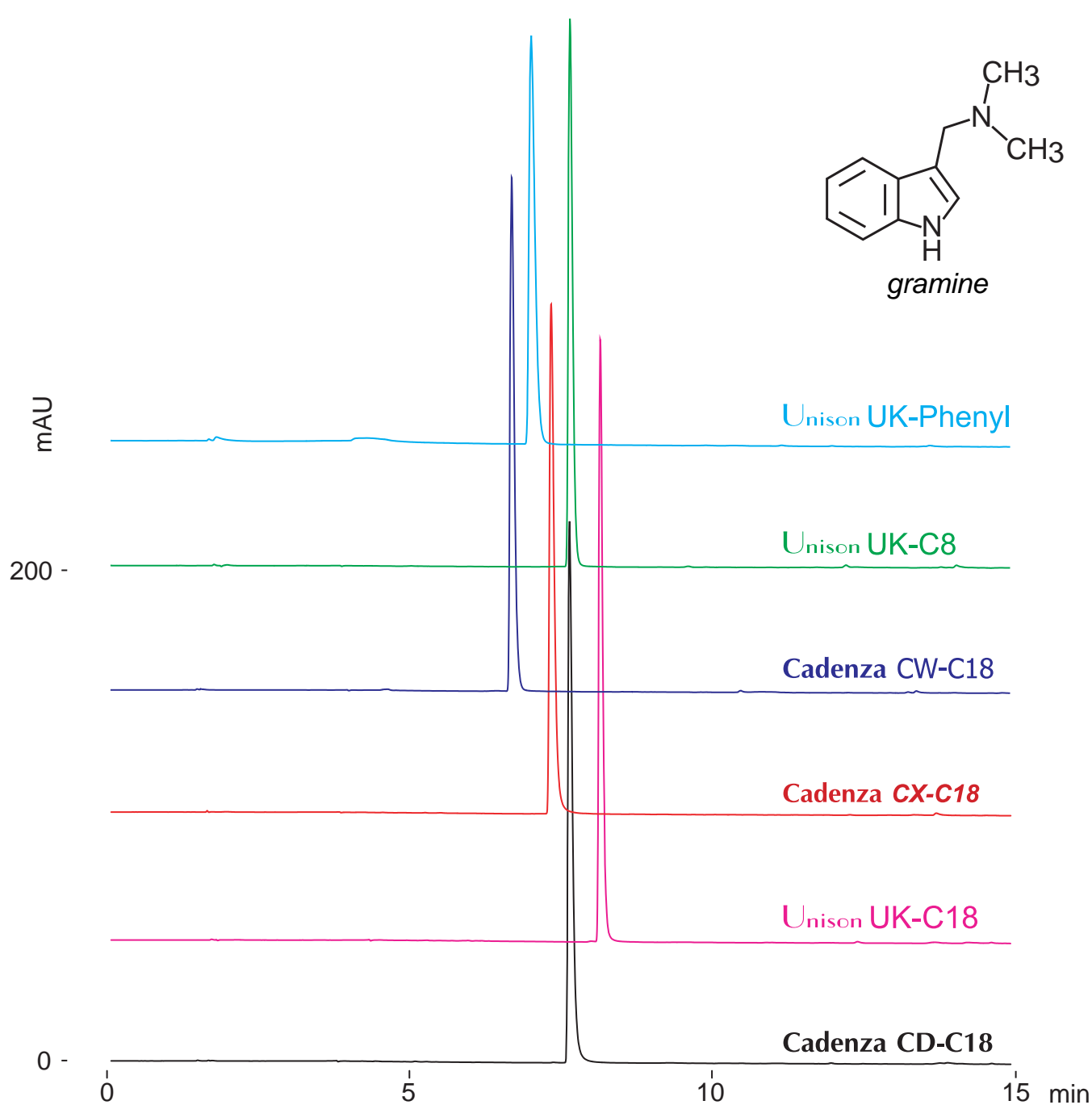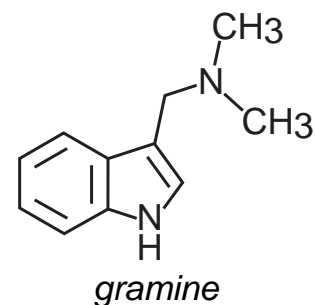

100 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 10 / 90

0%B (0-1min), 0-50%B (1-15min)

1.0mL/min (11MPa), 37°C, 260nm

10µL (0.12mg/mL, 100mM ammonium formate / acetonitrile = 67 / 33)