

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Flavokawain A (Chalcones)

フラボカワインA(カルコン類)

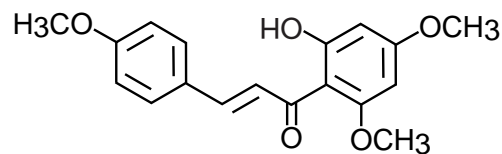

flavokawain A

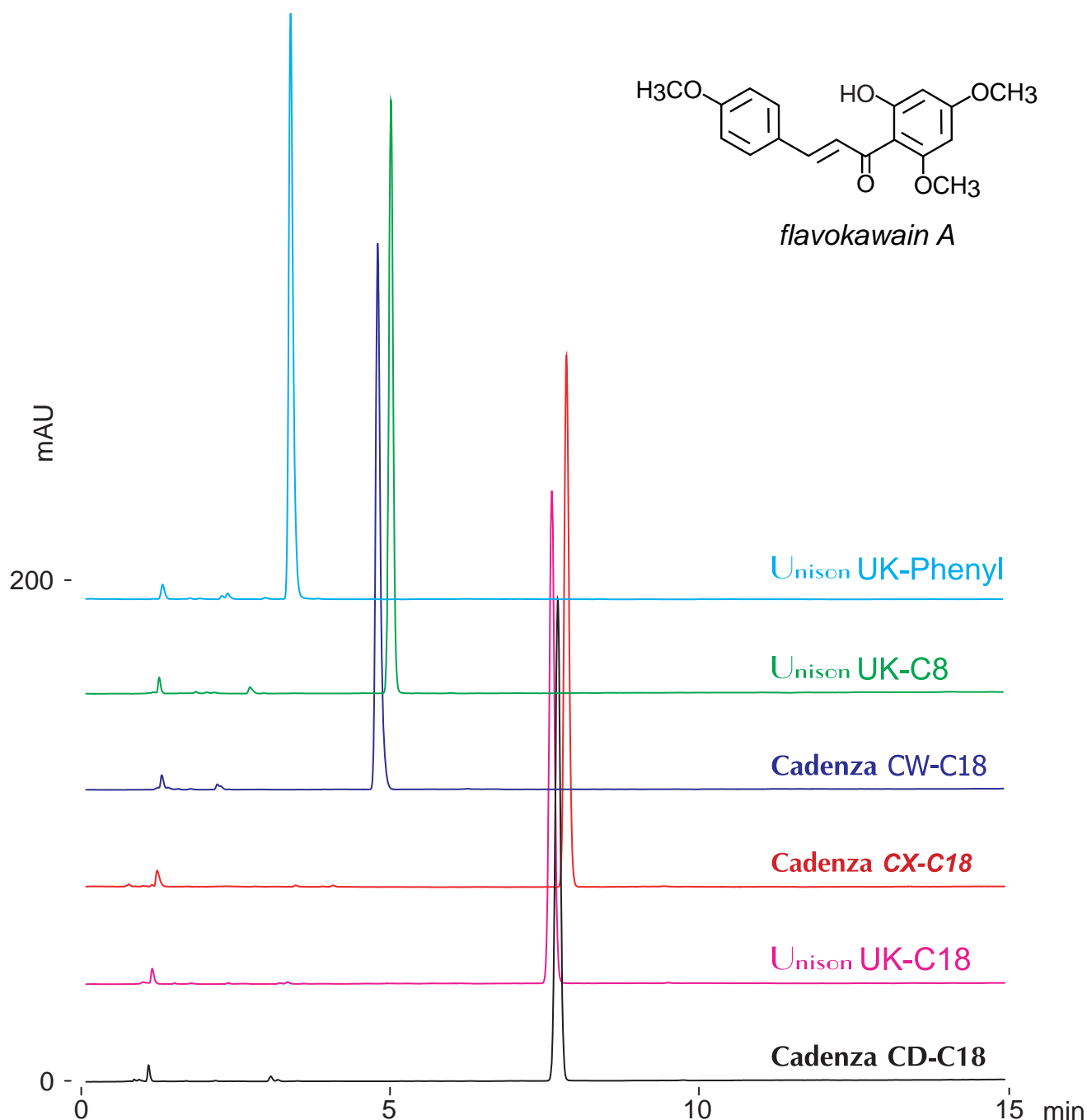

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

10-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.3mg/mL, water / acetonitrile / formic acid = 40 / 60 / 1)