

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18


Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Dicoumarol (Coumarin anticoagulant)

ジクマロール (クマリン系抗凝固薬)

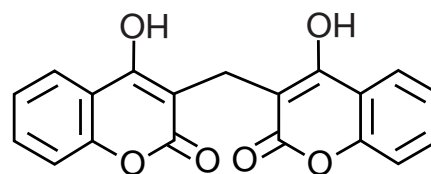

dicoumarol

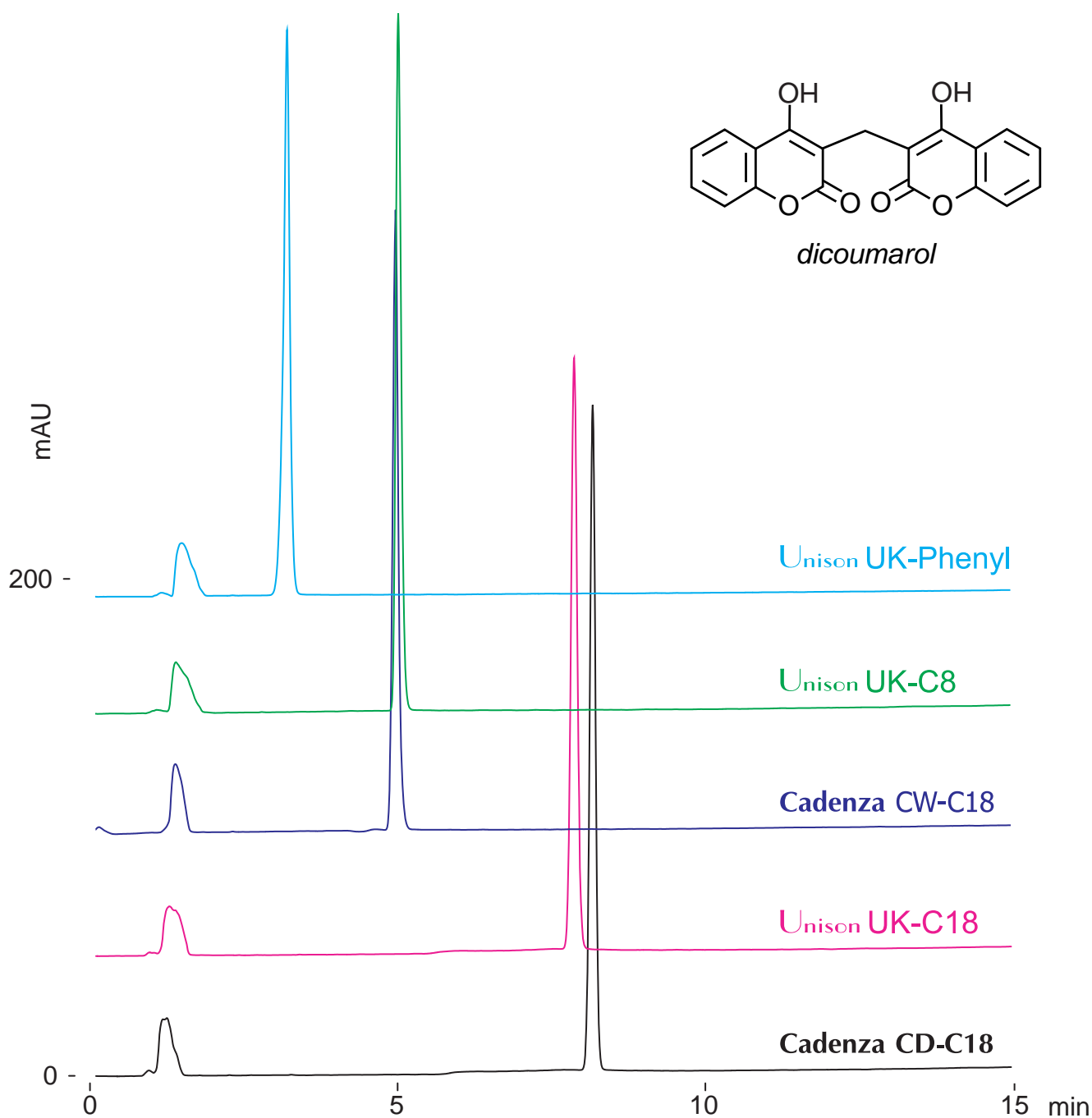

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.5

B: water / acetonitrile / formic acid = 10 / 90 / 0.5

10-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.2mg/mL, water / DMF / formic acid = 20 / 80 / 1)