


Unison UK-Phenyl Cadenza CX-C18
 Unison UK-C8 Cadenza CW-C18
 Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Methylcobalamin (Coenzyme form of Vitamin B12)
 メチルコバラミン (補酵素型ビタミン12)

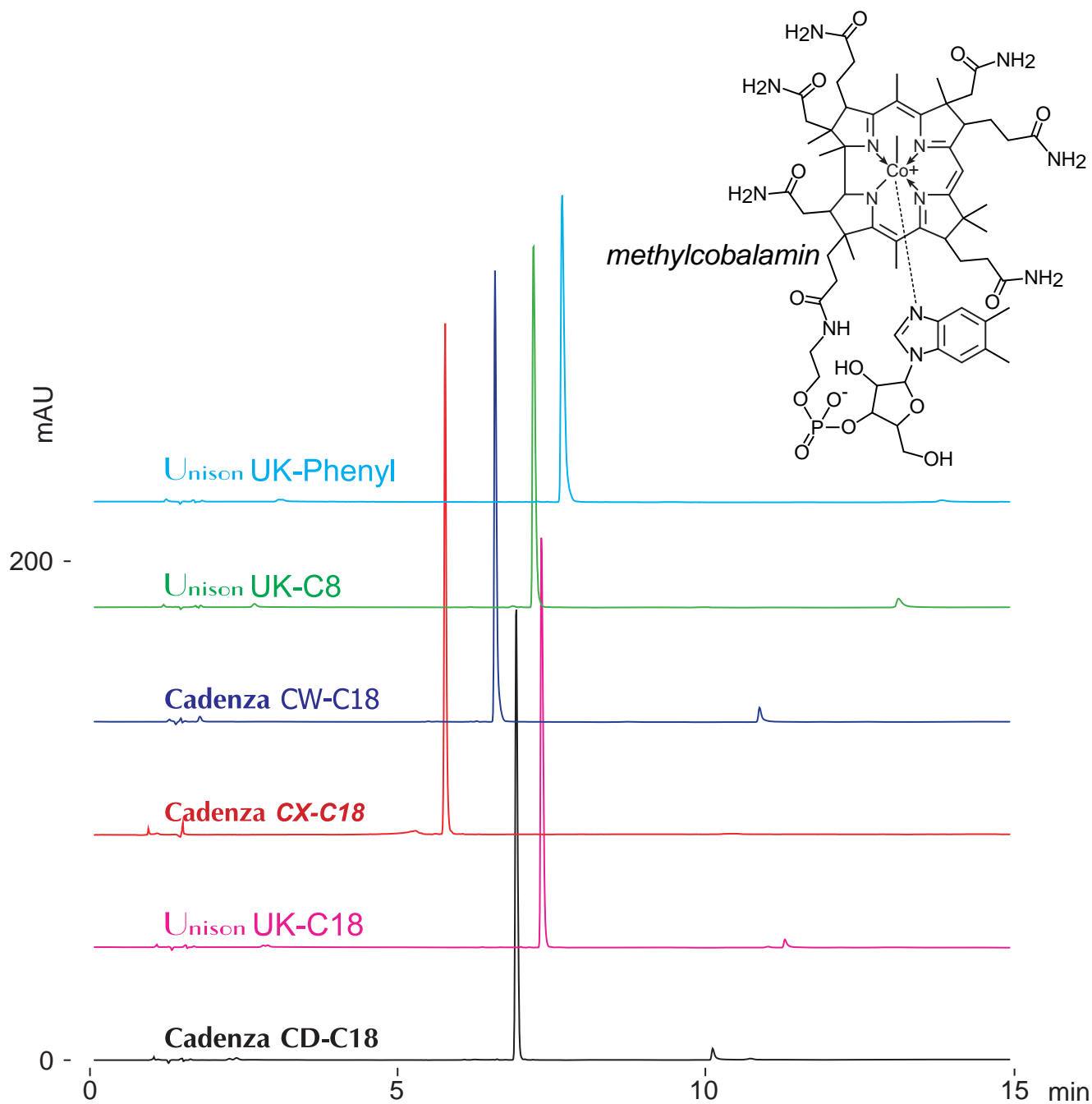

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

0-50%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.15mg/mL, 100mM ammonium formate)