

Unison UK-C18




100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Patulin

パツリンの LC-MS/MS 分析

Patulin STD 10 ppb

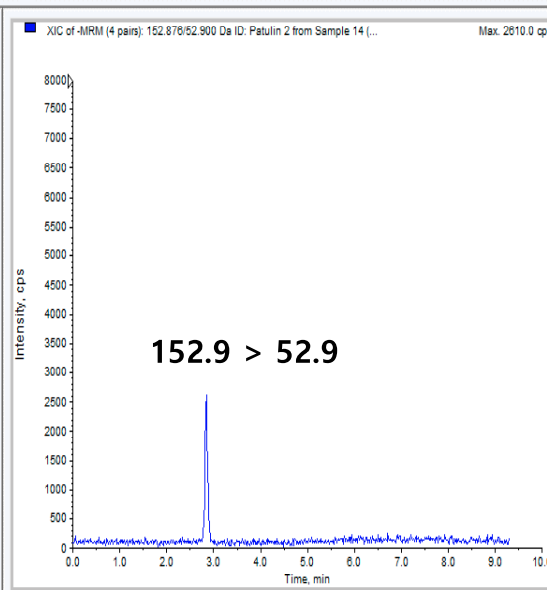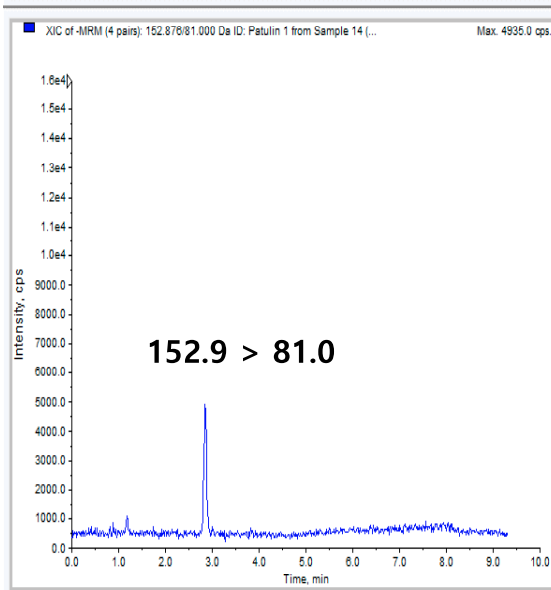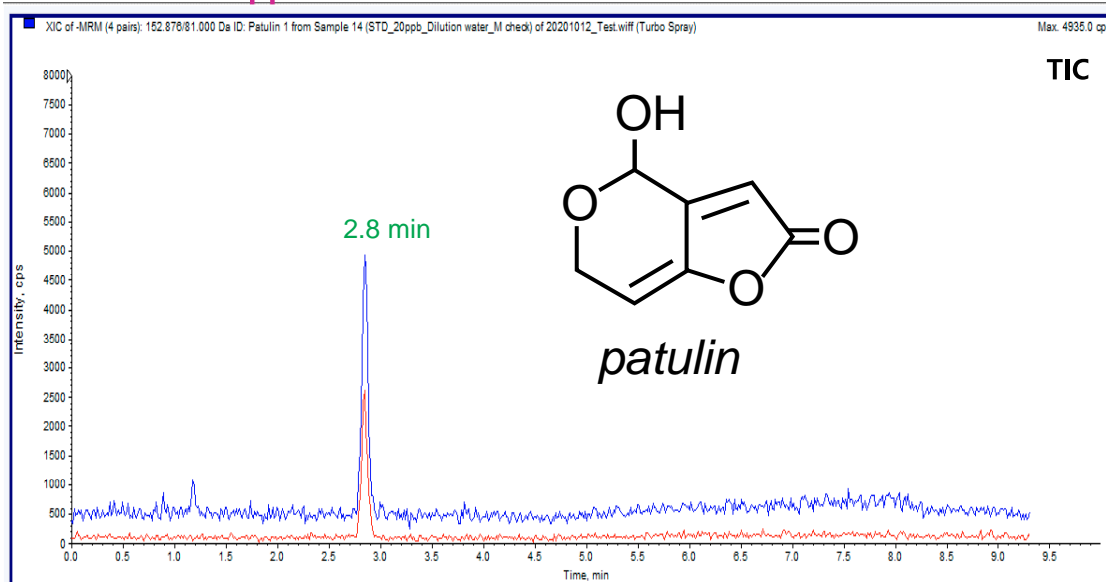

Unison UK-C18, 100 x 2 mm

A: 5mM ammonium formate

B: acetonitrile

10%B (0-1min), 10-90%B (1-4min), 90%B (4-6min),

90-10%B (6-6.1min), 10%B (6.1-10min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Negative)

5µL (10 ppb)

Courtesy of PEAKMAN SP, KOREA