


Unison UK-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Isonthiazolinone

イソチアゾリノン類のLC-MS/MS分析

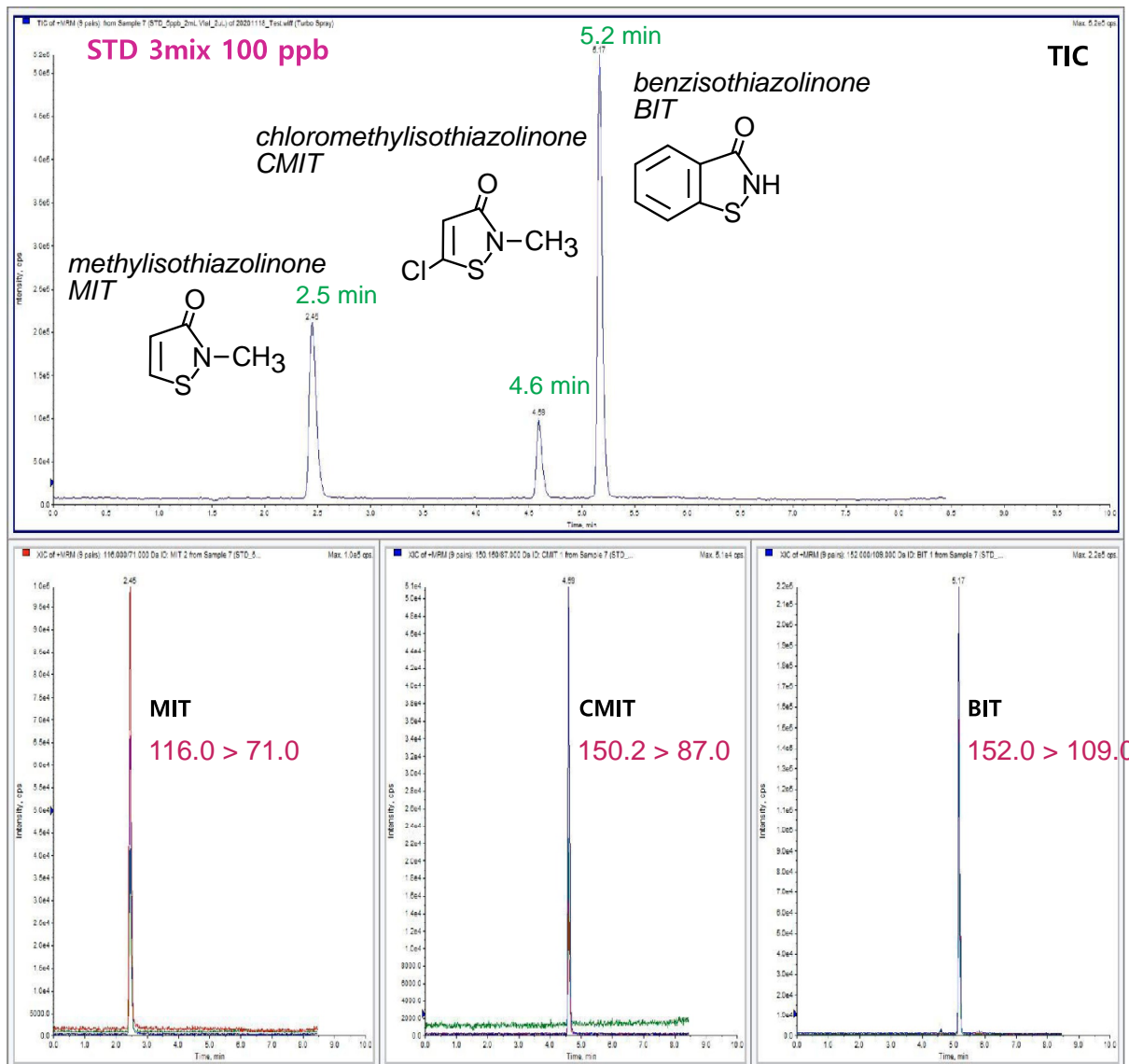

Unison UK-C18, 150 x 3 mm

A: 0.1% acetic acid in water

B: 0.1% acetic acid in methanol

30%B (0-1min), 30-90%B (1-4min), 90%B (4-6min),

90-30%B (6-6.1min), 30%B (6.1-10min)

0.4mL/min, 40°C, Triple Quad. MS (ESI, Positive)

2µL (100 ppb)

Courtesy of PEAKMAN SP, KOREA