


Unison UK-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Nitrosoamines

ニトロソアミンの LC-MS/MS 分析

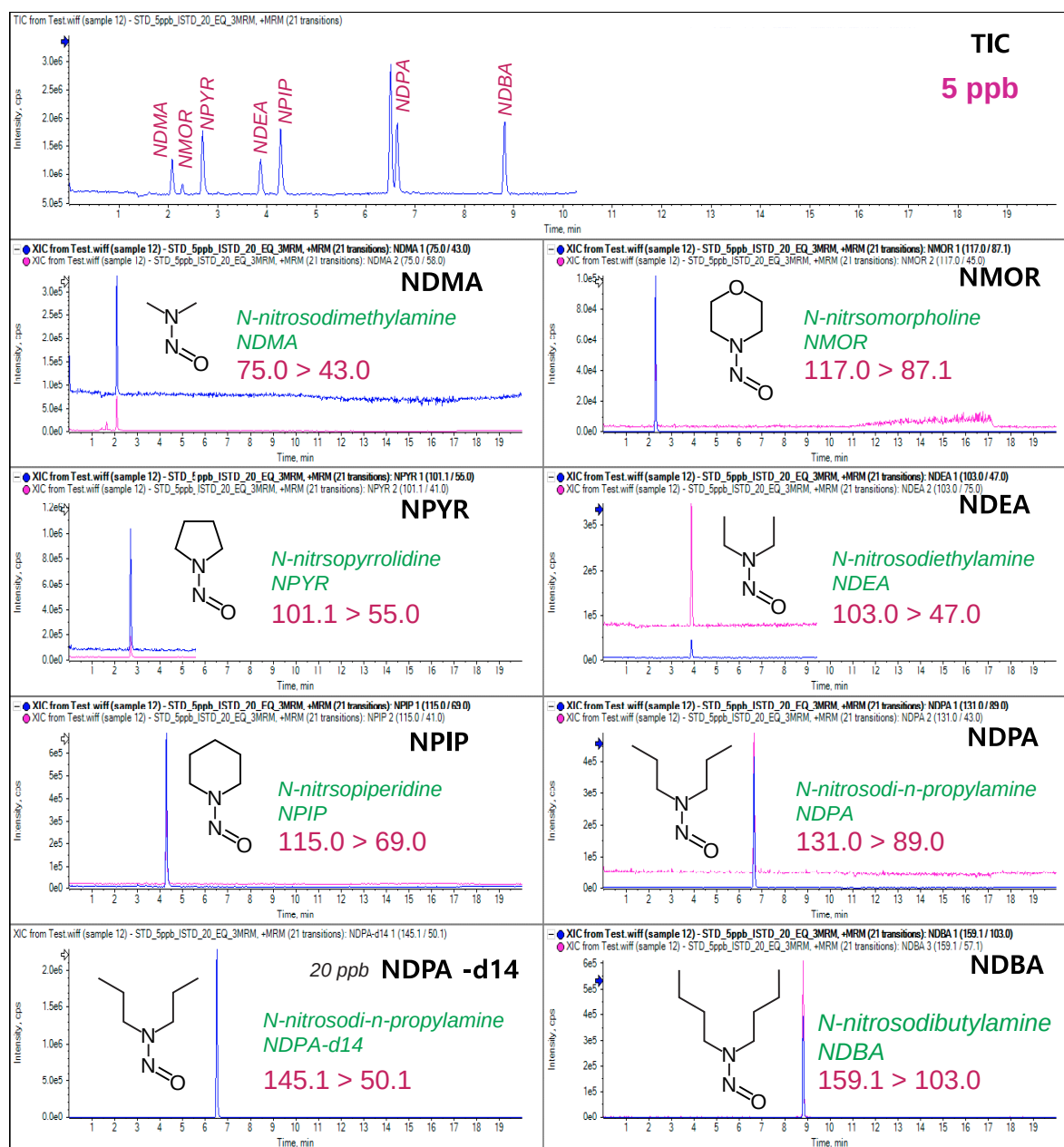

Unison UK-C18, 150 x 3 mm

A: 0.1% formic acid in water

B: 0.1% formic acid in methanol

30-100%B (0-10min), 100%B (10-15min),

100-30%B (15-15.1min), 30%B (15.1-20min)

0.5mL/min, 35°C, Triple Quad. MS (ESI, Negative)

20µL (5 ppb)

Courtesy of PEAKMAN SP, KOREA