

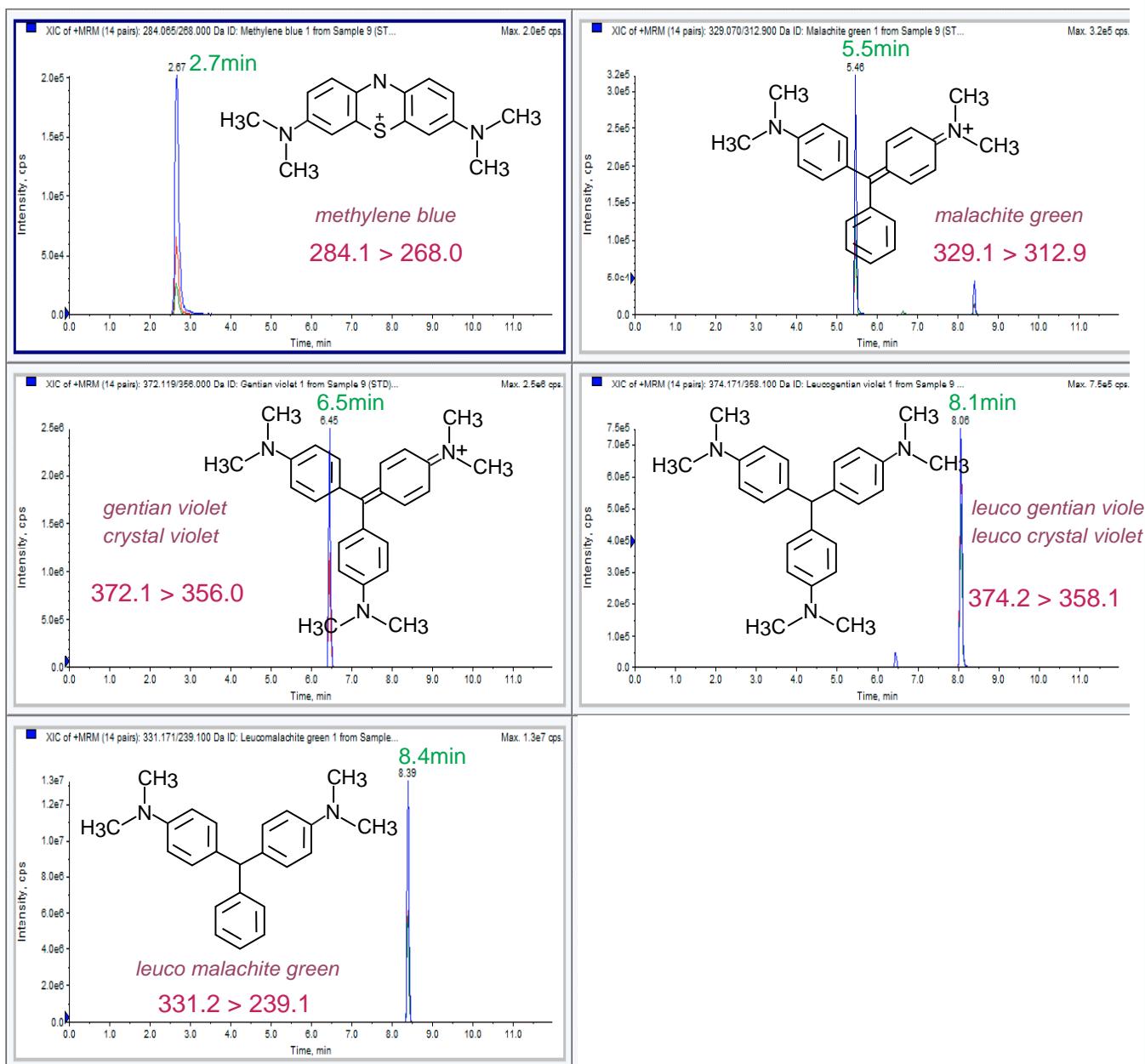
Unison UK-C18

100 x 2 mm

HPLC COLUMN APPLICATION

LC-MS/MS: Dye compounds

色素の LC-MS/MS 分析

*Unison UK-C18, 100 x 2 mm*

A: 0.1% formic acid + 10mM ammonium acetate in water

B: 0.1% formic acid in acetonitrile

30%B (0-1min), 30-90%B (1-7min), 90%B (7-9min),

90-30%B (9-9.1min), 30%B (9.1-12min)

0.3mL/min, 40°C, Triple Quad. MS (ESI, Positive)

2 μ L

Courtesy of PEAKMAN SP, KOREA