


Cadenza CD-C18 HT
150 x 2 mm
HPLC COLUMN APPLICATION

LC-MS/MS: Antibiotics (2)
  
抗生物質の LC-MS/MS 分析 (2)

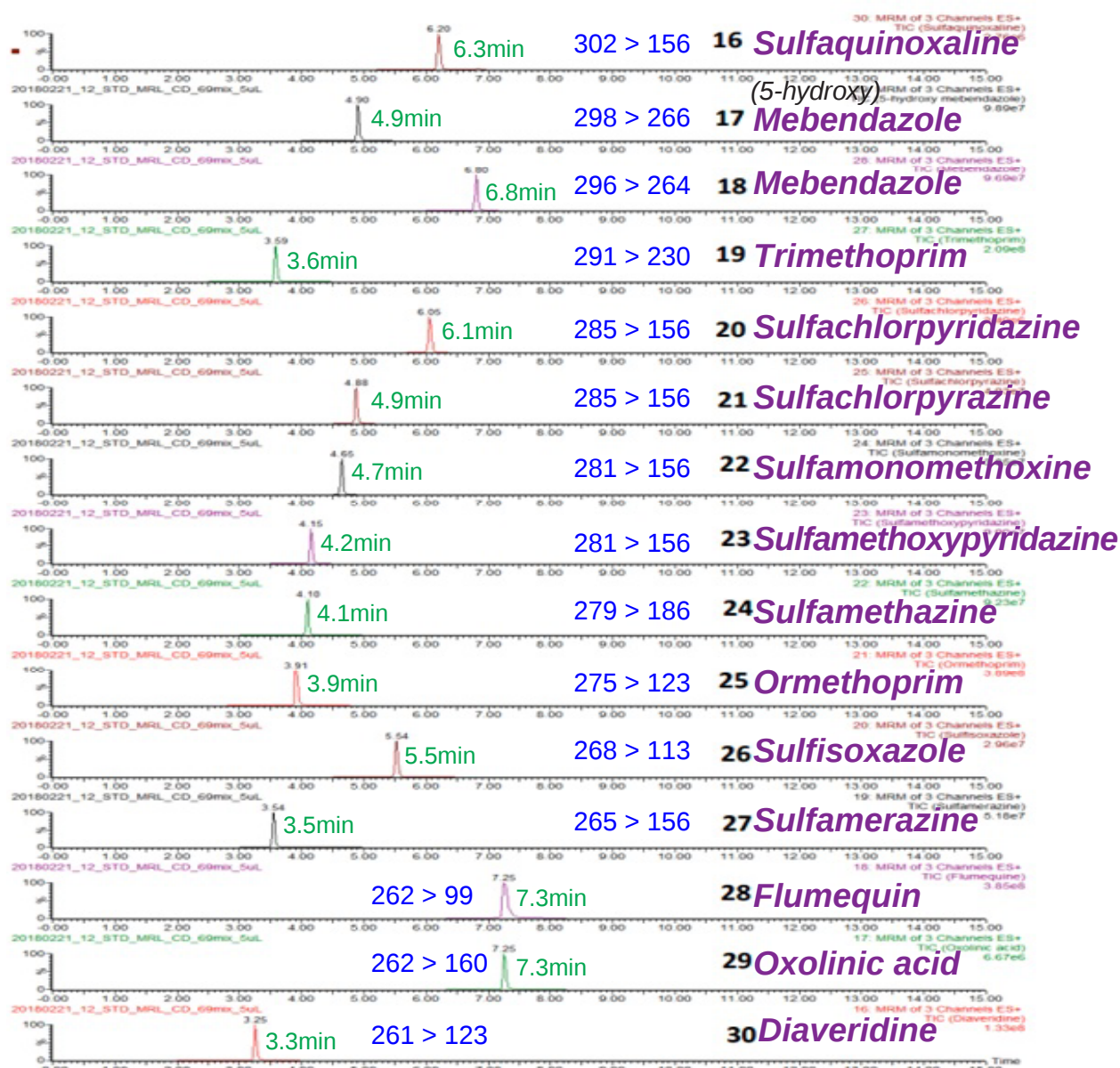

Cadenza CD-C18 HT, 150 x 2 mm  
A: 0.1% formic acid in water  
B: 0.1% formic acid + 2mM ammonium formate in acetonitrile  
10%B (0-1min), 10-40%B (1-6min), 40-95%B (6-10min),  
95%B (10-15min), 95-10%B (15-15.1min), 10%B (15.1-20min)  
0.4mL/min, 40°C, Triple Quad. MS (ESI, Positive/Negative)  
5µL

Courtesy of PEAKMAN SP, KOREA