



Nardis ND-NX

150 x 4.6 mm


HPLC COLUMN APPLICATION

Mercaptopurine (Antineoplastic)


メルカプトプリン (抗悪性腫瘍薬)

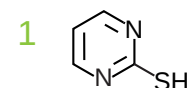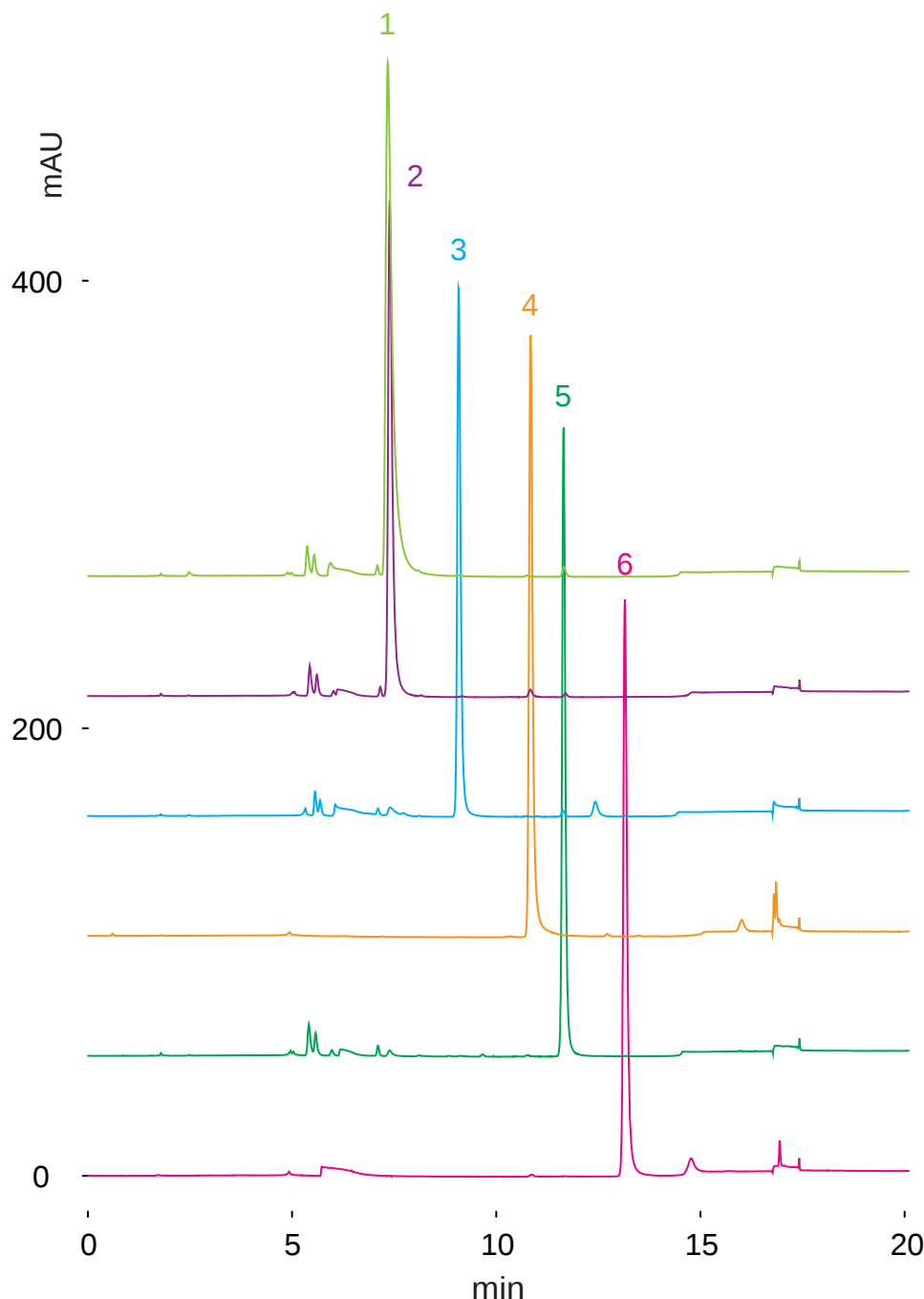


2-mercaptopyrimidine
(0.08mg/mL)

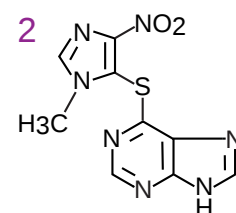


azathioprine
(0.1mg/mL)

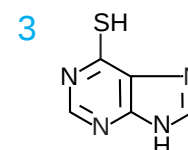


6-mercaptopurine
(0.4mg/mL)

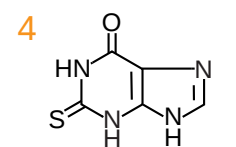

2-thioxanthine
(0.1mg/mL)

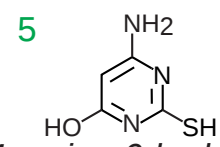

4-amino-6-hydroxy
-2-mercaptopyrimidine
(0.06mg/mL)

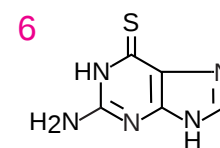

6-thioguanine
(0.5mg/mL)

Nardis ND-NX, 150 x 4.6 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium formate = 50 / 50

0 -40%B (0 -15 min), 100%B (15 -20 min)

1.0 mL/min (5MPa), 37°C, 280nm

10 µL (0.06-0.5mg/mL, diluted with water/HCOOH=100/1)