

Nardis ND-RX




Unison UK-Phenyl

250 x 4.6 mm

HPLC COLUMN APPLICATION

Cholane derivatives (Cholanic acids)

コラン誘導体（コラン酸）の分析

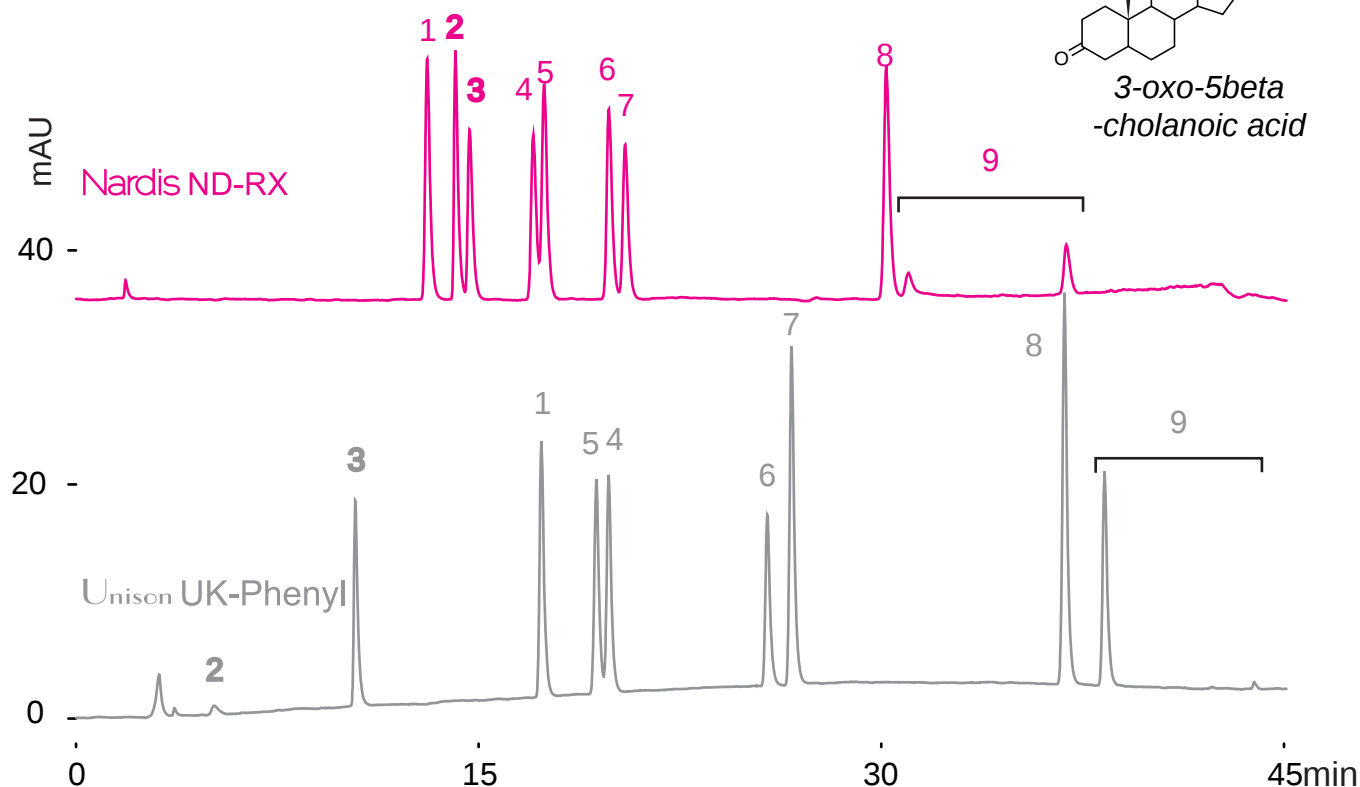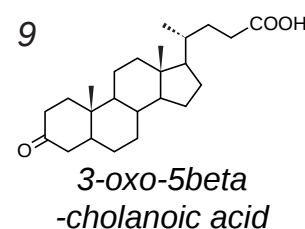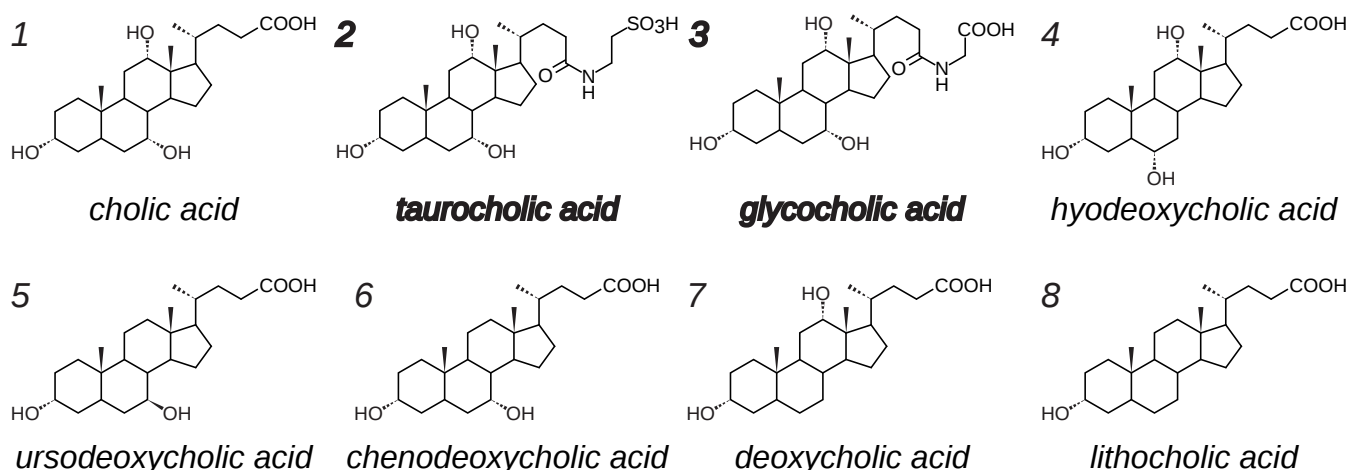

250 x 4.6 mm

A: water / acetonitrile / acetic acid = 70 / 30 / 0.5

B: 100mM ammonium acetate / acetonitrile = 50 / 50

0 -100%B (0 -40 min), 100%B (40 -45 min)

1.0 mL/min (17 MPa) , 37°C, ELSD (spray chamber 70°C, drift chamber 100°C)

10 µL (1mg/mL, methanol)