


Unison UK-Phenyl
Unison UK-C8
Nardis ND-RX
Cadenza CD-C18
Unison UK-C18

150 x 4.6 mm


HPLC COLUMN APPLICATION

Emtricitabine

エムトリシタビンの分析

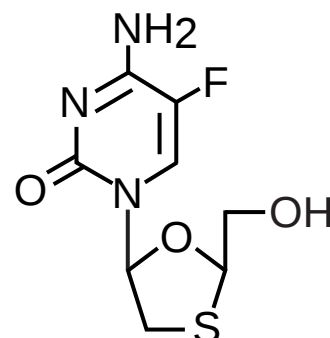

emtricitabine
(0.1mg/mL)

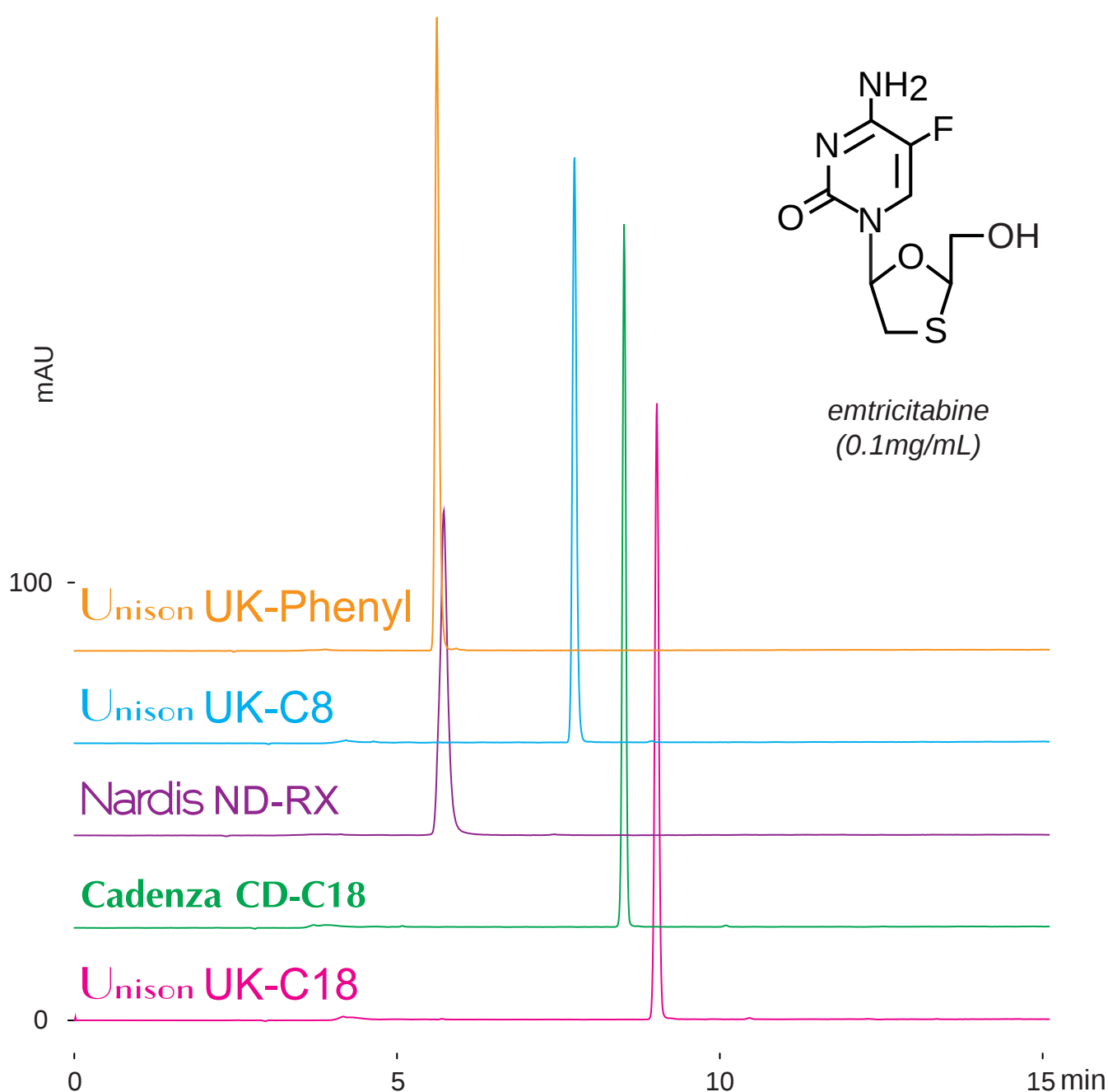

150 x 4.6 mm

A: water / acetic acid = 100 / 0.5

B: acetonitrile / acetic acid = 100 / 0.5

0-15%B (0-15min)

1.0mL/min (9-10MPa), 37°C, 260nm

10μL (0.1mg/mL, water)