


Unison UK-Phenyl
Unison UK-C8
Nardis ND-RX
Cadenza CD-C18
Unison UK-C18

150 x 4.6 mm


HPLC COLUMN APPLICATION

Brivudine

ブリブジンの分析

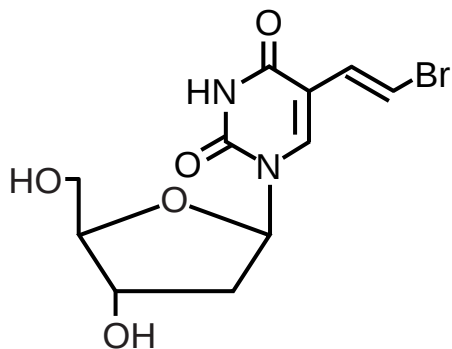

brivudine
((E)-5-(2-bromovinyl)
-2'-deoxyuridine)

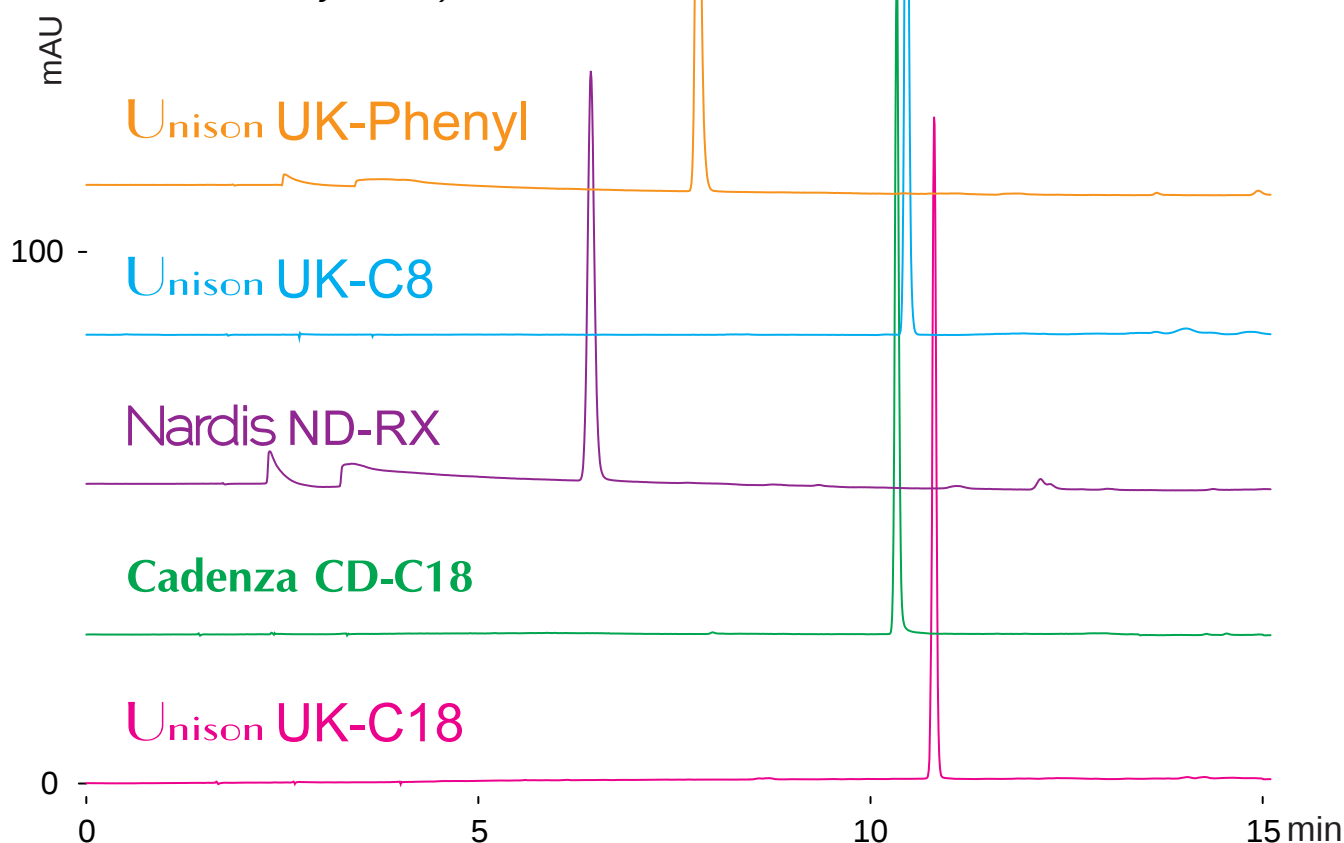

150 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-80%B (0-15min)

1.0mL/min (9-10MPa), 37°C, 260nm

10μL (50μg/mL, water/acetonitrile=50/50)