


Unison UK-Phenyl
Unison UK-C8
Nardis ND-RX
Cadenza CD-C18
Unison UK-C18

150 x 4.6 mm


HPLC COLUMN APPLICATION

Idoxuridine

イドクスウリジンの分析

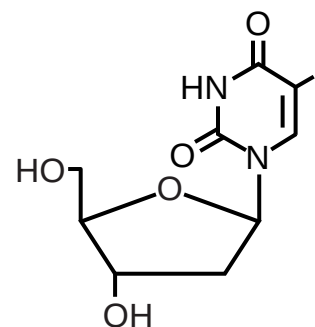

idoxuridine
(5-iodo-2'-deoxyuridine)

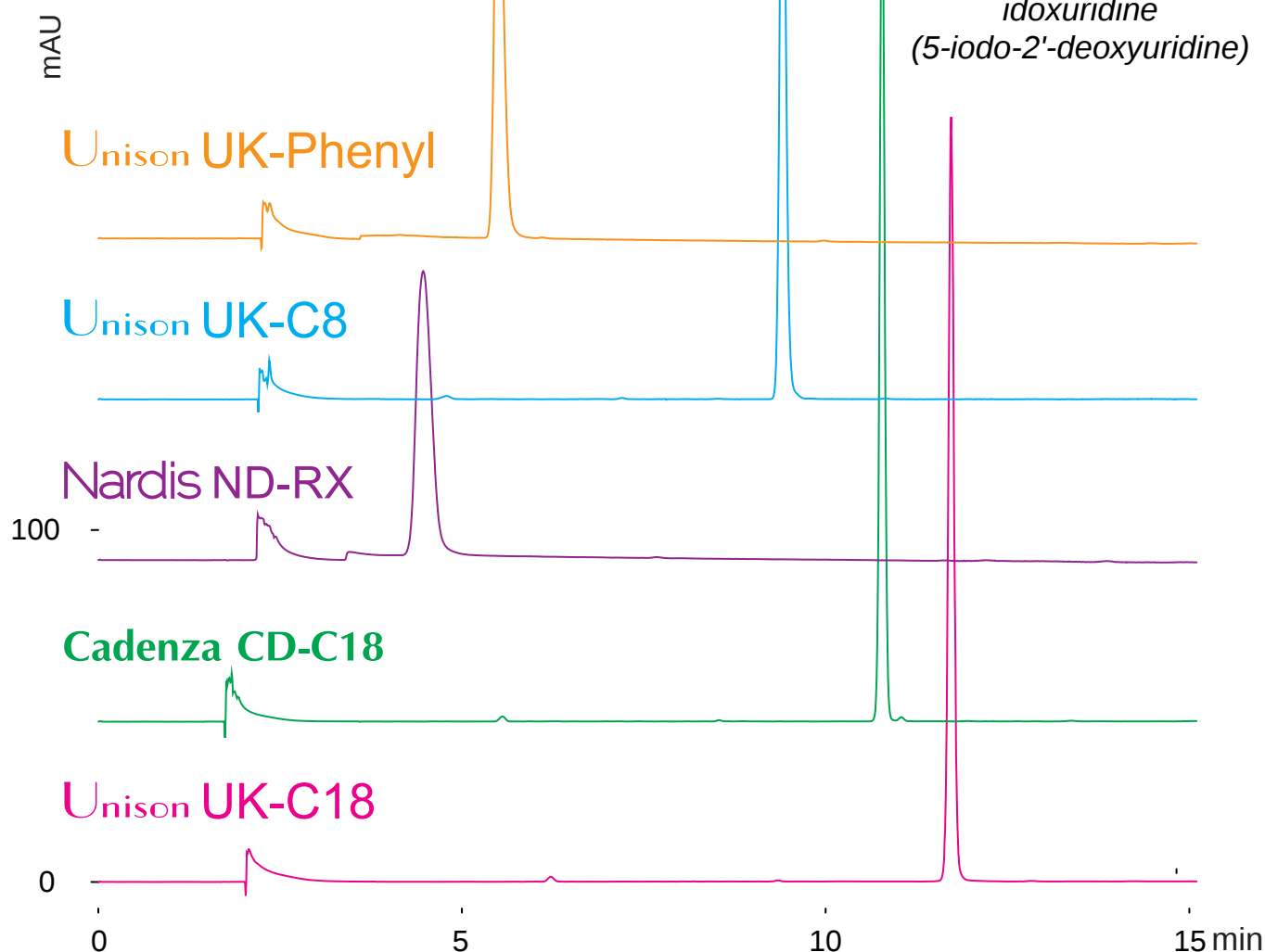

150 x 4.6 mm

A: 100mM ammonium formate

B: 100mM ammonium formate / acetonitrile = 50 / 50

0-25%B (0-15min)

1.0mL/min (9-10MPa), 37°C, 260nm

10μL (0.5mg/mL, dimethyl sulfoxide)