



Unison UK-Phenyl Cadenza CX-C18
Unison UK-C8 Cadenza CW-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Mevastatin (HMG-CoA reductase inhibitor)

メバスタチン (HMG-CoA還元酵素阻害剤)

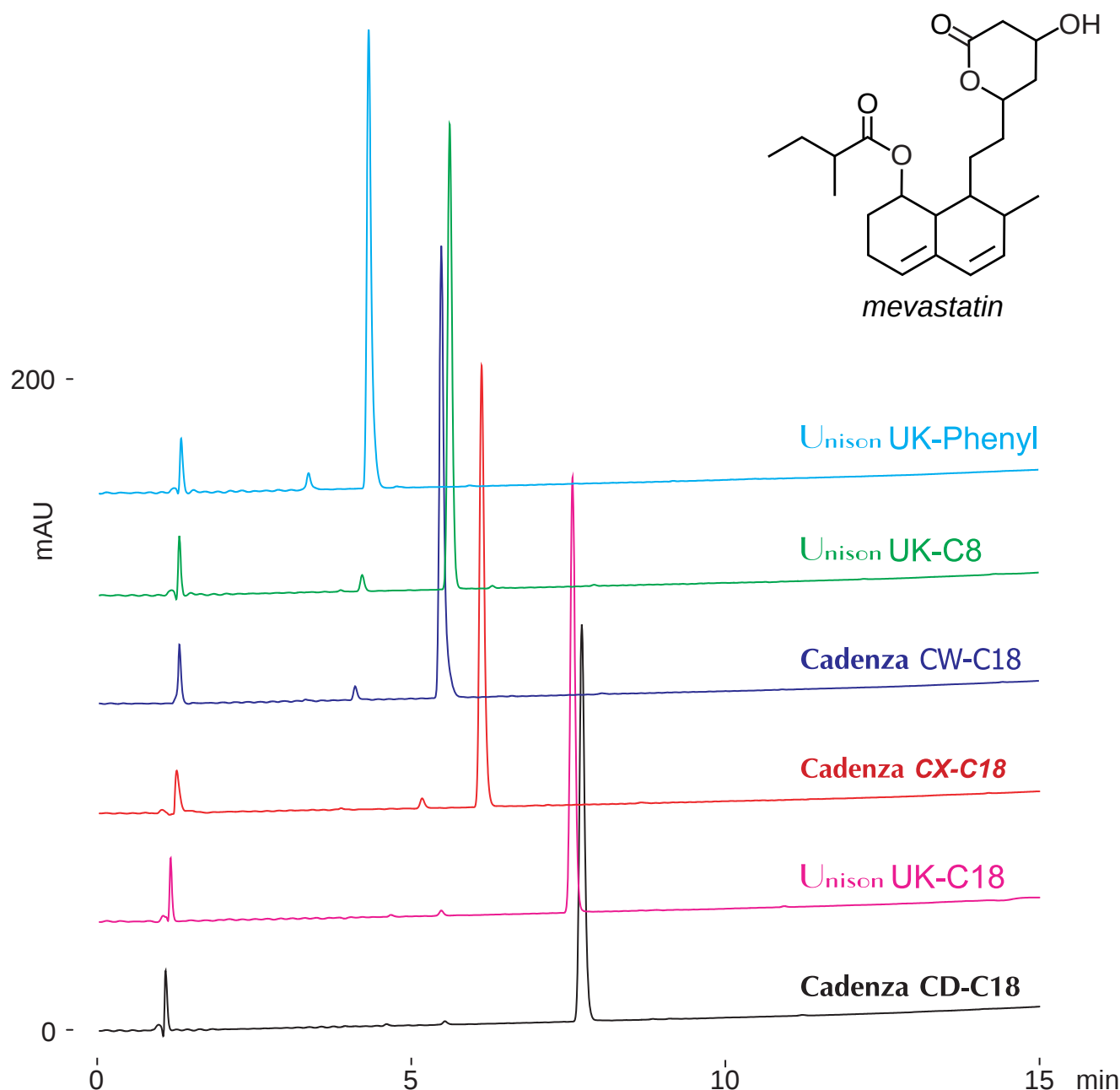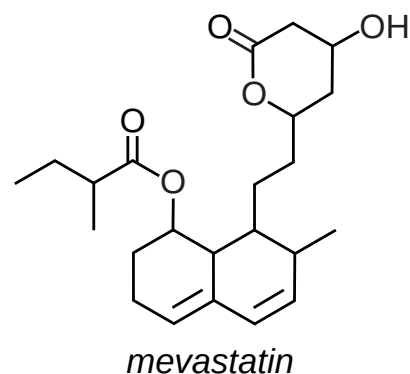

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: water / acetonitrile / formic acid = 10 / 90 / 0.1

45-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 250nm

10μL (0.15mg/mL, water / acetonitrile / formic acid = 40 / 60 / 1)