

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18

Cadenza CW-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

## Azulene sulfonate (Nonsteroidal anti-inflammatory drug)

アズレンスルホン酸 (非ステロイド性抗炎症薬)

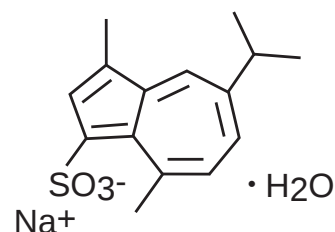

sodium azulene sulfonate hydrate

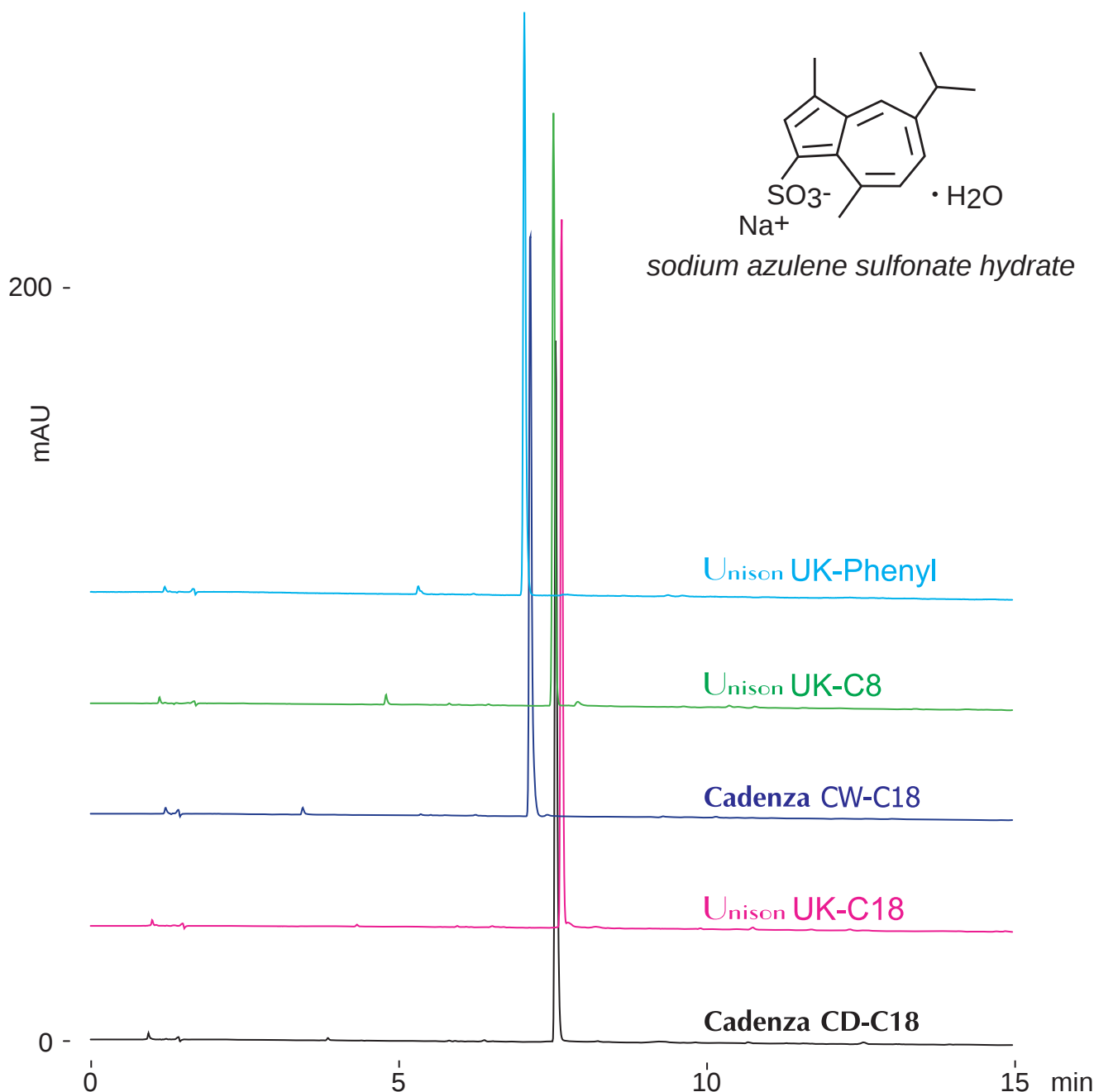

100 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 90 / 10

B: 100mM ammonium formate / acetonitrile = 10 / 90

0-85%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10μL (0.15mg/mL, 100mM ammonium formate)