


Nardis ND-NX

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Biological buffer componets

生化学用緩衝液成分の LC-MS 分析

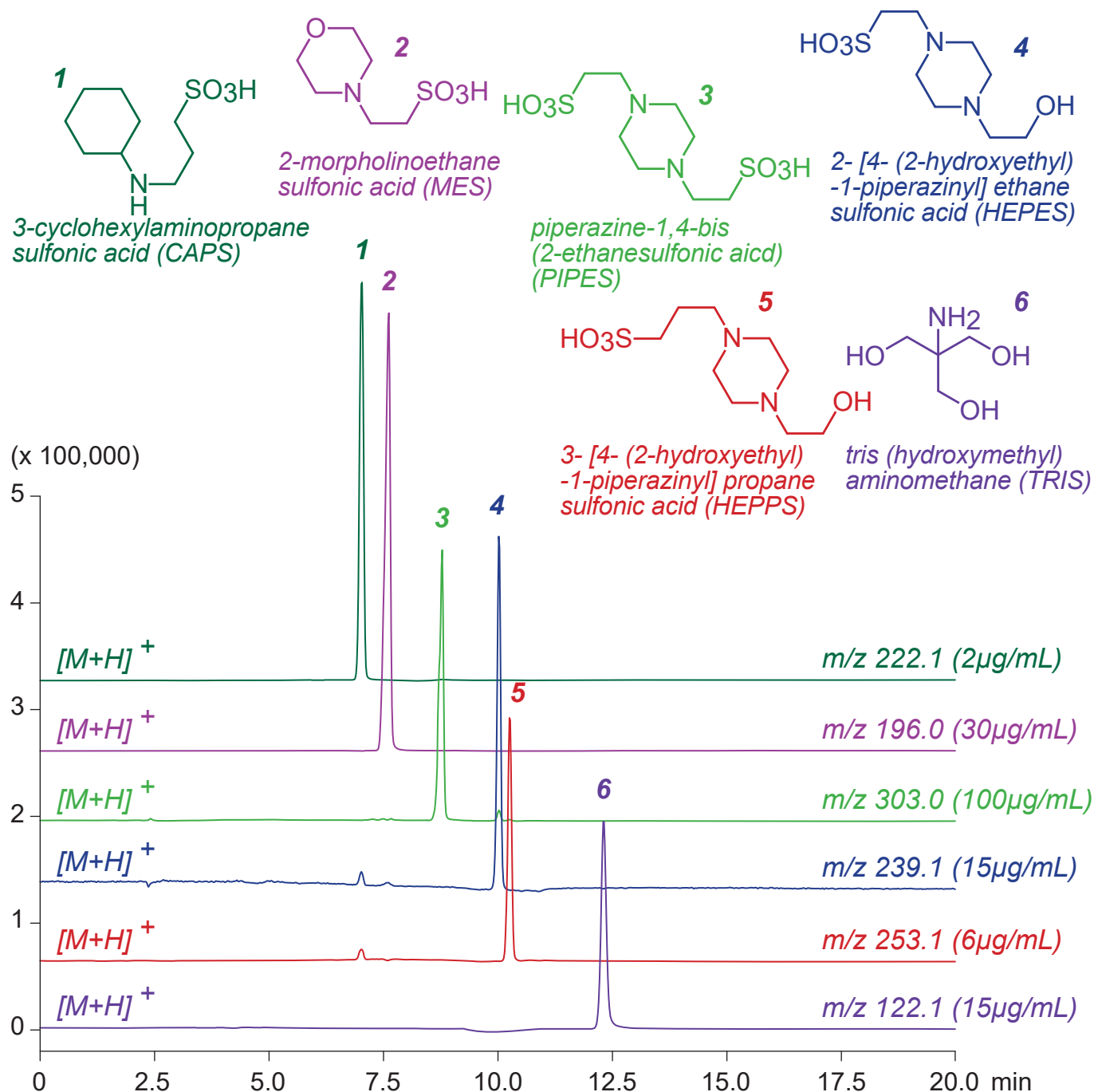

Nardis ND-NX, 150 x 3 mm

A: acetonitrile / water / formic acid = 95 / 5 / 0.1

B: acetonitrile / 100mM ammonium formate = 30 / 70

0-100%B (0-15min), 100%B (15-20min)

0.4mL/min (4MPa), 37°C, Single Quad. MS (ESI)

5 µL (2-100 µg/mL, 100mM ammonium formate)