

Unison UK-C8



Unison UK-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Organic light-emitting diode green dopant materials

有機エレクトロルミネッセンス緑色ドーパント材料

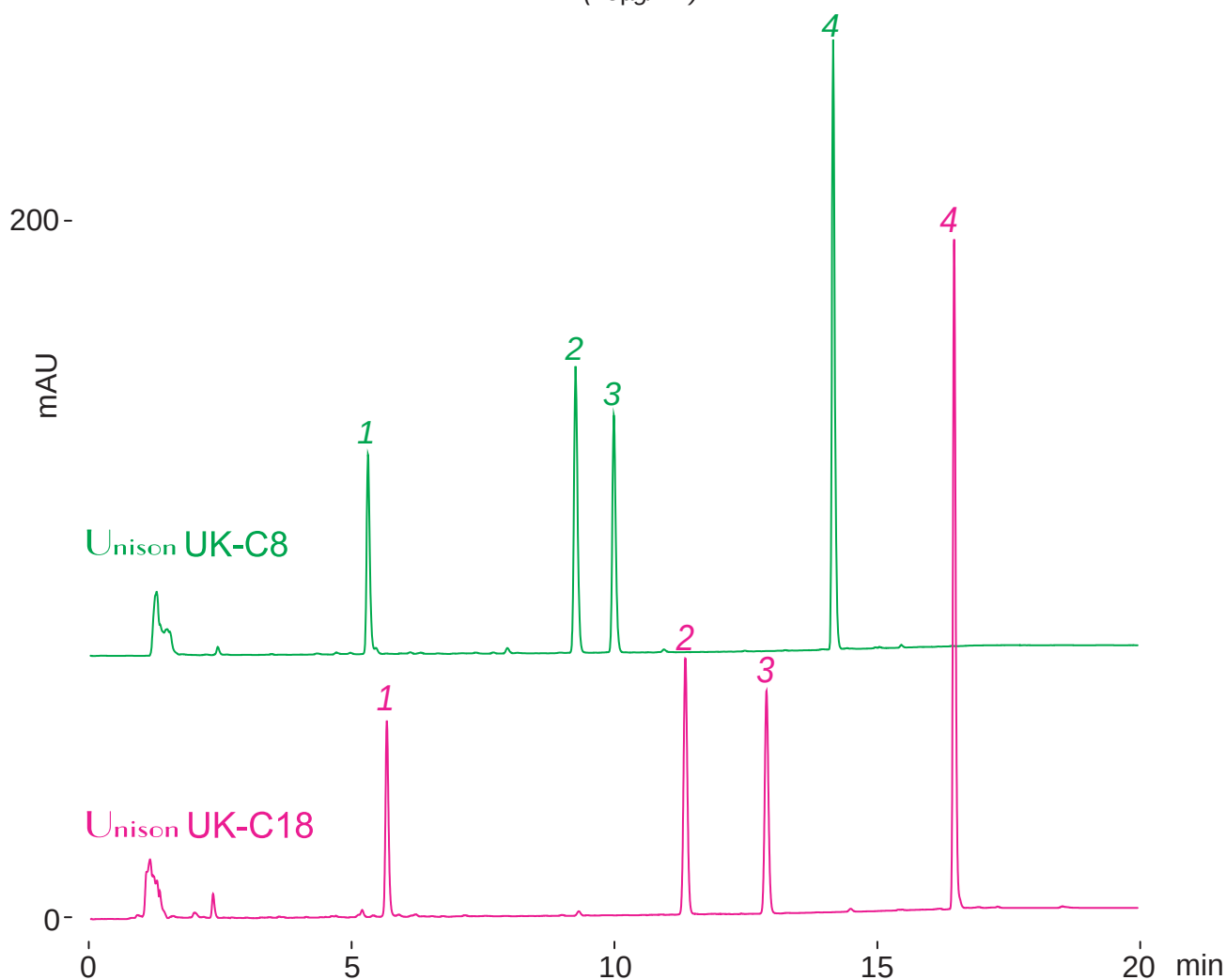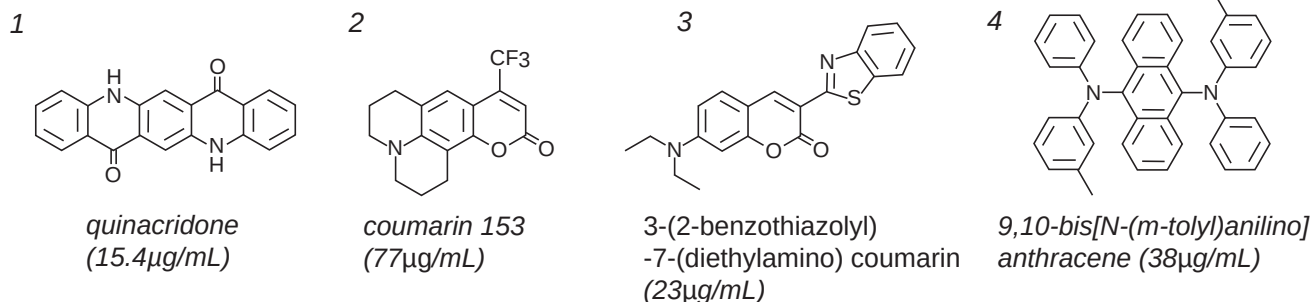

100 x 4.6 mm

A: water / acetonitrile / formic acid = 90 / 10 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

20-100%B (0-15min), 100%B (15-20min)

1.0mL/min (10MPa), 37°C, 260nm

10μL (DMSO / formic acid = 100 / 1)