


Unison UK-Phenyl Unison UK-C1
Unison UK-C8 Cadenza CX-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Tris (1,10-phenanthroline) ruthenium (II) Bis (hexafluorophosphate) (Organic light-emitting diode material)

トリス (1,10-フェナントロリン) ルテニウム (II) ビス (ヘキサフルオロホスファート)
(有機 E L 材料)

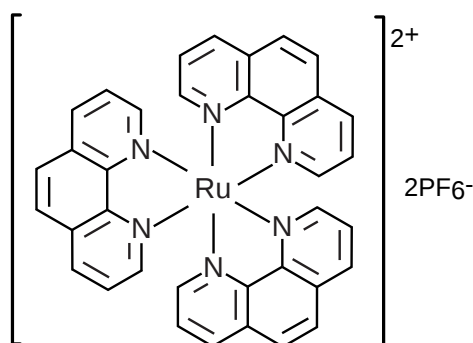

tris (1,10-phenanthroline) ruthenium (II)
bis (hexafluorophosphate)

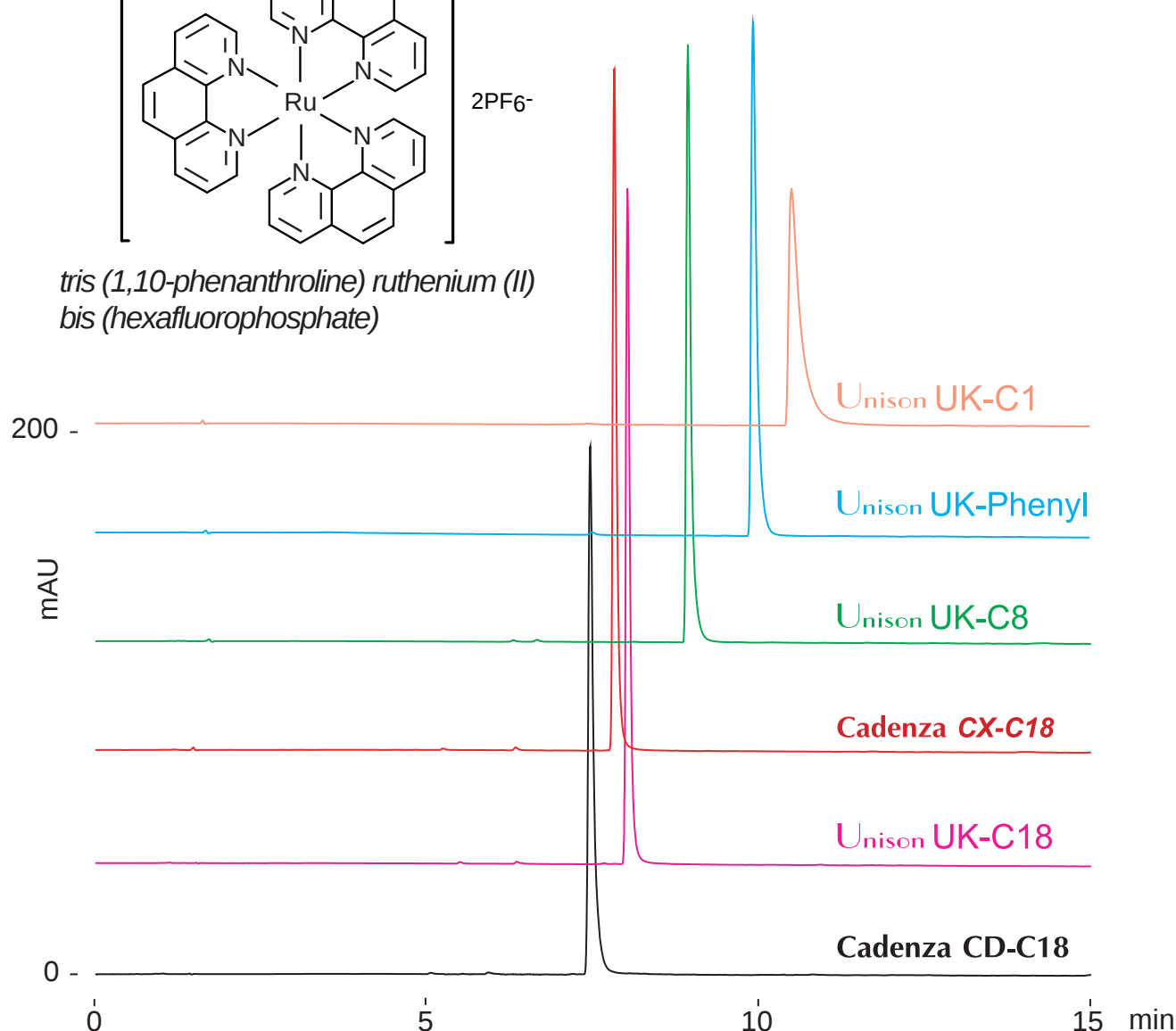

100 x 4.6 mm

A:100mM ammonium formate / acetonitrile = 90 / 10

B:100mM ammonium formate / acetonitrile = 10 / 90

0-50%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10μL (0.05mg/mL, 100mM ammonium formate)