

Cadenza CD-C18 Unison UK-C8


Unison UK-C18 Unison UK-C1

Cadenza CX-C18 Unison UK-Phenyl 150 x 4.6 mm

HPLC COLUMN APPLICATION

Hecogenin, diosgenin

ヘコゲニンとジオスゲニンの分析

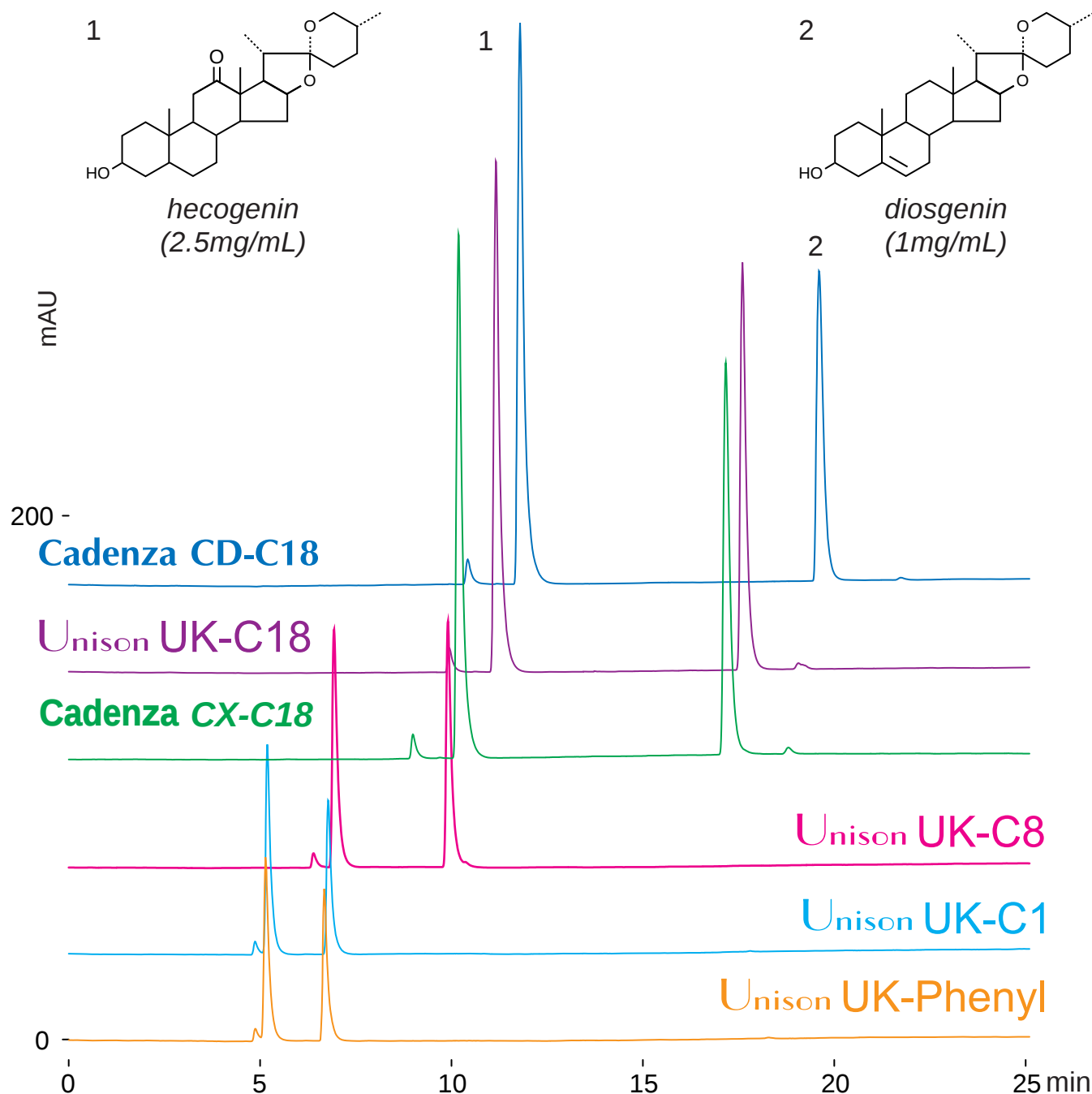

A: water / acetonitrile / acetic acid = 50 / 50 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

0-100%B (0-15min), 100%B (15-25min)

1.0mL/min (11-13MPa), 37°C

ELSD (spray chamber 70°C, drift chamber 100°C)

10μL (1-2.5mg/mL, ethanol)