


Unison UK-Phenyl Unison UK-C1
Unison UK-C8 Cadenza CX-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

N,N,N',N'-Tetrakis (p-tolyl) benzidine (Organic light-emitting diode material)

N,N,N',N'-テトラキス (p-トリル) ベンジジン
(有機EL材料)

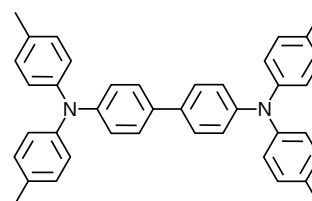

N,N,N',N'-tetrakis (p-tolyl) benzidine

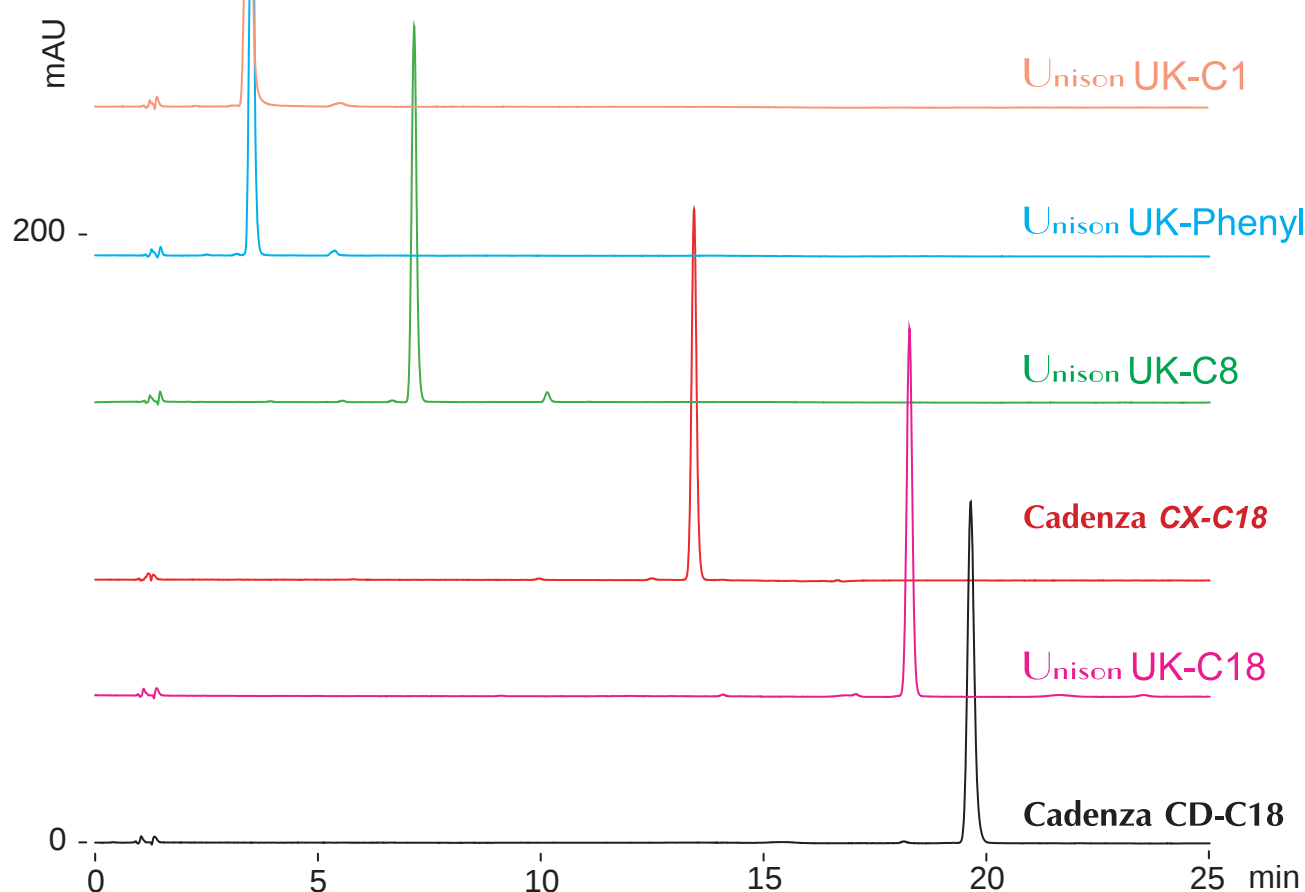

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

60-100%B (0-15min), 100%B (15-25min)

1.0mL/min (6MPa), 37°C, 350nm

10μL (0.05mg/mL, THF / ACN / HCOOH = 50 / 50 / 1)