


Unison UK-Phenyl Unison UK-C1
Unison UK-C8 Cadenza CX-C18
Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

N,N'-Bis [4- (diphenylamino) phenyl] -N,N'-diphenylbenzidine (Organic light-emitting diode material)

N,N'-ビス [4-(ジフェニルアミノ) フェニル] -N,N'-ジフェニルベンジジン
(有機EL材料)

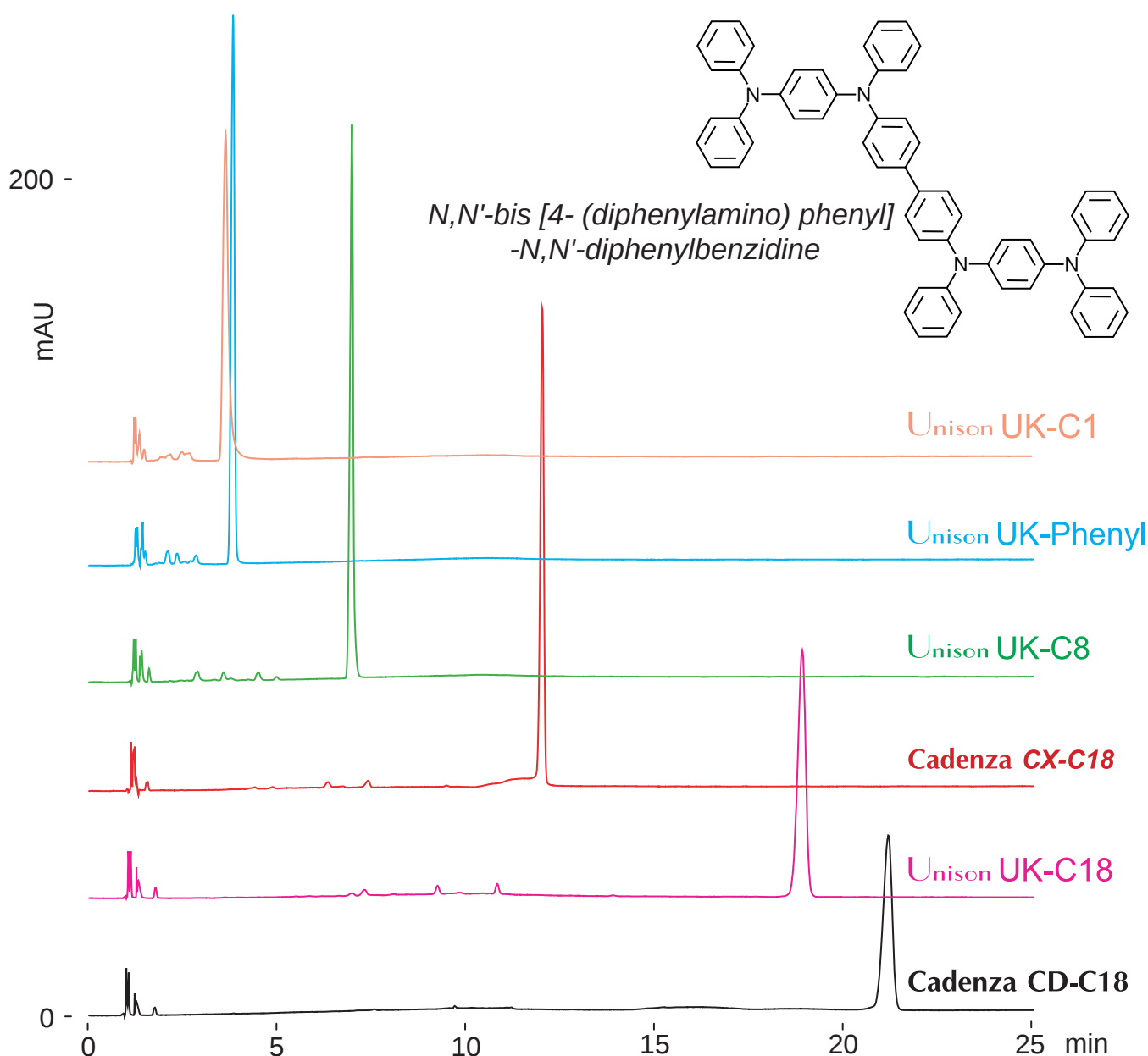

100 x 4.6 mm

A: water / acetonitrile / formic acid = 50 / 50 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

70-100%B (0-10min), 100%B (10-25min)

1.0mL/min (6MPa), 37°C, 260nm

10μL (0.2mg/mL, THF / ACN / HCOOH = 50 / 50 / 1)